


Synergy Neo 2

Multi-Mode Microplate Reader

User Manual

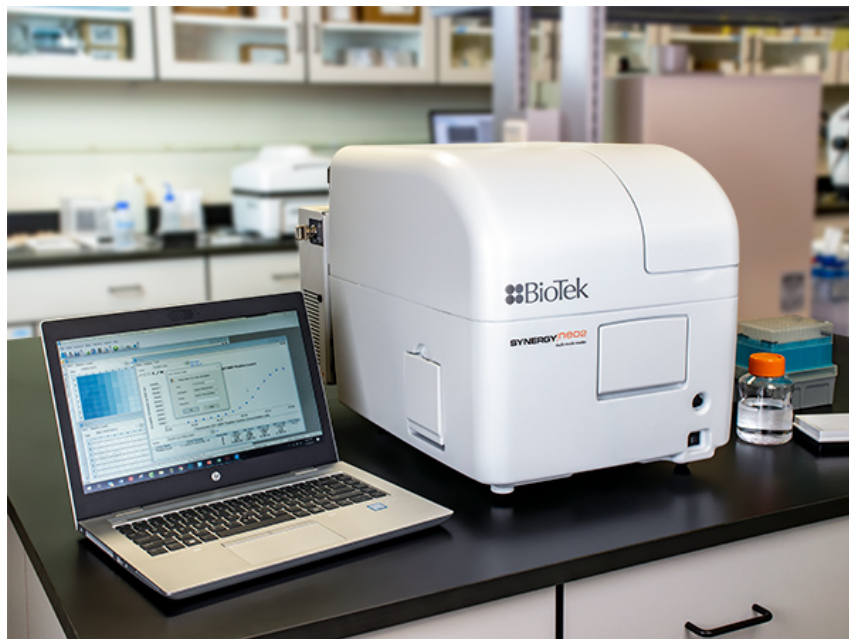

ERRATA NOTICE: This document contains references to BioTek. Please note that BioTek is now Agilent. For more information, go to www.agilent.com/lifesciences/biotek.



Synergy Neo 2

Multi-Mode Microplate Reader

User Manual

ERRATA NOTICE: This document contains references to BioTek. Please note that BioTek is now Agilent. For more information, go to www.agilent.com/lifesciences/biotek.

Document part number 1351000N, revision A

BioTek Instruments, Inc.

December 2021



Contents

Contents	ii
Notices	vi
Contact Information	vi
Intended Use Statement	vii
Quality Control	vii
Warranty and Product Registration	vii
Safety Notices	viii
Warnings and Precautions	viii
Symbols	xiii
Conformance to Standards	xxiv
EMC Information and Technical Description	xxvi
Ingress Protection Code	xxvi
Disposal	xxvi
Introduction	1
Product Description	2
Package Contents	3
Optional Accessories	4
Technical Support	7
Installation	8
Product Registration	9
Important Information	9
1: Unpack and Inspect the Reader	10
2: Unpack and Inspect the Dispenser	10
3: Unpack and Inspect the Gas Controller (if applicable)	11
4: Select an Appropriate Location	11
5: Remove the Shipping Hardware	12
6: Install the Power Supply	14
7: Install the Gas Controller (if applicable)	15
8: Install the Dispenser	15

9: Connect the Host Computer	19
10: Install Gen5 on the Host Computer	19
11: Turn on the Reader	20
12: Establish Communication	20
13: Set Dispenser Calibration Values	21
14: Run a System Test	22
15: Test the Injector System	24
Operational/Performance Qualification	25
Repackaging and Shipping Instructions	25
Getting Started	34
Modular Design	35
External Components	37
Internal Components	39
Gen5 Software	55
Recommendations for Optimum Performance	57
Periodic Maintenance	61
Periodic Maintenance Overview	62
Recommended Maintenance Schedule	62
Warnings and Precautions	63
Clean Exposed Surfaces	64
Inspect/Clean Filter Cubes	65
Flush/Purge the Fluid Path	65
Run a Dispense Protocol (Optional)	67
Empty/Clean the Tip Priming Trough	68
Clean the Priming Plate	68
Clean the Dispense Tubes and Injectors	68
As-Needed Maintenance	70
Decontamination	71
Procedure for Models without a Dispenser	72
Procedure for Models with a Dispenser	73
Dispenser Syringe Replacement	76

Instrument Qualification	79
Instrument Qualification Overview	80
IQ/OQ/PQ Description	80
Recommended Qualification Schedule	81
System Test	83
Plate Shaker Test	84
Absorbance Testing Overview	85
Absorbance Plate Test	86
Absorbance Liquid Tests	91
Fluorescence Testing	98
Luminescence Test	137
Dispense Module Tests	143
Alpha Detection Test	151
Specifications	157
General Specifications	158
Dispense/Read Specifications	159
Absorbance Specifications	161
Fluorescence Specifications (Mono-Based)	163
Fluorescence Specifications (Filter-Based)	164
Luminescence Specifications	166
Alpha Detection Specifications	166
Error Conditions	167
Error Codes Overview	168
Error Codes	168
Microplate Barcode Scanner	175
Overview	176
Install the Microplate Barcode Scanner	176
About Barcode Labels	179
Placing a Barcode Label	180
Microplate Scanner Diagnostic Utility	182

Safety Information	184
Safety Notices	185
In This Book	213

Notices

All Rights Reserved

© 2021, BioTek Instruments, Incorporated. No part of this publication may be reproduced, transcribed, or transmitted in any form, or by any means electronic or mechanical, including photocopying and recording, for any purpose other than the purchaser's use without written permission of BioTek Instruments, Inc.

Restrictions and Liabilities

Information in this document is subject to change and does not represent a commitment by BioTek Instruments, Inc. Changes made to the information in this document will be incorporated in new editions of the publication. No responsibility is assumed by BioTek for the use or reliability of software or equipment that is not supplied by BioTek or its affiliated dealers.

Contact Information



Agilent Technologies, Inc.
5301 Stevens Creek Blvd.
Santa Clara, CA 95051

Worldwide Sales and Support

www.agilent.com/en/contact-us/page

Technical Support and Service

www.agilent.com/en/support

bio.tac@agilent.com

Instrument service and repair is available worldwide at one of our international service centers and in the field at your location.

UK Responsible Person (UKRP)

Agilent LD UK Ltd
5500 Lakeside
Cheadle Royal Business Park
Cheadle, Cheshire SK8 3GR

Intended Use Statement

The Synergy Neo2 is a multi-mode microplate reader and intended to be used for the examination of specimens to analyze their characteristics and impact on a variety of analytes.

Quality Control

It is considered good laboratory practice to run laboratory samples according to instructions and specific recommendations included in the assay package insert for the test to be conducted. Failure to conduct Quality Control checks could result in erroneous test data.

Warranty and Product Registration

Review the warranty information that shipped with your product. Register to ensure you receive important information updates about the products you have purchased.

Safety Notices

Raadpleeg Bijlage D voor informatie in andere talen.

Reportez-vous à l'annexe D pour obtenir des informations dans d'autres langues.

Informationen in anderen Sprachen finden Sie in Anhang D.

Fare riferimento all'Appendice D per informazioni in altre lingue.

Consulte el Apéndice D para obtener información en otros idiomas.

Pay special attention to the following safety notices in all product documentation.

WARNING A WARNING notice denotes a hazard. It calls attention to an operating procedure, practice, or the like that, if not correctly performed or adhered to, could result in personal injury or death. Do not proceed beyond a WARNING notice until the indicated conditions are fully understood and met.

CAUTION A CAUTION notice denotes a hazard. It calls attention to an operating procedure, practice, or the like that, if not correctly performed or adhered to, could result in damage to the product or loss of important data. Do not proceed beyond a CAUTION notice until the indicated conditions are fully understood and met.

Warnings and Precautions

Electrical Hazards

WARNING **Internal Voltage.** Always turn off the power switch and unplug the power supply before cleaning the outer surface of the instrument.

WARNING **Power Rating.** The instrument's power supply or power cord must be connected to a power receptacle that provides voltage and current within the specified rating for the system. Use of an incompatible power receptacle may produce electrical shock and fire hazards.

WARNING **Electrical Grounding.** Never use a plug adapter to connect primary power to the external power supply. Use of an adapter disconnects the utility ground, creating a severe shock hazard. Always connect the power cord directly to an appropriate receptacle with a functional ground.

WARNING **Service.** Only qualified technical personnel should perform service procedures on internal components.

CAUTION **Power Supply.** Use only the power supply shipped with the instrument, and operate it within the range of line voltages listed on it.

Chemical/Environmental

WARNING **Potential Biohazards.** Some assays or specimens may pose a biohazard. Adequate safety precautions should be taken as outlined in the assay's package insert. Always wear safety glasses and appropriate protective equipment, such as chemical-resistant rubber gloves and apron.



WARNING **Liquids.** Avoid spilling liquids on the instrument; fluid seepage into internal components creates a potential for shock hazard or instrument damage. If a spill occurs while a program is running, stop the program and turn off the instrument. Wipe up all spills immediately. Do not operate the instrument if internal components have been exposed to fluid.

CAUTION **Liquids.** Do not immerse the instrument, spray it with liquid, or use a dripping-wet cloth on it. Do not allow water or other cleaning solution to run into the interior of the instrument. If this happens, contact Technical Support.

CAUTION **Environmental Conditions.** Do not expose the instrument to temperature extremes. For proper operation, temperature near the instrument should remain within the range in the **Specifications** section of this document. Performance may be adversely affected if temperatures fluctuate above or below this range

CAUTION **Sodium Hypochlorite.** Do not expose any part of the instrument to the recommended diluted sodium hypochlorite solution for more than 20 minutes. Prolonged contact may damage the instrument surfaces. Be certain to rinse and thoroughly wipe all surfaces.

CAUTION **Lubricants.** Do not apply lubricants to moving parts. Lubricant on components in the carrier compartment will attract dust and other particles, which may cause the instrument to produce an error.


Components

WARNING



Laser Beam. Serious eye injury may occur if you stare directly into the laser beam of the barcode scanner during operation of the scanner. This hazard is noted by the symbol shown here. Do not look directly into the laser beam during operation of the scanner.

WARNING

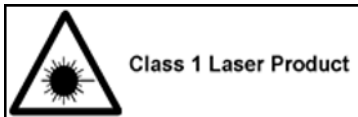

Class 1 Laser Product. Alpha laser “A” and TRF laser “T” models

See Report on laser safety on page xi.

WARNING



Pinch Hazard. Some areas of the external dispense module can present pinch hazards when the instrument is operating. Keep hands and fingers clear of these areas when the instrument is operating.

WARNING



Hot Surface. Models with incubation: The microplate carrier is hot when the incubator is turned on and for a period of time after the incubator is turned off. Do not touch the carrier or the inside surface of the door.

WARNING



Two-person lift. The instrument should be lifted by two people. The instrument with all available modules weighs up to 45.5 kg.

WARNING

Accessories. Only accessories that meet the manufacturer's specifications shall be used with the instrument.

CAUTION

Shipping Hardware. All shipping hardware must be removed before operating the instrument and reinstalled before repackaging the instrument for shipment.

CAUTION

Spare Parts. Only approved spare parts should be used for maintenance. The use of unapproved spare parts and accessories may result in a loss of warranty and potentially impair instrument performance or cause damage to the instrument.

CAUTION

Service. Only qualified technical personnel should perform service procedures on internal components.

Intended Product Use**WARNING**

Software Quality Control. The operator must follow the manufacturer's assay package insert when modifying software parameters and establishing reading methods. It is considered good laboratory practice to run laboratory samples according to instructions and specific recommendations included in the assay package insert for the test to be conducted. Failure to conduct quality control checks could result in erroneous test data.

WARNING

Data Reduction. No limits are applied to the raw measurement data. Data exported via computer control must be analyzed by the operator. The performance characteristics of the data reduction software have not been established with any laboratory diagnostic assay. Users must evaluate this instrument and PC-based software in conjunction with their specific assay (s). This evaluation must include the confirmation that performance characteristics for the specific assay(s) are met.

WARNING

Unspecified Use. Failure to operate equipment according to the guidelines and safeguards specified in the product user documentation could result in a hazardous condition.

CAUTION

Use of labware other than described in this document can result in positioning errors during program execution.

Report on laser safety**Synergy Neo2 T-models with the N2 laser**

Synergy Neo2 "T" models include a nitrogen laser MNL-100 manufactured by LTB. This laser is classified to class 3B and therefore the safety aspects are made according to laser class 3B.

Synergy Neo2 A-models with the Alpha Laser

Synergy Neo2 "A" models include an Alpha Laser, K6880E09FN-0.800W manufactured by BWT Beijing. This laser is classified to class 3B and therefore the safety aspects are made according to laser class.

Concerning class 3B lasers, the International Standard IEC 60825-1 on the safety of laser products states the following: “Class 3B: Lasers are normally hazardous when direct intrabeam exposure occurs (i.e., within the NOHD). Viewing diffuse reflections is normally safe.”

Regarding laser safety, it is extremely important to prevent the user from being exposed to laser beam, either directly or through reflection from a mirror surface. Although the laser used in Synergy Neo2 is a class-3B laser, the Synergy Neo2 instrument is a class 1 laser product. This is possible, because Synergy Neo2 has adequate laser shielding which prevents the user being exposed to the laser radiation.

The Synergy Neo2's interlock system prevents user exposure to the laser beam by disabling the laser whenever the top-filter door, plate carrier, or side (bottom-filter) hatch of the instrument is open. The state of the interlock is made visible to users via LEDs on the laser unit attached to the rear of the instrument. Red LEDs indicate interlock activity; green LEDs indicate safe laser usage.

When closed the three doors form a light-tight system. Light tightness is essential not only for laser safety, but also for the functional performance of the instrument. The laser itself has been installed on the instrument under a protective casing, and the laser beam is led by an optical fiber via the electronics compartment to the light-tight instrument compartment.

Symbols



Veiligheidssymbolen



Symboles de sécurité

Sicherheitssymbole

Simboli di sicurezza

Símbolos de seguridad

	<p>Caution, consult the instructions for use for important cautionary information such as warnings and precautions</p> <p>Voorzichtig, raadpleeg de gebruiksaanwijzing voor belangrijke voorzorgsinformatie zoals waarschuwingen en voorzorgsmaatregelen</p> <p>Attention, pour des informations de mise en garde importantes telles que des avertissements et des précautions, consultez le mode d'emploi.</p> <p>Achtung, lesen Sie die Gebrauchsanweisung für wichtige Vorsichtshinweise wie Warnungen und Sicherheitsvorkehrungen</p> <p>Attenzione, consultare le istruzioni per l'uso per importanti informazioni cautelative come avvertenze e precauzioni</p> <p>Precaución, consulte las instrucciones de uso para obtener información importante, como advertencias y precauciones</p>
	<p>Warning; Biological hazard</p> <p>Waarschuwing; biologisch gevaar</p> <p>Avertissement : Risque biologique</p> <p>Warnung; biologische Gefahr</p> <p>Avvertenza, rischio biologico</p> <p>Advertencia: peligro biológico</p>

	<p>Warning; Laser beam (If equipped with the optional external barcode scanner)</p> <p>Waarschuwing; Laserstraal (Als ik Uitgerust met de optionele externe barcode scanner)</p> <p>Avertissement; faisceau laser (Si j'EQUIPEE avec le scanner de codes à barres externe en option)</p> <p>Warnung; Laserstrahl (Wenn ich mit den optionalen externen Barcode-Scanner ausgerüstet)</p> <p>Avvertimento; Raggio laser (Se equipaggiata con il barcode scanner esterno opzionale)</p> <p>Advertencia; rayo láser (Si equipado con el escáner de código de barras externo opcional)</p>
	<p>Warning; Class 1 Laser Product (Alpha laser "A" and TRF laser "T" models)</p> <p>Waarschuwing; Klasse 1 Laser Product (Alpha Laser "A" en TRF laser "R" -modellen)</p> <p>Avertissement; Laser de classe 1 (Alpha Laser "A" et TRF laser modèles "R")</p> <p>Warnung; Laser Klasse 1 (Alpha Laser "A" und TRF-Laser "R" -Modelle)</p> <p>Avvertimento; Classe prodotto 1 Laser (Alpha Laser "A" e la FR laser modelli "R")</p> <p>Advertencia; Producto láser de clase 1 (Alpha Laser "A" y TRF láser modelos "R")</p>



Complies with FDA performance standards for laser products except for conformance with IEC 60825-1 Ed.3., as described in Laser Notice No.56, dated May 8, 2019. (Alpha laser "A" and TRF laser "T" models)

Voldoet aan FDA-prestatienormen voor laserproducten behalve conform IEC 60825-1 Ed.3., Van beschreven in Laser Notice No.56, gedateerd 8 mei 2019. (Alpha Laser "A" en TRF laser "R" -modellen)

Conforme aux normes de performance de la FDA pour les produits laser sauf pour la conformité à la norme IEC 60825-1. Ed.3, De décrits dans Laser Notice No.56, datée du 8 mai 2019. (Alpha Laser "A" modèles et TRF laser "R")

Erfüllt Performance FDA-Standards für Laserprodukte außer für die Konformität mit IEC 60825-1 Ed.3., In Laser Notice No.56 beschrieben Aus, vom 8. Mai 2019 (Alpha Laser "A" und TRF-Laser "R" -Modelle)

Conforme agli standard di prestazione FDA per i prodotti laser tranne che per la conformità alla norma IEC 60825-1. Ed.3, Da descritti nella Laser Notice No.56, in data 8 maggio 2019. (e laser modelli Alpha Laser "A" della Fondazione "R")

Conforme con las normas de rendimiento de la FDA para productos láser, excepto para la conformidad con la norma IEC 60825-1. Ed.3, desde que se describen en Laser Notice No.56, de fecha 8 de mayo de 2019. (láser y modelos Alfa láser "A" TRF "R")



Warning; Pinch hazard



Waarschuwing; beknellingsgevaar

Avertissement : risque de pincement

Warnung; Quetschgefahr

Avvertenza, rischio di pizzicamento

Advertencia: peligro de atrapamiento

 <div data-bbox="440 226 667 426" style="background-color: black; color: white; padding: 5px;"> <p>⚠ Caution</p> <p>Carton Exceeds 50lbs (22.5kg) When handling, two or more people are required</p> </div>	<p>Caution; carton exceeds 50lbs (22.5kg). When handling, two or more people are required.</p> <p>Voorzichtig; de doos weegt meer dan 22,5 kg. Bij het hanteren zijn twee of meer personen nodig.</p> <p>Attention : le carton dépasse 50 lb (22,5 kg). Lors de la manipulation, deux personnes ou plus sont nécessaires.</p> <p>Vorsicht; Karton überschreitet 22,5 kg (50 lbs). Zur Handhabung sind zwei oder mehr Personen erforderlich.</p> <p>Attenzione, la confezione pesa più di 22,5 kg. Per la movimentazione sono necessarie due o più persone.</p> <p>Precaución: la caja pesa más de 22,5 kg (50 lb). Para manipularla hacen falta dos o más personas.</p>
	<p>Warning; Hot surface</p> <p>Waarschuwing; heet oppervlak</p> <p>Avertissement : surface chaude</p> <p>Warnung; heiße Oberfläche</p> <p>Avvertenza, superficie molto calda</p> <p>Advertencia: superficie caliente</p>



Disposal Notice: Dispose of the instrument according to Directive 2012/19/EU, “on waste electrical and electronic equipment (WEEE)” or local ordinances



Kennisgeving van verwijdering: Verwijder het instrument volgens Richtlijn 2012/19/EU betreffende afgedankte elektrische en elektronische apparatuur (AEEA) of lokale verordeningen


Avis concernant la mise au rebut : mettez l'instrument au rebut conformément à la directive 2012/19/EU portant sur les déchets d'équipement électrique et électronique (DEEE) ou aux dispositions locales.

Entsorgungshinweis: Entsorgen Sie das Gerät gemäß der Richtlinie 2012/19/EU „für Elektro- und Elektronik-Altgeräte (WEEE)“ bzw. den Landesvorschriften.

Avviso per lo smaltimento: smaltire lo strumento in base alla Direttiva 2012/19/EU, sui "rifiuti di apparecchiature elettriche ed elettroniche (WEEE)" o le ordinanze locali

Aviso de eliminación: elimine el instrumento de conformidad con la Directiva 2012/19/UE sobre residuos de aparatos eléctricos y electrónicos (RAEE) o las ordenanzas locales

	<p>CE Marking – Indicates compliance with the requirements of the Directive 2014/30/EU on Electromagnetic Compatibility and the Directive 2014/35/EU on Low Voltage</p> <p>CE-markering – Geeft aan dat wordt voldaan aan de vereisten van Richtlijn 2014/30/EU inzake elektromagnetische compatibiliteit en Richtlijn 2014/35/EU inzake laagspanning</p> <p>Marquage CE – Indique la conformité aux exigences de la directive 2014/30/UE sur la compatibilité électromagnétique et de la directive 2014/35/UE sur la basse tension</p> <p>CE-Kennzeichnung – Zeigt die Einhaltung der Anforderungen der Richtlinie 2014/30/EU über elektromagnetische Verträglichkeit und der Richtlinie 2014/35/EU über Niederspannung</p> <p>Marcatura CE – Indica la conformità ai requisiti della Direttiva 2014/30/UE sulla Compatibilità Elettromagnetica e della Direttiva 2014/35/UE sulla Bassa Tensione</p> <p>Marcado CE: indica el cumplimiento de los requisitos de la Directiva 2014/30 / UE sobre compatibilidad electromagnética y la Directiva 2014/35 / UE sobre baja tensión.</p>
	<p>Date of manufacture</p> <p>Productiedatum</p> <p>Date de fabrication</p> <p>Herstellungsdatum</p> <p>Data di produzione</p> <p>Fecha de fabricación</p>

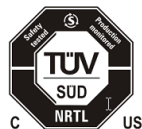

TÜV SÜD Certification Mark – Type tested; production monitored

TÜV SÜD certificeringsmerk - type getest; productie bewaakt

TÜV SÜD Marque de certification – Type testé ; production contrôlée

TÜV SÜD-Prüfzeichen – Typ geprüft; Produktion überwacht

Marchio di certificazione TÜV SÜD: tipo testato, produzione monitorata

Marca de certificación TÜV SÜD: tipo probado, producción controlada



This product complies with environmental protection use period as defined in People's Republic of China Electronic Industry Standard SJ/T11364-2006. Toxic or hazardous substances will not leak or mutate under normal operating conditions for 40 years.

Dit product voldoet aan de milieubeschermingsgebruiksperiode zoals gedefinieerd in de Electronic Industry Standard SJ/T11364-2006 van de Volksrepubliek China. Giftige of gevaarlijke stoffen zullen onder normale bedrijfsomstandigheden gedurende 40 jaar niet lekken of muteren.

Ce produit est conforme à la période d'utilisation dans le cadre de la protection de l'environnement telle que définie par la norme de l'industrie électronique de la République populaire de Chine SJ/T11364-2006. Les substances toxiques ou dangereuses ne fuiront pas ou ne subiront pas de mutation dans des conditions de fonctionnement normales pendant 40 ans.

Dieses Produkt entspricht der Umweltschutz-Nutzungsdauer gemäß der Definition im Electronic Industry Standard SJ/T11364-2006 der Volksrepublik China. Giftige oder gefährliche Stoffe werden unter normalen Betriebsbedingungen 40 Jahre lang nicht austreten oder mutieren.

Questo prodotto è conforme al periodo di utilizzo della protezione ambientale come definito nello Standard del settore elettronico della Repubblica Popolare Cinese SJ/T11364-2006. Le sostanze tossiche o pericolose non fuoriescono o non subiscono mutazioni in condizioni operative normali per 40 anni.

Este producto cumple con el periodo de uso de protección ambiental según el estándar SJ/T11364-2006 de la República Popular China para la industria electrónica. Las sustancias tóxicas o peligrosas no se filtrarán ni mutarán en condiciones de funcionamiento normales durante 40 años.



UK Conformity Assessed marking is a certification mark that indicates conformity with the applicable requirements for products sold within Great Britain.

De 'UK Conformity Assessed'-markering is een certificeringsmerk dat aangeeft dat producten die in Groot-Brittannië worden verkocht, voldoen aan de toepasselijke eisen.

Le marquage UK Conformity Assessed est une marque de certification qui indique la conformité aux exigences applicables aux produits vendus en Grande-Bretagne.

Die Kennzeichnung „UK Conformity Assessed“ ist ein Zertifizierungszeichen, das die Konformität mit den geltenden Anforderungen für in Großbritannien verkaufte Produkte anzeigt.

Il marchio UKCA (conformità valutata del Regno Unito) è un marchio di certificazione che indica la conformità ai requisiti applicabili per i prodotti venduti in Gran Bretagna.

El mercado UKCA (UK Conformity Assessed) es una marca de certificación que indica la conformidad con los requisitos aplicables para los productos vendidos en Gran Bretaña.




EAC-MED is a certification mark to indicate products that conform to all the safety and quality requirements of the Eurasian Customs Union. It means that the EAC-MED marked products meet all requirements of the corresponding technical regulations and have passed all conformity assessment procedures.

EAC-MED is een certificeringsmerk om producten aan te duiden die voldoen aan alle veiligheids- en kwaliteitseisen van de Euraziatische douane-unie. Dit betekent dat de producten met een EAC-MED-markering aan alle eisen van de desbetreffende technische voorschriften voldoen en alle conformiteitsbeoordelingsprocedures hebben doorlopen.

EAC-MED est une marque de certification qui indique la conformité des produits à toutes les exigences de sécurité et de qualité de l'Union douanière eurasiatique. Cela signifie que les produits marqués EAC-MED satisfont à toutes les exigences des réglementations techniques correspondantes et ont passé toutes les procédures d'évaluation de la conformité.

EAC-MED ist ein Zertifizierungszeichen zur Kennzeichnung von Produkten, die allen Sicherheits- und Qualitätsanforderungen der Eurasischen Zollunion entsprechen. Das bedeutet, dass die EAC-MED-gekennzeichneten Produkte alle Anforderungen der entsprechenden technischen Bestimmungen erfüllen und alle Konformitätsbewertungsverfahren bestanden haben.

EAC-MED è un marchio di certificazione che indica prodotti conformi a tutti i requisiti di sicurezza e qualità dell'Unione doganale eurasiatica. Ciò significa che i prodotti con marchio EAC-MED soddisfano tutti i requisiti dei regolamenti tecnici corrispondenti e hanno superato tutte le procedure di valutazione della conformità.

	<p>EAC-MED es una marca de certificación para indicar productos que cumplen con todos los requisitos de seguridad y calidad de la Unión Aduanera Euroasiática. Significa que los productos con la marca EAC MED cumplen todos los requisitos de los reglamentos técnicos correspondientes y han superado todos los procedimientos de evaluación de conformidad.</p>
	<p>Product complies with Australian Communications Requirements EESS - The Regulatory Compliance Mark (RCM) ACMA Labeling Requirements</p> <p>Product voldoet aan de Australische communicatie-eisen EESS - De markering voor naleving van de regelgeving (RCM) ACMA-etiketteringsvoorschriften</p> <p>Le produit est conforme aux exigences australiennes en matière de communication EESS - Marque réglementaire de conformité (RCM) Exigences en matière d'étiquetage ACMA</p> <p>Das Produkt entspricht den australischen Kommunikationsanforderungen. EESS – Kennzeichnung „Regulatory Compliance Mark“ (RCM) ACMA-Kennzeichnungsanforderungen</p> <p>Il prodotto è conforme ai requisiti Australian Communications Requirements EESS: marchio di conformità alle normative Requisiti di etichettatura ACMA</p> <p>El producto cumple con los requisitos de comunicaciones de Australia. EESS: marcado RCM (Regulatory Compliance Mark) de cumplimiento de la normativa. Requisitos de etiquetado de ACMA</p>



Korea Certification (KC) mark signifies Korea product compliance mark for safety and EMC/Radio/SAR of electrical and electronic equipment. The EMC requirements are applied to Agilent products.

Korea Certification (KC)-merkteken staat voor Korea-productconformiteitsmerk voor veiligheid en EMC/Radio/SAR van elektrische en elektronische apparatuur. De EMC-eisen worden toegepast op Agilent-producten.

La marque de certification coréenne (KC) signifie la marque de conformité des produits coréens pour la sécurité et l'EMC/Radio/SAR des équipements électriques et électroniques. Les exigences CEM s'appliquent aux produits Agilent.

Das Korea-Zertifizierungszeichen (KC) bezeichnet das koreanische Produktkonformitätszeichen für Sicherheit und EMV/Funk/SAR von elektrischen und elektronischen Geräten. Die EMV-Anforderungen gelten für Agilent-Produkte.

Il marchio Korea Certification (KC) indica il marchio di conformità del prodotto Corea per la sicurezza e EMC/Radio/SAR di apparecchiature elettriche ed elettroniche. I requisiti EMC vengono applicati ai prodotti Agilent.

La marca de certificación de Corea (KC) significa la marca de cumplimiento de productos de Corea para la seguridad y EMC / Radio / SAR de equipos eléctricos y electrónicos. Los requisitos de EMC se aplican a los productos Agilent.

Conformance to Standards

The Synergy Neo 2 meets the requirements of the following standards:

2014/35/EU – Low Voltage Directive

2014/30/EU – EMC Directive

2011/65/EU (with exemptions) and (EU) 2015/863 – RoHS Directives

2012/19/EU – WEEE Directive as amended by (EU) 2018/849

2006/42/EC of the European Parliament and of the Council of 17 May 2006 on machinery

Standard	Description
IEC QC 080000	IEC Quality Assessment System for Electronic Components (IECQ System) - Hazardous Substance Process Management (HSPM) System Requirements
UL 61010-1	UL Standard for Safety Electrical Equipment For Measurement, Control, and Laboratory Use; Part 1: General Requirements
EN 61010-1	Safety Requirements for Electrical Equipment For Measurement, Control, and Laboratory Use – Part 1: General Requirements
EN 61010-2-010	Safety Requirements for Electrical Equipment For Measurement, Control, and Laboratory Use – Part 2-010: Particular requirements for laboratory equipment for the heating of materials
EN 61010-2-081	“Particular requirements for automatic and semi-automatic laboratory equipment for analysis and other purposes.”
EN 61010-2-010	“Particular requirements for laboratory equipment for the heating of materials.”
EN 60825-1	"Safety of laser products. Part 1: Equipment classification and requirements."
CAN/CSA C22.2 No. 61010-1	Safety Requirements for Electrical Equipment For Measurement, Control, and Laboratory Use – Part 1: General Requirements
CAN/CSA C22.2 No. 61010-2-010	Safety Requirements for Electrical Equipment For Measurement, Control, and Laboratory Use – Part 2-010: Particular requirements for laboratory equipment for the heating of materials

EMC Information and Technical Description

The Synergy Neo 2 conforms to:

Emissions:

EN55011/CISPR 11, Class A

CFR Title 47 FCC Part 15 Subpart B, Class A

ICES-001, Issue 5, Class A (CAN ICES-001(A)/NMB-001(A))

ACMA AS/NZS CISPR 11, Class A

Immunity:

EN/IEC 61326-1 and 61326-2-6

ELECTRICAL EQUIPMENT FOR MEASUREMENT, CONTROL AND LABORATORY USE

PART 1: GENERAL REQUIREMENTS FOR (NON IVD) LISTED PRODUCTS

Ingress Protection Code

IP 20. Protected against solid foreign objects of 12.5 mm diameter and greater. No protection against water.

Disposal

Dispose of the instrument according to Directive 2012/19/EU, “on waste electrical and electronic equipment (WEEE)” or local ordinances.

Chapter 1

Introduction

This chapter introduces the Synergy Neo 2 Multi-Mode Microplate Reader, describes its hardware and software features, and provides contact information for technical assistance.

Product Description	2
Package Contents	3
Optional Accessories	4
Technical Support	7

Product Description

The Synergy Neo 2 is a multi-mode microplate reader. Depending on the model, Synergy Neo 2 detection modes include fluorescence intensity (FI), fluorescence polarization (FP), time-resolved fluorescence (TRF) — using a xenon flash or nitrogen laser, luminescence, UV-visible absorbance, and alpha. The reader is modular, and upgrade options are available; contact BioTek Customer Care for more information.

The reader is computer-controlled using Gen5 software for all operations, including data reduction and analysis. The Synergy Neo 2 is robot accessible and compatible with the BioTek Synergy Neo Stacker and BioStack 4; and with BioSpa 8 Automated Incubator. Gen5 supports OLE automation to facilitate the Synergy Neo 2's integration into an automated system.

The Synergy Neo 2 can perform reads using barcoded filter cubes or a monochromator.¹ The filter-based system can perform top and bottom fluorescence, luminescence, and alpha reads. Filter fluorescence uses a xenon flash light source or TRF laser in equipped models, along with interference filters and dichroic mirrors for wavelength specificity and up to three photomultiplier tubes (PMTs). To run a fluorescence polarization protocol, the filter cube must contain polarizing filters. Luminescence is measured through an empty filter position in the filter cube, (filters can be used if light filtering is necessary) and a dedicated broadband luminescence fiber.

The monochromator-based system, which has both top and bottom probes, is used for absorbance, fluorescence, and luminescence spectral reads. The xenon lamp allows for both UV and visible light measurements. The monochromator provides wavelength selection from 230–999 nm in 1 nm increments. Available absorbance and fluorescence read methods are endpoint, area scan, spectral scanning, and pathlength correction. You can use the monochromator optics for luminescence spectral scanning.

The alpha detection method can be used for endpoint reads using the top filter system in non-synchronized plate mode.

The Synergy Neo 2 has 4-Zone temperature control from 4°C over ambient to 70°C. (For alpha detection mode, the range is 3°C over ambient to 30°C.) Internal plate shaking, with both linear and orbital modes, is supported to ensure that reagents are properly mixed prior to reading.

The Synergy Neo 2 supports the reading of 6-, 12-, 24-, 48-, 96-, 384, and 1536-well microplates with 128 x 86 mm geometry as well as the Take3 and Take3 Trio Multi-Volume Plates.

Use of microplates other than those listed can result in positioning errors during program execution.

¹This dual light path capability is protected by U.S. patent number 8,218,141.

Models with injectors support dual-reagent dispensing to 6-, 12-, 24-, 48-, 96-, and 384-well microplates. An external dispense module pumps fluid from the supply bottles to the two injectors located inside the instrument.

See **Appendix A** for performance and technical specifications.

Package Contents

Package contents and part numbers are subject to change.

Item	Part #
<i>Synergy Neo 2 User Manual</i>	1351000N
Power cord set (specific to installation environment):	
Europe (Schuko)	75010
USA/International	75011
United Kingdom	75012
Australia/New Zealand	75013
USB cable	75108
Filter cube rack (except NEO2MON models)	1030541
Models with an external dispense module (packed separately), with the following accessories:	
Injector	8040541
Inlet tubes (2) from supply bottles to syringe drives	7082121
250- μ L syringes (2)	7083000
Syringe thumbscrews	19511
Priming plate	8042202
Injector tip priming trough	8042068
Dispense module communication cable	75107
Dispense module front cover	8042313
Dispense module box	8040579
Supply bottles (2, 30 mL)	7122609

Item	Part #
Supply bottle holders (2)	8042193
Injector tip cleaning stylus and plastic storage bag	2872304
Strap reagent racks (6)	7212035
Optional accessories per the sales order, unless shipped separately.	

Optional Accessories

Accessory availability and part numbers are subject to change.

Item	Part #
7-filter Absorbance Test Plate for absorbance measurement testing	7260522
Absorbance Test Plate for absorbance measurement testing at 340 nm	7260551
Synergy Neo 2 Product Qualification (IQ-OQ-PQ) package	1350526N
Microplate Barcode Scanner	1030008
Take3 Micro-Volume Plate	TAKE3
Take3 Trio Micro-Volume Plate	TAKE3TRIO
RS-232 serial cable	75034
PCR Tube Adapter Plates	6002072 6002076
BioCell Quartz Vessel	7272051
BioCell Adapter Plate	7270512
Replacement Shipping Materials	1033027 (instrument shipping box for all except T-models); 1843002 for T-models 1030009 Shipping hardware kit

Item	Part #
Luminometer Reference Microplate (includes microplate carrier adapter PN #8042263 for Synergy Neo2)	8030015
Gas controller, CO ₂ control only	1210012
Gas controller, CO ₂ and O ₂ control	1210013
Filter cubes	Contact BioTek
Fluorescence Test Plate	1400501

The Synergy Neo 2 is compatible with the Synergy Neo Stacker and BioStack 4. The stacker rapidly and systematically transfers microplates to and from the instrument's microplate carrier.

Materials for Conducting Liquid Tests

Manufacturer part numbers are subject to change.

Item	Part Number
Absorbance Liquid Tests	
BioTek Wetting Agent Solution	7773002
BioTek QC Check Solution #1 (25 mL)	7120779
BioTek QC Check Solution #1 (125 mL)	7120782
Phosphate-Buffered Saline (PBS) tablets (pH 7.2–7.6)	Sigma #P4417
β-NADH Powder (β-Nicotinamide Adenine Dinucleotide, reduced form)	98233 or Sigma #N6785-10VL
Fluorescence Liquid Tests	
Test Kits	
Kit with microplates and test solutions for conducting Corners/Sensitivity/Linearity (FI) tests using Sodium Fluorescein and Methylumbelliferone, and Time-Resolved Fluorescence (TRF) tests using Europium	7160010 (contains 7160013, 7160012, and 7160011 described below)

Item	Part Number
Kit for FI tests using Sodium Fluorescein	7160013
Kit for FI tests using Methylumbelliferone	7160012
Kit for TRF tests using Europium	7160011
Kit for Fluorescence Polarization (FP) test	7160014 or Invitrogen #P3088
Individual Materials	
Sodium Fluorescein Powder, 1 mg vial	98155
Methylumbelliferone, 10 mg vial	98156
Carbonate-Bicarbonate Buffer (CBB) capsules	Sigma #3041
Phosphate-Buffered Saline (PBS) tablets, pH 7.2–7.6	Sigma #P4417
Sodium Borate, pH 9.18	Fisher Scientific #159532 or equivalent
Injection System Tests	
BioTek Green Test Dye	7773003

Technical Support

See also "Contact Information" on page vi

Please be prepared to provide the following information:

- Your name and company information, along with a daytime phone or fax number, and/or an e-mail address
- The product name, model, and serial number
- The instrument software part number and basecode version (available via Gen5 for the Synergy Neo 2 by selecting **System > Instrument Control > Information**)
- The version of Gen5. From the main Gen5 screen, select **Help > About Gen5**.
- A copy of the Synergy Neo 2 system test results.
- For troubleshooting assistance or instruments needing repair, the specific steps that produce your problem and any error codes displayed in Gen5 (see also **Appendix B, Error Conditions**)

If you need to return an instrument to BioTek for service or repair, please contact Technical Support for instructions. Repackage the instrument according to the instructions at the end of **Chapter 2**

Installation

This chapter includes instructions for unpacking and setting up the Synergy Neo 2 and, if applicable, the external dispenser, the barcode reader, and the Synergy Neo Stacker. Instructions are also included for preparing the reader and dispenser for shipment.

Product Registration	9
Important Information	9
1: Unpack and Inspect the Reader	10
2: Unpack and Inspect the Dispenser	10
3: Unpack and Inspect the Gas Controller (if applicable)	11
4: Select an Appropriate Location	11
5: Remove the Shipping Hardware	12
6: Install the Power Supply	14
7: Install the Gas Controller (if applicable)	15
8: Install the Dispenser	15
9: Connect the Host Computer	19
10: Install Gen5 on the Host Computer	19
11: Turn on the Reader	20
12: Establish Communication	20
13: Set Dispenser Calibration Values	21
14: Run a System Test	22
15: Test the Injector System	24
Operational/Performance Qualification	25
Repackaging and Shipping Instructions	25

Product Registration

Register with BioTek to ensure that you receive important information and updates about the products you have purchased.

Important Information

WARNING

Two-person lift. The instrument should be lifted by two people. The instrument with all available modules weighs up to 45.5 kg.



CAUTION

Shipping Hardware. All shipping hardware must be removed before operating the instrument and reinstalled before repackaging the instrument for shipment.

- This chapter contains installation and setup tasks for the Synergy Neo 2 and accessories. Perform the tasks in the order presented.
- Save all packaging materials. Be sure to use packaging materials supplied by the manufacturer when shipping the reader. Using other forms of commercially available packaging, or failing to follow the repackaging instructions, may void your warranty.
- During the unpacking process, inspect the packaging, reader, and accessories for shipping damage. If the reader is damaged, notify the carrier and your BioTek representative. Keep the shipping boxes and the packaging materials for the carrier's inspection.

1: Unpack and Inspect the Reader

Save all packaging materials. If you need to ship the reader to BioTek for repair or replacement, you must use the BioTek-supplied materials. Using other forms of commercially available packaging, or failing to follow the repackaging instructions, may **void your warranty**.

During the unpacking process, inspect the packaging, reader, and accessories for shipping damage. If the reader is damaged, notify the carrier and your BioTek representative. Keep the shipping boxes and the packaging materials for the carrier's inspection. BioTek will arrange for repair or replacement immediately.

1. Open the shipping box, remove the instrument from the box, and place it on a level, stable surface.
2. Place the packaging materials back into the shipping box for reuse if the instrument needs to be shipped again.

2: Unpack and Inspect the Dispenser

Save all packaging materials. If you need to ship the reader to BioTek for repair or replacement, you must use the BioTek-supplied materials. Using other forms of commercially available packaging, or failing to follow the repackaging instructions, may **void your warranty**.

During the unpacking process, inspect the packaging, reader, and accessories for shipping damage. If the reader is damaged, notify the carrier and your BioTek representative. Keep the shipping boxes and the packaging materials for the carrier's inspection. BioTek will arrange for repair or replacement immediately.

If applicable, perform these steps to unpack the dispenser.

1. Open the shipping box. Remove the accessories box and foam insert that contains the injector tubing and bottle holders.
2. Lift out the dispenser and place it on a level surface.
3. Open the accessories box and remove its contents. The accessories should include the dispenser-related items listed under **Package Contents** in Chapter 1.
4. Place all packaging materials into the shipping box for reuse if the dispenser needs to be shipped.

3: Unpack and Inspect the Gas Controller (if applicable)

If applicable, perform these steps to unpack the gas controller.

1. Open the shipping box.
2. Lift out the accessories (power supply, tubing, and manual), and set them aside.
3. Lift out the gas controller, and place it on a level surface.
4. Place all packaging materials into the shipping box for reuse if the gas controller needs to be shipped.

4: Select an Appropriate Location

Install the reader on a level, stable surface in an area where temperatures between 18°C (64°F) and 40°C (104°F) can be maintained.

For alpha laser reads, the temperature should be no warmer than 30°C (86°F).

Leave at least six inches of space between the instrument's rear panel and any other object. This space ensures proper air flow in and out of the instrument.

If you intend to use the bottom filter access door on the left side of the instrument, ensure that you select a location that allows for easy access to the door.

The reader is sensitive to extreme environmental conditions. Avoid the following:


- **Excessive humidity.** Condensation directly on the sensitive electronic circuits can cause the instrument to fail internal self-checks. The humidity must be in the range of 10–85%, non-condensing.
- **Excessive ambient light.** Bright light may affect the reader's optics and readings, reducing its linear range.
- **Dust.** Readings may be affected by extraneous particles (such as dust) in the microplate wells. A clean work area is necessary to ensure accurate readings.


BioStack/Barcode Scanner: If you are installing the barcode scanner "Microplate Barcode Scanner" on page 175 or the Synergy Neo Stacker or BioStack 4 for operation with the Synergy Neo 2, you may wish to set them up now. This is a good time, for example, to seat the instruments in their aligning plates. (Refer to the *Synergy Neo Stacker User Manual* or *BioStack User Manual* for more information.)

5: Remove the Shipping Hardware

CAUTION **Shipping Hardware.** All shipping hardware must be removed before operating the instrument and reinstalled before repackaging the instrument for shipment.

1. Locate the shipping hardware: plate carrier shipping brackets, and, if equipped, top filter chamber shipping brackets, the bottom filter chamber shipping bracket, and T-model set screw with yellow tag.

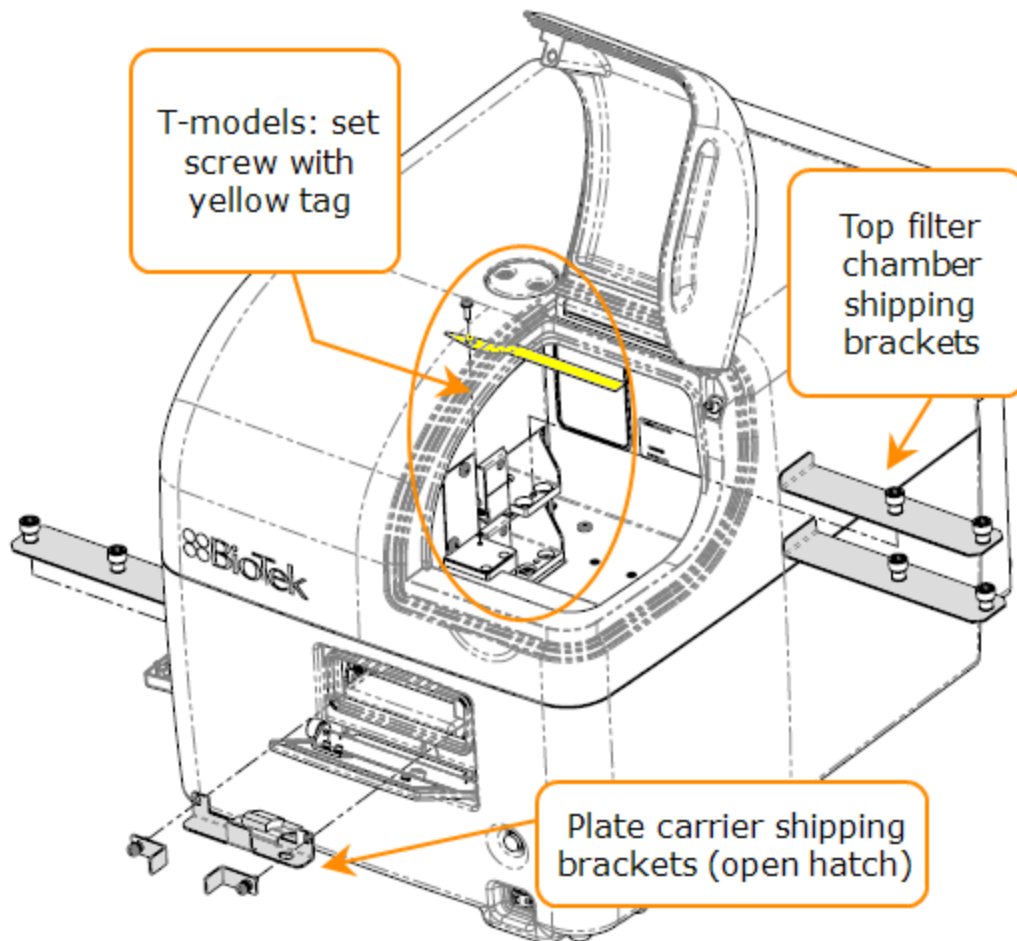

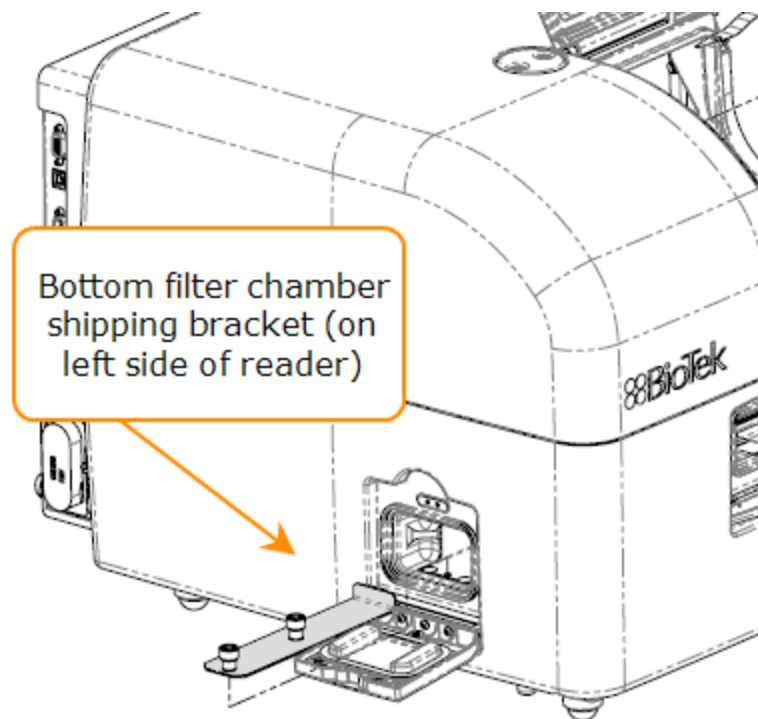

2. Using a #2 Phillips screwdriver, remove the screws holding the carrier shipping bracket. Push the carrier off the rear of the bracket, and then lift and withdraw the bracket from the instrument.

If equipped:

3. Remove the top filter chamber shipping bracket(s).
4. T-Models: remove the set screw and its yellow tag from the switching block.
5. Remove the bottom filter chamber shipping bracket.



6. Store the shipping brackets and hardware with the original packaging for reuse in case you need to ship the instrument.
7. T-Models: Remove safety keys taped to the top of the TRF laser, insert a key and turn to ON.

The filter cubes are shipped in a protective case. The **Getting Started** chapter explains when and how to install filter cubes.

6: Install the Power Supply

WARNING **Power Rating.** The instrument must be connected to a power receptacle that provides voltage and current within the specified rating for the system. Use of an incompatible power receptacle may produce electrical shock and fire hazards

WARNING **Electrical Grounding.** Never use a plug adapter to connect primary power to the instrument. Use of an adapter disconnects the utility ground, creating a severe shock hazard. Always connect the system power cord directly to an appropriate receptacle with a functional ground.

CAUTION **Power Supply.** Use only the power supply shipped with the instrument, and operate it within the range of line voltages listed on it.

Do **not** plug the power supply into a power receptacle until after the power supply is connected to the instrument.

1. Locate the power inlet on the left side of the reader.
2. Examine the power supply's plug. It has three small holes and a small groove that line up with the prongs and a tab inside the power inlet.
3. Insert the plug into the power inlet and finger-tighten the collar of the plug. Plug the power supply's cord into an appropriate power receptacle.

7: Install the Gas Controller (if applicable)

The gas controller is an external module that enables the user to control CO₂ and O₂ concentrations inside the attached instrument's reading chamber. If you purchased the module for operation with the Synergy Neo 2, refer to the *Gas Controller User Guide* for installation instructions.

8: Install the Dispenser

Place the dispense module on top of the reader or on top of the gas controller (if equipped). Do not place the dispenser next to the reader.



Synergy Neo2 with dispenser, and with dispenser and gas controller

1. Open the plastic bag containing the injector tube and tips. Remove the clear plastic shrouds from the tubes.

2. Remove the two inlet tubes from their plastic canisters.
3. Identify the two syringe valves on the dispense module. Each is labeled with a left-pointing arrow.


When installing the inlet and outlet tubes, do not use any tools. Finger-tighten only!

4. Screw the fitting of one inlet tube into the right side of the Syringe 1 valve.
5. Identify the #1 outlet tube, and screw it into the left side of the Syringe 1 valve.
6. Repeat these steps to attach the inlet and outlet tubing for Syringe 2.


It is critical that the tubing is installed in the correct ports. Otherwise, injected fluid may miss the intended well.

7. Remove the round tubing feed-through cover from the top of the reader (2 screws). Store the cover and screws with the shipping hardware in case the reader needs to be shipped again.
8. Thread the injector tip holder, with outlet tubing connected to both ports, through the hole in the top of the reader.
9. Open the reader's top door, and, holding the injector tip holder by the tab, insert the injector tips into the appropriate holes inside the reader, behind the filter cube chamber (s).

Removing any filter cubes that are installed in the filter cube chambers may be necessary to reach the injector tip holes.

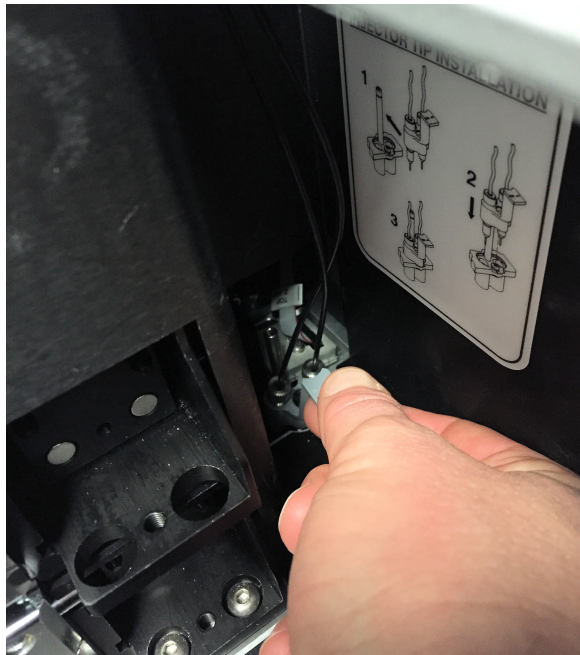

Injector tip holder in its socket

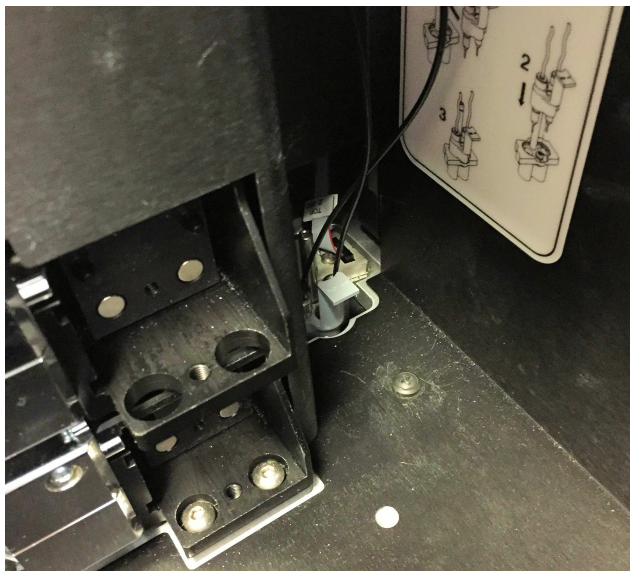


A magnet located between the injector tips helps to guide the tips into place and secures them in the reader.

10. Place the tubing feed-through cover over the hole in the top of the reader, and finger-tighten the thumbscrews to secure it.
11. Remove the two syringes from their protective boxes. They are identical and interchangeable.

12. Install both syringes.

- Hold the syringe vertically with the threaded end at the top.
- Screw the top of the syringe into the bottom of the syringe valve. Finger-tighten only.
- Carefully pull down the bottom of the syringe until it rests inside the hole in the bracket.
- Pass a thumbscrew up through this hole and thread it into the bottom of the syringe. Hold the syringe to prevent it from rotating while tightening the thumbscrew. Finger-tighten only.

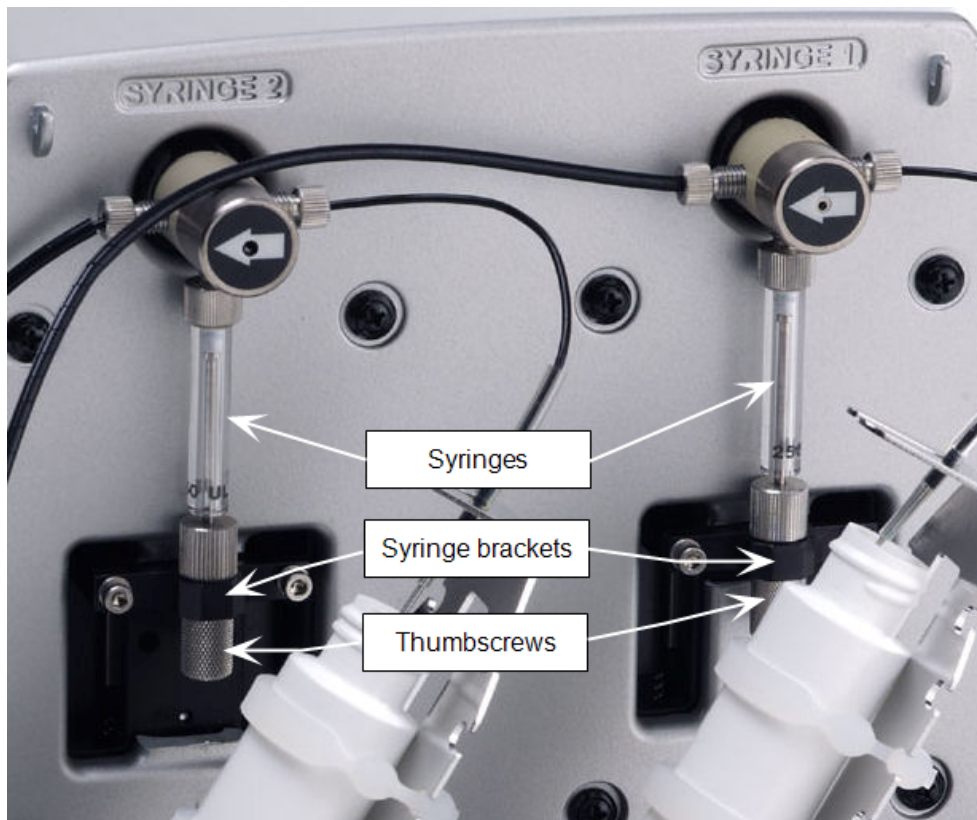

13. Locate the dispenser cable. Plug one end into the port on the back of the dispenser. Plug the other end into the “Dispenser Port” on the left side of the reader.
14. Locate the injector tip-cleaning stylus, packaged in a small cylinder. Attach the cylinder to the back of the dispenser for storage.

Perform a visual inspection or a Performance Qualification test after reconnecting the tubes.

9: Connect the Host Computer

The Synergy Neo 2 is equipped with two communication ports: USB and Serial (RS-232). Both ports are located on the left side of the reader.

- A USB cable is included in the accessories box; the RS-232 cable is available as an optional accessory.
- Connect one end to the appropriate port on the reader and the other end to the appropriate port on the host computer.

10: Install Gen5 on the Host Computer

The Synergy Neo 2 is controlled by Gen5 software running on a host computer. There is a certain sequence of events that **must** be followed to ensure that the software is properly installed and configured. Please follow the instructions provided with the Gen5 to install the software and the USB Driver.

Before Installing Gen5

- Gen5 v3.11 and higher require Windows 10 (Professional edition, 64-bit).
- If you will use Gen5 Quick Export or PowerExport, Microsoft Excel 2010 through 2016 (32- or 64-bit edition) is required.
- You must have administrator privileges to install Gen5. Log in to Windows as “Administrator” or consult your IT department for assistance

Procedure

1. Insert the Gen5 software USB flash drive into an available port on the computer. Using File Explorer, navigate to the flash drive and locate the **setup** file.
2. Double-click **setup** and follow the instructions to install the software. Choose **Typical** when prompted for the install type. Click **Finish** when prompted.
3. Enter the Gen5 Serial Number (located on the Gen5 storage case) when prompted, and click **Continue**.

11: Turn on the Reader

1. If Gen5 is open, close it now.

T-Models: Make sure the laser safety key is in its **ON** position.

2. The reader's power switch is located on the lower-right corner of the front panel. Turn the reader on. The reader performs a system test. When the test is completed, the reader extends the microplate carrier.

The carrier eject button, located above the reader's power switch, can be used to extend/retract the microplate carrier.

When the reader is idle, the light on the carrier eject button is green. When the reader is busy (when running an assay or a self-test) the light is red. If an error occurs, the light is red and blinking.

12: Establish Communication

If using the USB cable, refer to the instructions that shipped with the USB Driver Software on the Gen5 software media to install the necessary drivers.

1. Start Gen5 and log in if prompted. The default System Administrator password is **admin**.
2. Choose **Standard** mode in Gen5 v.3.06 and higher to complete instrument installation and setup. **Simple** mode can be used to operate the reader.
3. From the Gen5 main screen, select **System > Instrument Configuration** and click **Add Reader**.
4. Set the Reader Type to **Synergy Neo 2**, and click **OK**.
5. Perform one of the following steps:
 - Select **Plug & Play** and click on the Reader Type you are using.

A Synergy Neo 2 must be connected via USB to the computer and turned on to appear in the Available Plug & Play Readers list.

- Set the Com Port to the computer's COM port to which the reader is connected.

Com Port: When using the RS-232 cable you must define the COM port. With either method, com port information can be found via the Windows Control Panel, under Ports in the Hardware/Device Manager area of System Properties (e.g., USB Serial Port (COM5)).

6. Click **Test Communication**. Gen5 attempts to communicate with the reader. If the communication attempt is successful, return to Gen5's main screen.

Communication Errors

If the communication attempt is not successful, try the following:

- Is the reader connected to the power supply and turned on?
- Is the communication cable firmly attached to both the reader and the computer?
- Did you select the correct Reader Type in Gen5?
- Did you install the USB driver software, if applicable?
- Was the reader allowed to fully complete its power-up self-test before initiating communication with Gen5?

If you remain unable to get Gen5 and the reader to communicate with each other, contact Technical Support .

13: Set Dispenser Calibration Values

Applies only to models equipped with injectors.

Before using the external dispense module with the Synergy Neo2, you must set its calibration values in Gen5.

The calibration values for both dispensers (#1 and #2) are printed on labels affixed to the rear of the dispense module. Each label lists six target calibration values (e.g., 200, 80, 40) with their actual measured values (e.g., 199.3, 79.7, 39.9). You will enter the **measured** calibration values into Gen5.


1. If you have not already done so, power on the instrument, and establish communication.
2. In Gen5, go to **System > Instrument Configuration**, select your instrument, and click **View/Modify**.

3. Click **Setup**, and then select the **Dispenser 1 tab**.
4. On the keyboard, press CTRL+SHIFT+M to enter maintenance mode for the Dispenser 1 window.
5. Enter the syringe calibration values from the corresponding label on the rear of the dispenser box.
6. Click **Send Volumes**, and then click **Get Volumes** to verify that the entered values were sent to the instrument.
7. Select the **Dispenser 2 tab**, and repeat steps 4 through 6 for Dispenser 2.

14: Run a System Test

Running a system test will confirm that the reader is set up and running properly, or will provide an error code if a problem is detected.

"T" models: Make sure the laser safety key is installed and turned to ON; otherwise, the self-test will fail.

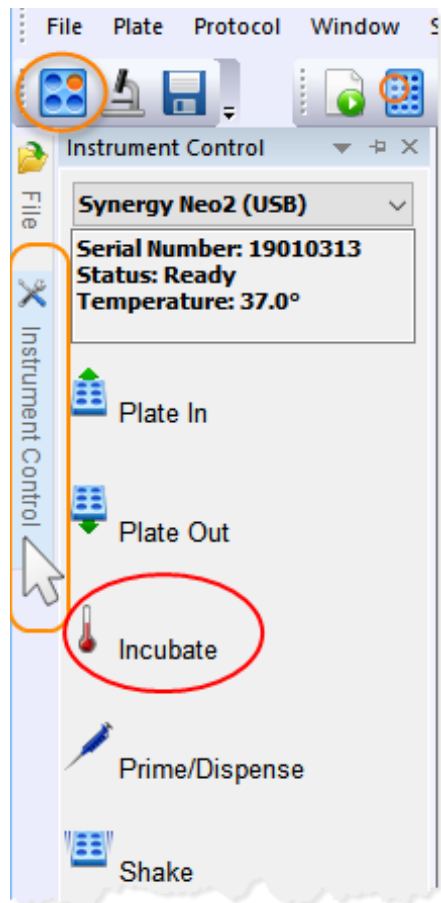

1. Turn on the incubator:
 - Open the **Instrument Control** panel by clicking the tab on the left border. At the top of the Instrument Control panel, make sure the correct instrument is selected (if you have more than one instrument).
 - Click **Incubate** and enter a Requested temperature of at least 37°C and click **On**.

Wait until the incubator temperature reaches the set point before continuing.

2. From Gen5's main view, select **System > Diagnostics > Run System Test**. If prompted to select a reader, select **Synergy Neo 2** and click **OK**.
3. When the test is completed, a dialog requesting additional information appears. Enter the information and click **OK**.

If a message appears that a pending system test is waiting from the initial power-up self-test, view the pending system test and repeat steps 2 and 3.


4. The results report appears. Near the top of the report in the Test Results area, the text should read "SYSTEM TEST PASS."
 - You may wish to print the report and store it with your records. The report can also be saved as a .txt file.
 - Gen5 stores system test information in its database; you can retrieve it at any time.
 - Sign and date the report, and store it with your test documentation.

If an error code is returned, refer to **Appendix B** and look up the code. If the problem is something you can fix, do so now and run another system test. If the problem is something you cannot fix, or if the test continues to fail, contact Technical Support.

5. Turn off the incubator.

15: Test the Injector System

1. If necessary, press the carrier eject button to extend the microplate carrier.
2. Place the tip priming trough in the rear pocket of the carrier.
3. Place the priming plate on the carrier.

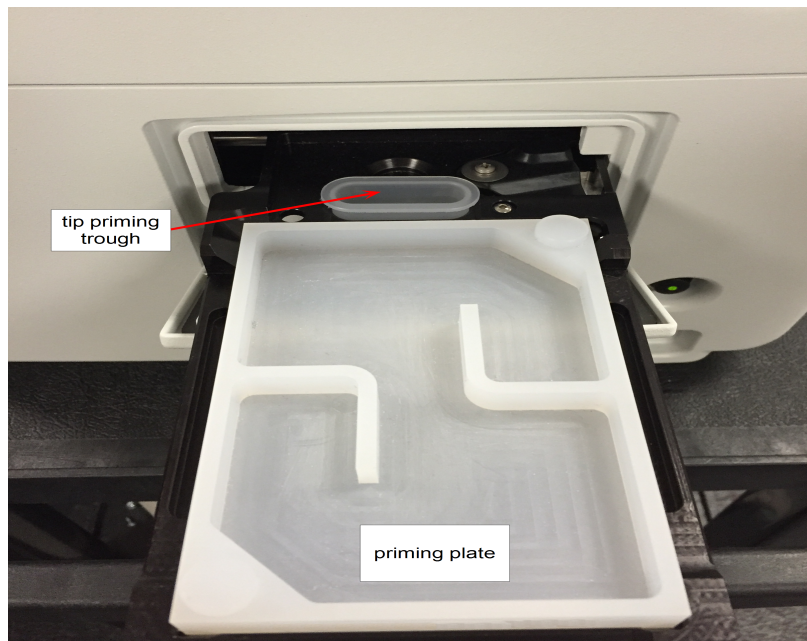

4. Fill the two reagent bottles with distilled or deionized water. Place the bottles in their holders, and place the holders directly in front of the syringes. Insert the inlet tubes into the bottles.
5. Select **Prime/Dispense** from the Instrument Control panel or from the Gen5 main screen, select **System > Instrument Control > Synergy Neo 2**.
6. Click the **Prime** tab.
7. With Dispenser set to **1**, set the Volume to **5000 µL** and click **Prime**.

The syringe should move down and up repeatedly, drawing fluid from the bottle. The fluid should pump through the tubing and dispense into the priming plate. Examine the fittings; no leaks should be detected. If leaks are detected, tighten all fittings and repeat the prime. If leaks are still detected, contact **Technical Support**.

8. When the prime finishes, set Volume to **2000 µL** and click **Purge** to clear the fluid lines.
9. Set Dispenser to **2** and repeat steps 7 and 8.
10. When finished, remove and empty the priming plate.
11. Return to the Gen5 main screen.

The installation and setup process is complete.

Operational/Performance Qualification

Your Synergy Neo 2 was fully tested prior to shipment and should operate properly following the successful completion of the installation and setup procedures described in this chapter.

If you suspect that problems occurred during shipment, if you received the reader back following service or repair, or if regulatory requirements dictate that Operational/Performance Qualification is necessary, refer to the **Instrument Qualification Process** chapter.

A Product Qualification & Maintenance (IQ/OQ/PQ) package for the Synergy Neo 2 is available for purchase (PN 1350526N).


Repackaging and Shipping Instructions


Important! Please read all of the information provided below before preparing the Synergy Neo 2 for shipment.

- Contact Technical Support before returning equipment for service.
- Decontamination prior to shipment is required by the U.S. Department of Transportation regulations.
- If the Synergy Neo 2 has been exposed to potentially hazardous material, decontaminate it to minimize the risk to all who come in contact with the instrument during shipping, handling, and servicing. The **Maintenance** chapter contains decontamination instructions.


- Remove the microplate and tip prime trough (if equipped) from the carrier before shipment. Spilled fluids can contaminate the optics and damage the instrument.
- Install the shipping hardware (see next section).
- The instrument's packaging design is subject to change. If the instructions in this document do not apply to the packaging materials you are using, contact Technical Support for guidance.
- Be sure to use packaging materials supplied by the manufacturer. Other forms of commercially available packaging are not recommended and can void the warranty. You can order replacement shipping materials (PN 1033027 , except "T" models order PN 1843002) if necessary. The shipping box, accessories box, foam caps, and so on are included as a whole set under these part numbers and cannot be ordered separately.
- If you have misplaced the shipping hardware, order part number 1030009 which includes the following items:
 - Carrier shipping bracket [x1] (PN 1030539)
 - Carrier shipping bracket screws [x2] (PN 1032190)
 - Filter chamber shipping bracket [x3] (PN 1030527)
 - T-Models: Set screw for switching block (PN 01735) and its yellow tag (7091006)

Replace the Shipping Hardware

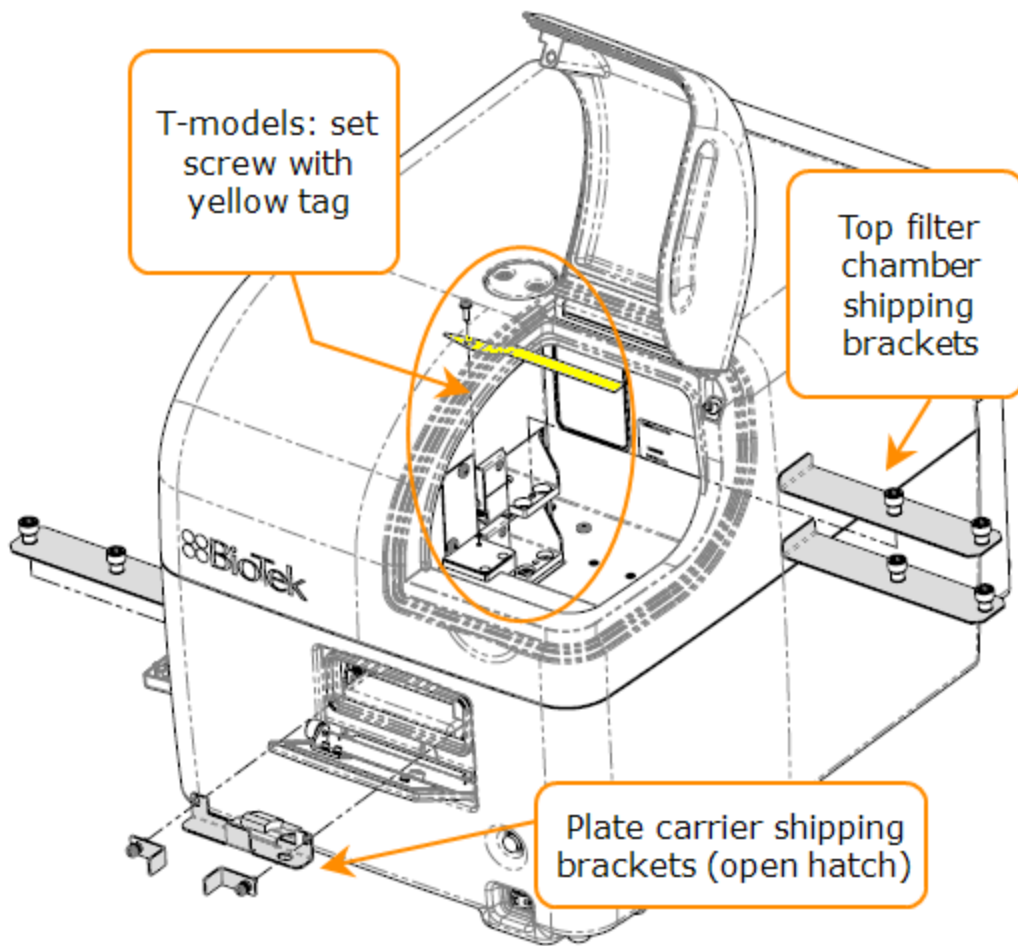

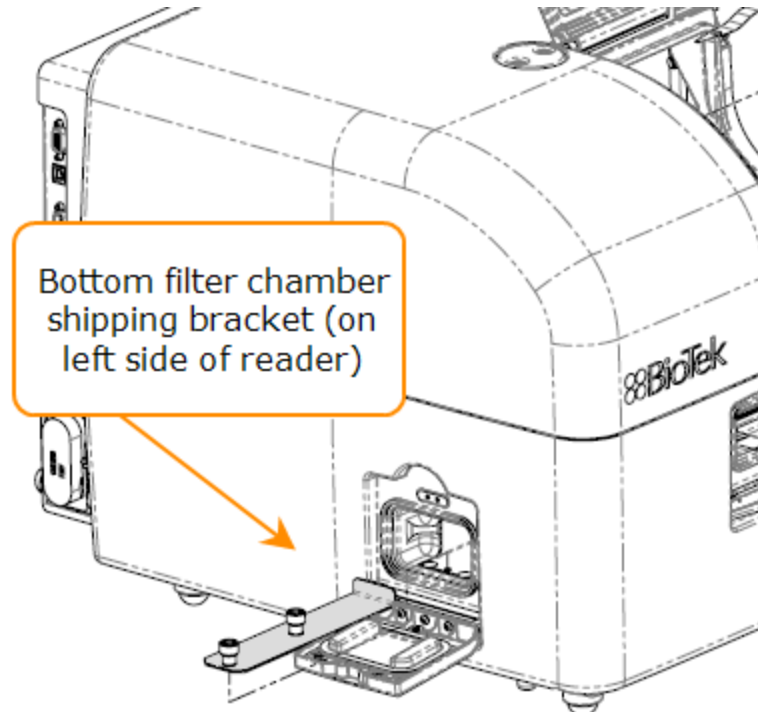

1. Open the carrier door and push the carrier into the reader about 2 inches.
2. Place the big plate carrier shipping bracket (without the hardware) in the read chamber so that the two holes on the horizontal surface of the bracket slip over the two studs sticking up from the read chamber floor.
3. Pull the carrier out toward the front of the instrument until it rests against the foam piece on the sloped part of the shipping bracket.
4. Holding the carrier against the shipping bracket, place one of the screw brackets against one of the holes at either end of the carrier bracket, so that the screw is pointing toward the hole in the carrier shipping bracket.
5. Insert the screw so that it passes through both shipping brackets and catches the threads in the carrier itself.
6. Repeats steps 4 and 5 to install the second screw bracket.
7. Remove all filter cubes from the instrument. Place them into plastic bags, and then place them into the filter cube case.

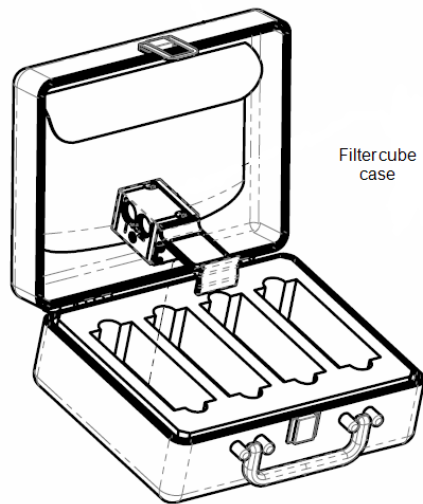


8. Install the top filter chamber shipping brackets using the inner mounting screws.
9. T-Models: Install the set screw and its yellow tag into the switching block.
10. If equipped, install the bottom filter chamber shipping bracket using the outer mounting screw.

Repack the Instrument


1. Contact Technical Support for shipping instructions before returning equipment for service.
2. Decontaminate the reader and, if attached, the dispense module, according to the instructions provided in the **As-Needed Maintenance** chapter.
3. If you will also be shipping the dispense module, see **Preparing the Dispenser for Shipment** on page 31.

If you are not shipping the dispenser, purge all fluids from the lines, and then disconnect it from the reader.

4. Place the accessories in the accessories box, then seal the accessories box with tape.

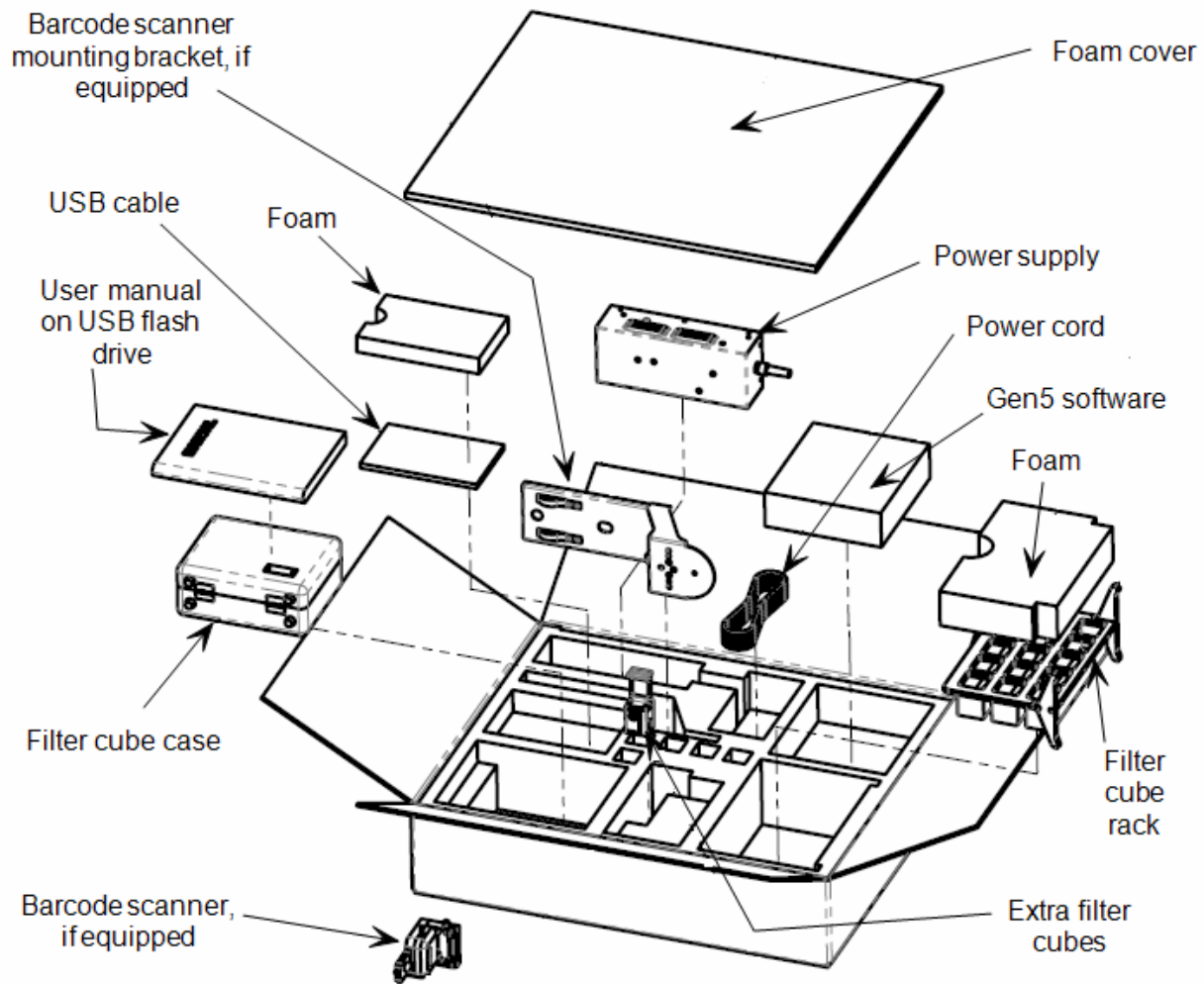

5. Place the instrument in the large plastic bag, then place it in the interior box, surrounded by foam boards. Then, place the interior shipping box, surrounded by the foam corners, into the external shipping box. Refer to the next figure.

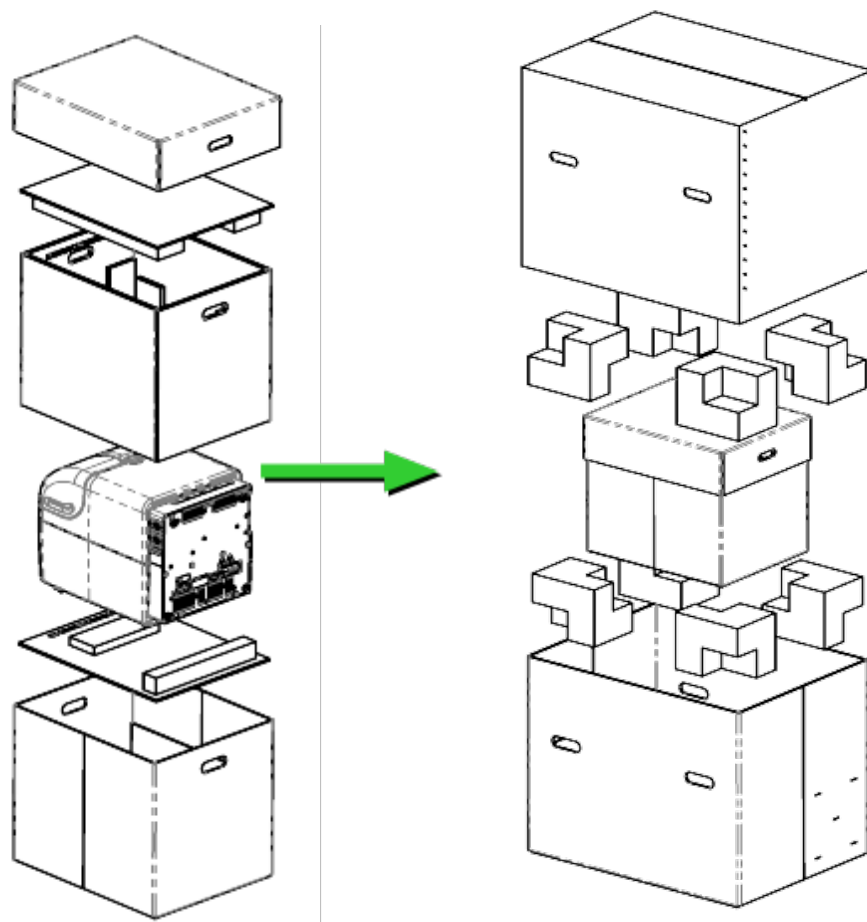

Preparing the Dispenser for Shipment

1. If you have not already done so, contact Technical Support for shipping instructions.
2. Decontaminate the dispenser according to the instructions in the **As-Needed Maintenance** chapter. Be sure to purge the dispenser of all fluid when finished.
3. With the reader on, start Gen5 and select **System > Instrument Control > Synergy Neo 2**.
4. Perform this step twice, for both dispensers: Click the **Prime** tab and set the dispenser number (1 or 2). Click **Maintenance**. The syringe bracket lowers itself. Remove the thumbscrew from underneath the bracket. Carefully unscrew the top of the syringe from the syringe valve. Lift out the syringe and store it in its original box.
5. Fully detach the dispenser from the reader. (The screws are stored in the plastic bag attached to the back of the dispenser.) Set the dispenser aside for the moment.
6. Remove the tip priming trough and store it in the dispenser accessories bag.

7. Remove the two inlet tubes from the syringe valves and store them in their plastic canisters.
8. Remove the two outlet tubes from the syringe valves. Attach the clear plastic shrouds to the fittings of the outlet tubes. Place the tubes in a plastic bag.
9. Remove the front cover from the dispenser.
10. Insert the bottom foam end cap in the dispenser accessories shipping box and place the accessories in the insert.
11. Insert the bottom foam end cap in the shipping box, and place the dispenser inside the end cap.
12. Insert the foam insert that holds the reagent bottle holders and injector tubing into the shipping box and place the bottle holders and tubing in it.
13. Slide the dispenser accessories box into the shipping box.
14. Insert the top foam end cap. Close and seal the outer box with tape.

Getting Started

This chapter describes some of the Synergy Neo 2's external and internal components, and provides an introduction to using Gen5 software to control the instrument.

Modular Design	35
External Components	37
Internal Components	39
Gen5 Software	55
Recommendations for Optimum Performance	57

Modular Design

The Synergy Neo 2 is a multi-mode microplate reader, with a design that allows you to initially purchase only the detection capabilities you need and then upgrade later as your requirements expand. Please contact BioTek to learn more about your upgrade options.

Gen5 software is used to control the reader. If the reader is connected and turned on, Gen5 will present only those options that apply to your reader model. For example, if your model is not equipped with the alpha detection system, alpha will not be available as a detection method.

The module letters form the part number for each Synergy Neo 2 model, which is shown on a label on the reader.

These configurations include dual top PMTs for fluorescence and luminescence (Lum).

Part Number	Monochromator-based			Filter-based			Alpha laser	TRF laser
	UV-Vis absorbance	top/bottom FI	Lum	top FI/FP/TRF	bottom FI/FP/TRF	Lum		
N2-SN	•			•		•		
N2M-SN	•	•	•	•		•		
N2MB-SN	•	•	•	•	•	•		
N2MA-SN	•	•	•	•		•	•	
N2MAB-SN	•	•	•	•	•	•	•	
N2MT-SN	•	•	•	•		•		•
N2MAT-SN	•	•	•	•		•	•	•
N2MABT-SN	•	•	•	•	•	•	•	•
N2MBT-SN	•	•	•	•	•	•		•
N2T-SN	•			•		•		•
N2AT-SN	•			•		•	•	•
N2ABT-SN	•			•	•	•	•	•


These configurations include a single top PMT for fluorescence and luminescence.

Part Number	Monochromator-based			Filter-based			Alpha laser
	UV-Vis absorbance	top/bottom FI	Lum	top FI/FP/TRF	bottom FI/FP/TRF	Lum	
N2S-SN	•			•		•	
N2MON-SN	•	•	•				
N2SM-SN	•	•	•	•		•	
N2SMB-SN	•	•	•	•	•	•	
N2SA-SN	•			•		•	•
N2SMA-SN	•	•	•	•		•	•
N2SMAB-SN	•	•	•	•	•	•	•

External Components



1	Upper and lower top filter cubes access door
2	Light-blocking microplate carrier access door
3	Microplate carrier eject button
4	Power switch
5	Bottom filter cube access door (if equipped)
6	Entry port for the dispense outlet tubes and injectors (if equipped)
7	not shown -"T" models, TRF laser attached to back of instrument (page 53)

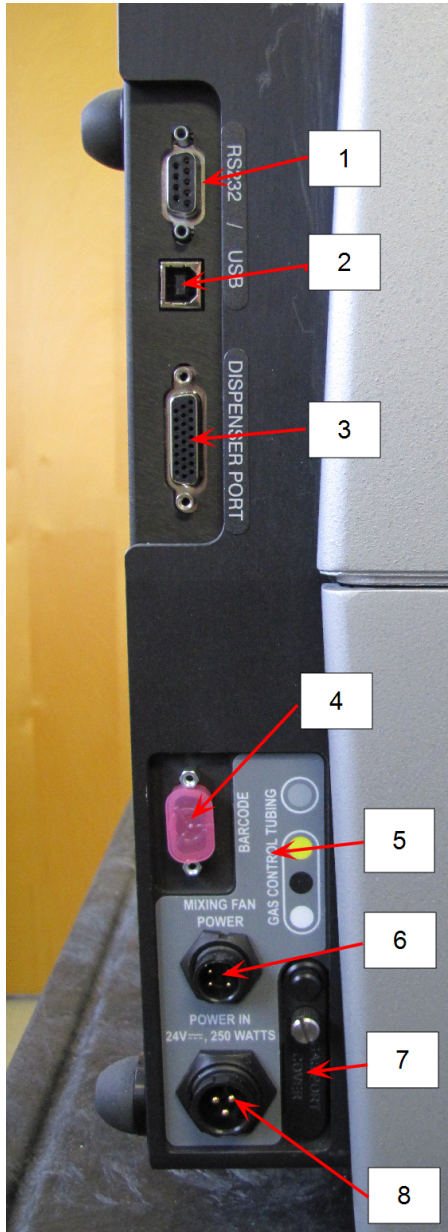

- 1: RS232 cable port
- 2: USB port
- 3: Dispenser port
- 4: Barcode reader cable port (with cover installed)
- 5: Legend for gas controller ports
- 6: Mixing fan power port
- 7: Gas controller ports (with cover installed)
- 8: Instrument power inlet

Note: TRF Laser is not shown.

Internal Components

Component	Description	Page
Barcoded Filter Cubes	The barcoded filter cubes contain excitation and emission filters, mirrors, and polarizing filters. Preconfigured cubes are available for purchase.	page 39
Shaking System	The instrument supports six different shake rates from each of three different shake types: linear, slow orbital, and fast orbital.	page 47
Injector System	Optional dispense module accessory: Syringes and tubing may require replacement over time. The tubing and injectors require cleaning at regular intervals.	page 51
Alpha Laser	A-models only: Equipped with an alpha laser for AlphaScreen, Alpha LISA, and other assays.	page 53
TRF Laser	T-models with TRF Laser only: Equipped with a 337 nm nitrogen laser for TR-FRET/HTRF/TRF applications. Gen5 version 3.06 or higher is required.	page 53

Filter Cubes


The Synergy Neo 2 is equipped with up to three filter cubes to perform filter-based fluorescence (using xenon flash or TRF laser), luminescence, and alpha reads. Each cube consists of excitation and emission filters, dichroic mirrors, polarizing filters, and/or unique apertures. When placed into the instrument, the unique component information for each cube is automatically identified via an internal barcode scanner.

Important: If two or more custom filter cubes (with ID 255) are configured in Gen5, it is the customer's responsibility to install the correct filter cube when performing a read.

The top filter cubes are accessed through the hinged door on the front of the reader. The bottom filter cube (if equipped) is accessed through a door on the left side of the reader. See the photo under **External Components** on page 37.

Each filter cube is labeled with an ID number (e.g., "107" or "61") and a short description of its contents and purpose (e.g., "360/460, 485/528" or "FP 485/530, LUM"). Each cube also has a barcode label that represents this information.

Be careful not to scratch the surface of the barcode label (see below). If the label is damaged, the Synergy Neo 2 may have trouble reading it correctly, requiring you to manually enter the filter cube information in Gen5.

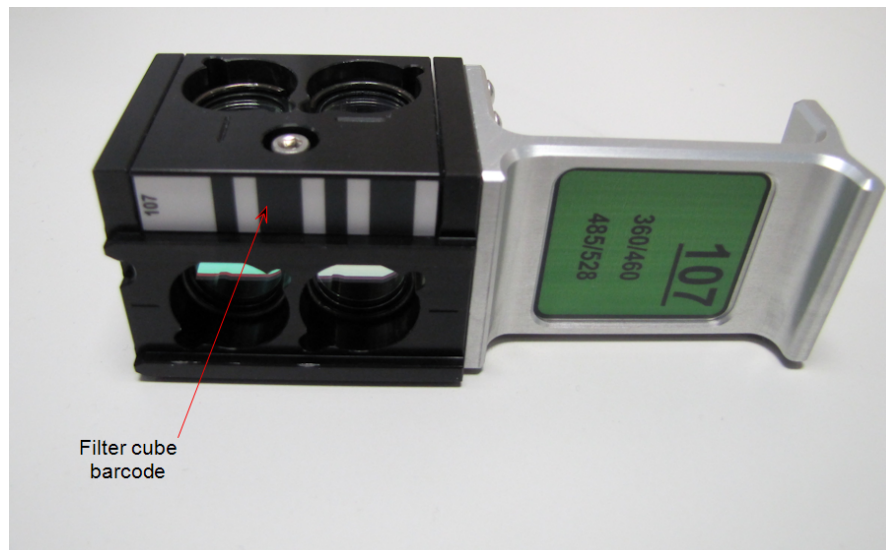


Filter Cubes and Gen5


When the instrument reads a new filter cube barcode label, information such as cube position type, filter and mirror specifications, and PMT location is passed to Gen5.

You can view the filter cube's ID number by selecting **System > Information Configuration > Synergy Neo > View/Modify > Setup**. If a location shows "?", either no filter cube is installed or the system could not read the barcode label. You can manually enter the filter cube ID number using the Override feature.

You can use the Gen5 Filter Cube Library (**System > Optics Library > Filter Cubes**) to define all cubes that you have available at your site and to import filter cube definitions, if necessary. Read step options in a protocol or experiment will be available for selection based on the filter cubes that you identify as being "On Site."

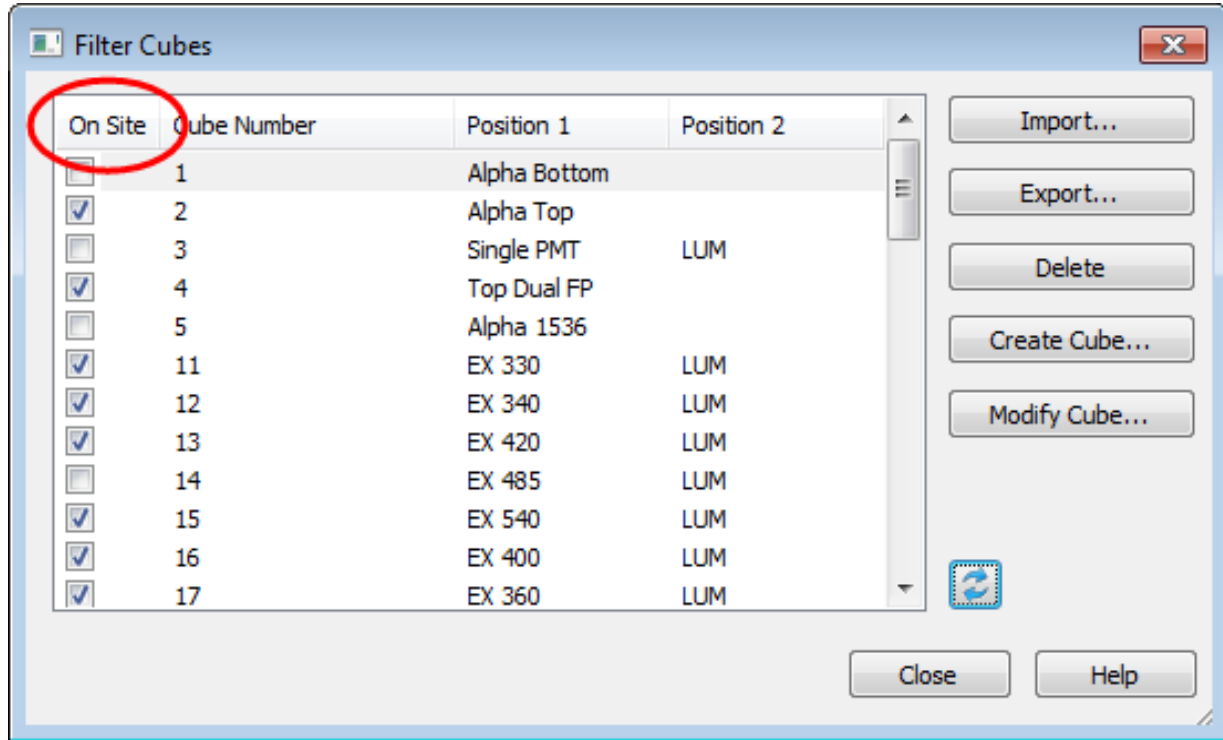

Gen5 also uses the filter cube information to prompt for filter cube installation at runtime, if necessary.

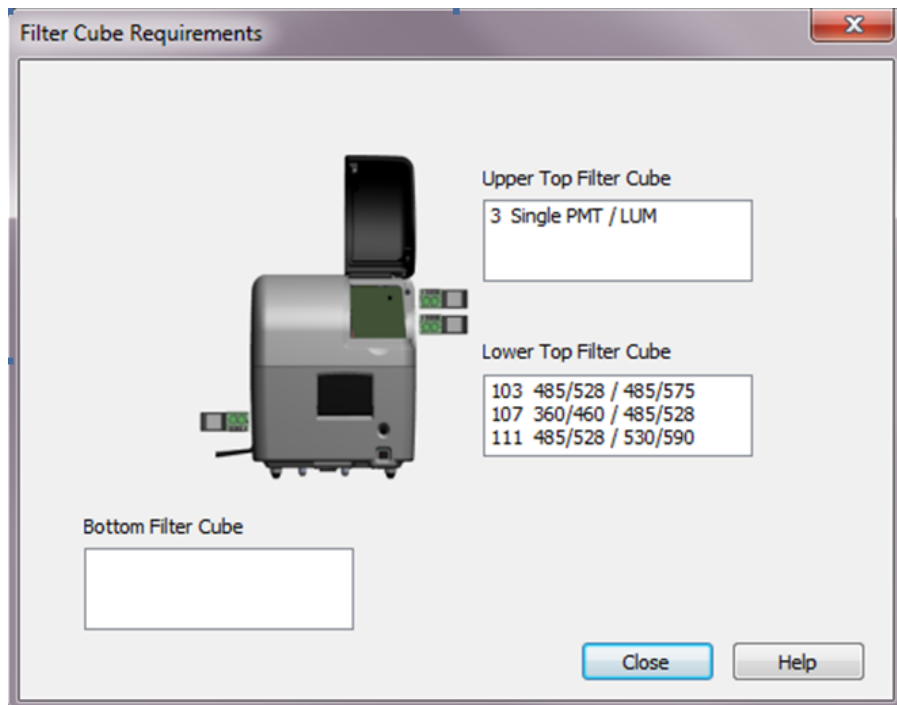

Installing or Removing a Filter Cube

Filter cubes can be installed or removed when the instrument is turned on or off, with one exception:

Do not open the filter cube access doors while a plate is being read (unless Gen5 prompts you to install a particular filter cube). Doing so may affect measurements and result in invalid data.

Open the top or side access door, as appropriate. Gently slide the filter cube into (or out of) its chamber, and then close the door. The barcode label will be scanned upon door closure.

Cleaning Filter Cubes

Do not disassemble filter cubes.

Refer to **Inspect/Clean Filter Cubes** in **Periodic Maintenance**.

Filter Cubes Available for Purchase

Preconfigured barcoded filter cubes are available for purchase. Custom filter cubes are also available. Contact BioTek with any questions.

Custom Filter Cubes



Custom barcoded filter cubes are assembled at BioTek. The information in this section is for reference purposes.

Custom barcoded filter cubes support single PMT measurements. For dual simultaneous-measurement filter cubes (e.g., FP, TR-FRET), please contact BioTek.



The Synergy Neo 2 supports the use of custom filter cubes. Each custom cube has the same ID: 255. It is the user's responsibility to ensure that the correct custom filter cube is installed before performing a read.

The following drawings shows the orientation of the filters inside the cubes.



Fluorescence Intensity Filter Cube: Use Position 1 or 2

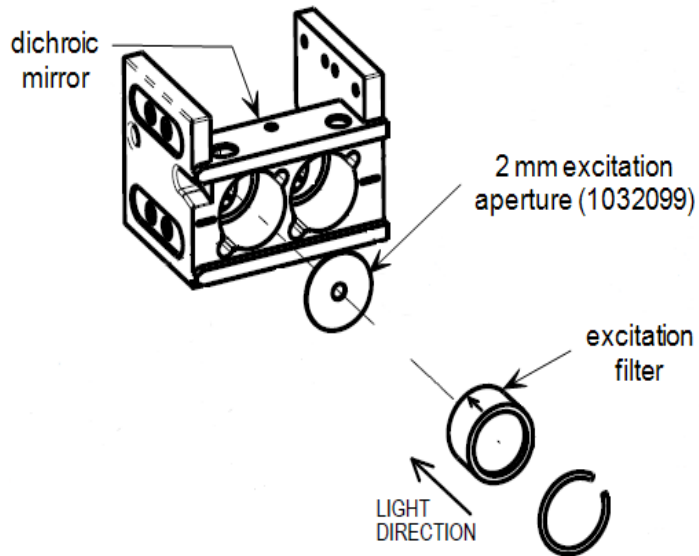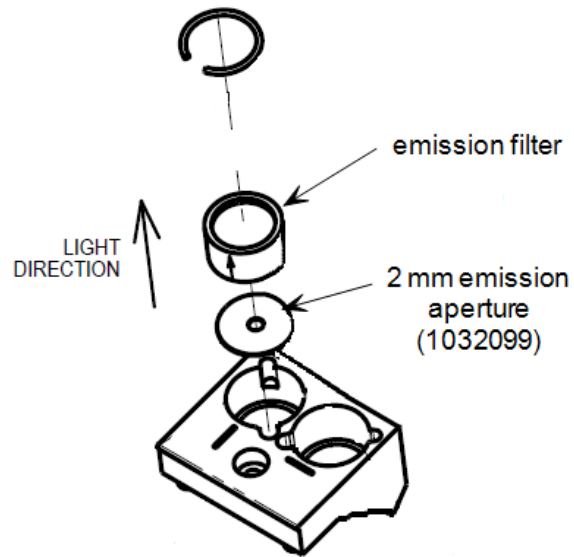



TRF Filter Cube: Use Position 1 or 2

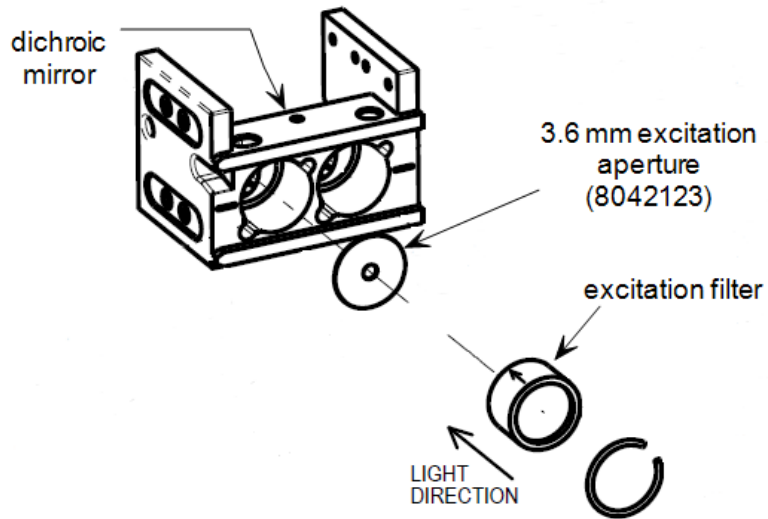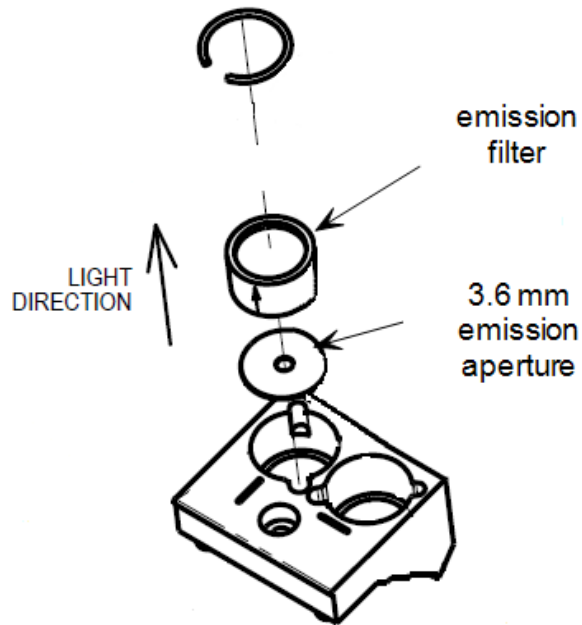

FP Filter Cube: Use both Positions 1 and 2

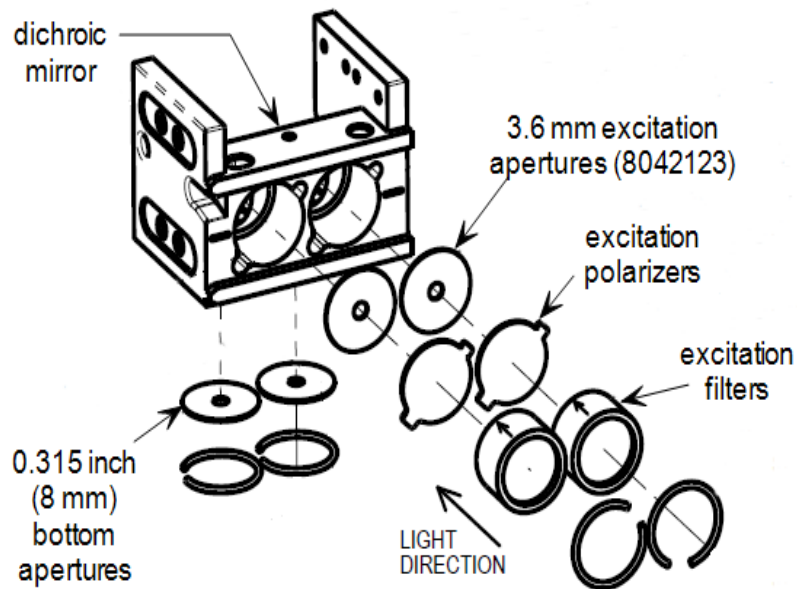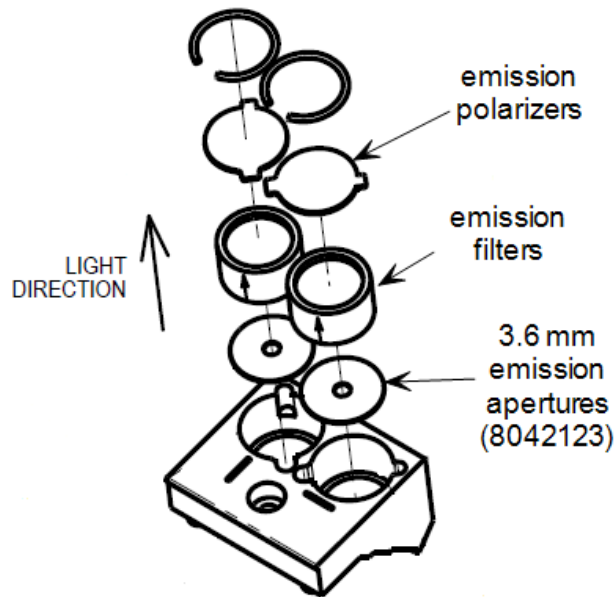

Luminescence Filter Cube: Use Position 1 or 2

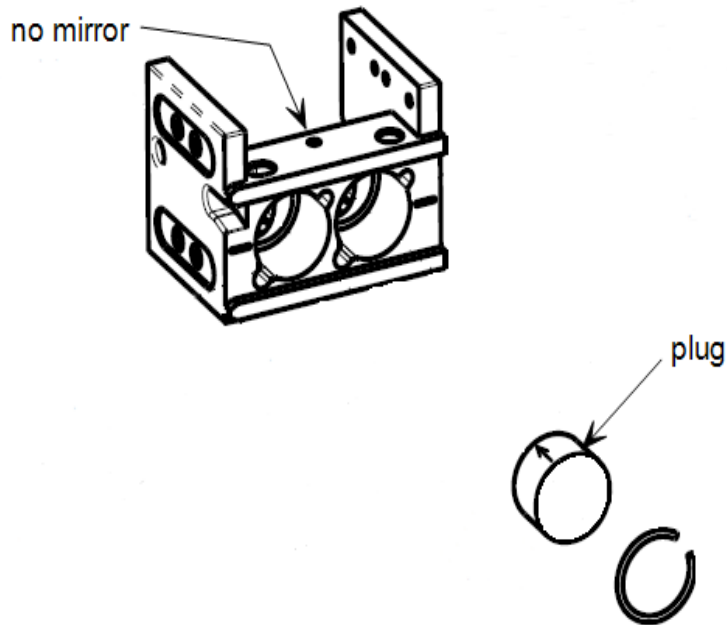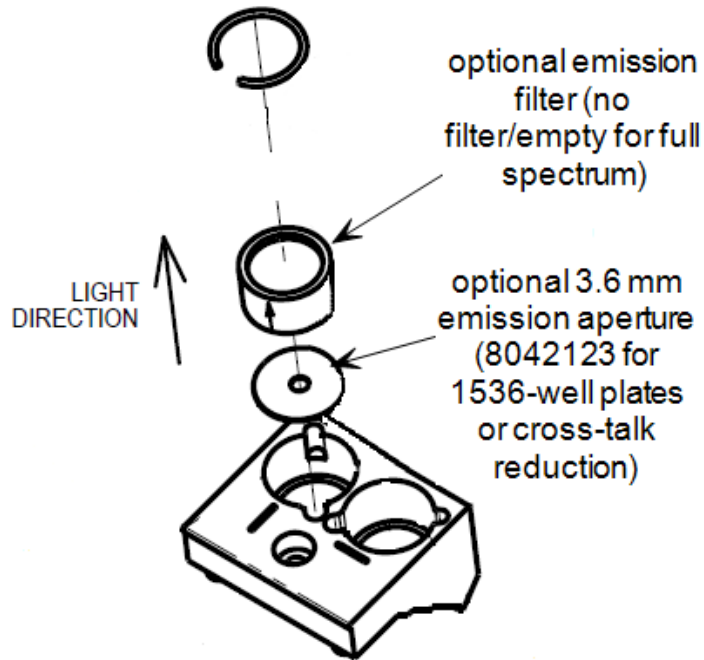

Creating or Modifying Custom Barcoded Filter Cubes in Gen5

Each custom filter cube has the same barcode ID: "255." You must enter the custom filter cube's information in Gen5.

If you define two or more filter cubes with ID 255, it is your responsibility to install the correct filter cube when performing a read. To facilitate this, it is recommended to match the **Cube name** field in Gen5 with the name you write on the cube's label.

To create a new cube:

1. Click **System > Optics Library > Filter Cubes > Create Cube**.
2. Enter a unique Cube Name.
3. Select the detection type for Position 1 (away from the cube's handle), and then click **Edit**.
4. Define information as applicable to the selected detection type. Click **Help** for assistance. Click **OK** with finished.
5. Repeat steps 3 and 4 for Position 2 (closest to the cube's handle), if used.

To modify an existing cube:

1. Click **System > Optics Library > Filter Cubes**.
2. Select the desired cube and click **Modify Cube**.
3. Edit information as necessary.

Shaking System

The user can choose six different shake rates from each of three different shake types: linear, slow orbital, and fast orbital. Each shake rate corresponds to a specific carrier displacement for one or both axes, ranging from 1 mm to 6 mm. A linear shake simply moves the carrier's y-axis back and forth in a line, whereas an orbital shake moves both carrier axes to scribe a circle. Each orbital type may also be selected in a double format, in which the carrier scribes a figure-eight pattern instead.

Carrier Shake Definitions

Shake Type	Displacement (mm/steps)	Period (msecs)	Frequency (Hz)	RPM	Ramp Profile
Linear	1/6	54.8	18.3	1096	13

Shake Type	Displacement (mm/steps)	Period (msecs)	Frequency (Hz)	RPM	Ramp Profile
	2/12	82.0	12.2	731	14
	3/18	105.8	9.5	567	15
	4/24	121.6	8.2	493	16
	5/30	146.3	6.8	410	17
	6/36	166.9	6.0	360	18
Slow orbital	1/6	107.3	9.3	559	19
	2/12	164.4	6.1	365	20
	3/18	212.5	4.7	282	21
	4/24	253.5	3.9	237	22
	5/30	292.1	3.4	205	23
	6/36	334.2	3.0	180	24
Fast orbital	1/6	74.3	13.4	807	25
	2/12	109.6	9.1	548	26
	3/18	141.3	7.1	425	27
	4/24	169.0	5.9	355	28
	5/30	195.4	5.1	307	29
	6/36	222.8	4.5	269	30

Each shake consists of a series of repeated ramped moves for one or both carrier axes. The sinusoidally defined ramps are used for acceleration and deceleration only; the plateau section of the ramp is skipped. When an axis reaches the end of the ramp profile, it simply repeats the ramp in the opposite direction (except for double-orbital shakes, where the y-axis repeats a ramp twice before changing directions).

The circular pattern of the orbital shake is achieved by starting the carrier's y-axis move sequence first. When that axis reaches peak speed at the top of its ramp, the carrier's x-axis begins its own series of repeated moves using the same profile. The axes remain in synch but out of phase, with the maximum speed of one always coinciding with the minimum speed of the other, until the specified shake time has expired.

Maximum Shaking Amplitude Based on Assay Volume and Plate Type

The Synergy Neo 2 offers a broad range of shaking speeds and amplitudes. If the wells of a microplate are almost full, shaking can result in spillage inside the instrument. The following table is designed to help avoid this issue. Select your microplate type and assay volume, and you will get the acceptable shake amplitude on your instrument.

Sample Volume	Linear	Orbital Slow	Orbital Fast	Double Orbital Slow	Double Orbital Fast
6-well Microplate					
0–3 mL	1–6 mm	1–6 mm	1–6 mm	1–6 mm	1–6 mm
3–4 mL	1–6 mm	1–6 mm	1–6 mm	1–6 mm	1 mm
4–5 mL	1 mm	1–6 mm	1–6 mm	1–6 mm	1 mm
5–7 mL	1 mm	1–6 mm	1–3 mm	1 mm	No
7–8 mL	1 mm	No	1–3 mm	No	No
9 mL	1 mm	No	No	No	No

Sample Volume	Linear	Orbital Slow	Orbital Fast	Double Orbital Slow	Double Orbital Fast
12-well Microplate					
0–1.5 mL	1–6 mm	1–6 mm	1–6 mm	1–6 mm	1–6 mm
1.5–2.5 mL	1–2 mm	1–6 mm	1–6 mm	1–6 mm	1 mm
2.5–3 mL	1–2 mm	1–6 mm	No	No	1 mm
3–4 mL	No	1–6 mm	No	No	No

Sample Volume	Linear	Orbital Slow	Orbital Fast	Double Orbital Slow	Double Orbital Fast
24-well Microplate					
0–0.75 mL	1–6 mm	1–6 mm	1–6 mm	1–6 mm	1–6 mm
0.75–1 mL	1 mm	1–6 mm	1–6 mm	1–6 mm	1–6 mm
1–1.5 mL	1 mm	1–6 mm	1–6 mm	1–6 mm	No
1.5–2 mL	1 mm	1–6 mm	No	1–6 mm	No

Sample Volume	Linear	Orbital Slow	Orbital Fast	Double Orbital Slow	Double Orbital Fast
48-well Microplate					
0–0.5 mL	1–6 mm	1–6 mm	1–6 mm	1–6 mm	1–6 mm
0.5–1 mL	1 mm	1–6 mm	1–6 mm	1–6 mm	1–6 mm
1–1.3 mL	No	1–6 mm	1–6 mm	1–6 mm	No

Sample Volume	Linear	Orbital Slow	Orbital Fast	Double Orbital Slow	Double Orbital Fast
96-well Microplate					
0–0.25 mL	1–6 mm	1–6 mm	1–6 mm	1–6 mm	1–6 mm

Injector System


The tubing and injectors should be cleaned at least every three months. See **Periodic Maintenance**, for instructions.

Inspect the injector system daily for leaks, preferably immediately after priming and whenever tubing changes have been made.

If a syringe is leaking, it may need to be replaced. See **As Needed Maintenance**, for instructions.

Dispense Module

The dispense module sits on top of the reader and pumps fluid from the reagent bottles to the injectors located inside the instrument. Fluid is injected into one well at a time. The injectors support plate types from 6- to 384-wells.

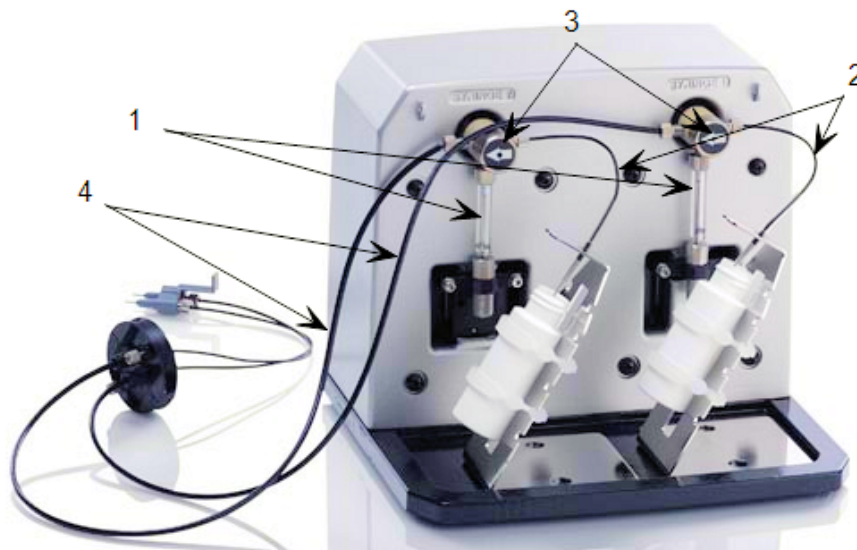

- 1 Two 250 μ L syringes draw fluid from the supply bottles.
- 2 Inlet tubes transport fluid from the supply bottles to the syringes. These tubes are short pieces of opaque PTFE (Teflon) tubing connected to stainless-steel probes on one end and threaded fittings on the other end.

- 3 Solenoid valves allow the fluid drawn from the supply bottles by the syringe pumps to flow into the outlet tubes.
- 4 Outlet tubes transport fluid from the syringes into the instrument, through the tubing ports on the Synergy Neo 2's top cover. The outlet tubes are opaque PTFE tubes with threaded fittings on each end.

Avoid continuous contact with harsh chemicals. Rinse the fluid path with deionized water after contact with any strong acid, base, or solvent.

See the **Periodic Maintenance** chapter for cleaning instructions.

Priming the Injector System

Before running a dispense assay, prime the system with the reagent or dispensing fluid. In addition, tip priming can be performed at the start of an assay and, sometimes, just before each dispense to a well. The tip prime compensates for any fluid loss at the injector tip due to evaporation since the last dispense. All priming activities are controlled via Gen5.

If the injector system is not primed adequately, air bubbles can get trapped in the system and affect injection volumes. Air bubbles in the system can also result in fluid spraying or scattering inside the reader.

Both types of primes require a fluid reservoir to be present on the microplate carrier.


- The priming plate is placed on the microplate carrier for a Prime operation (to prime the dispense system with fluid).
- The tip priming trough is placed in the rear pocket of the carrier, and is used for performing the Tip Prime before dispensing. The trough holds up to 1.5 mL of liquid and must be periodically emptied and cleaned by the user.

Do not perform tip priming when using tall plates. Generally, plates with fewer than 96 wells are too tall for error-free tip priming; and, tip priming is rarely required for these larger-volume plates.

The priming tray should be empty before priming and should contain fluid after priming.

Alpha Laser

WARNING

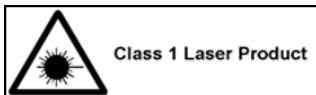

Class 1 Laser Product. Alpha laser "A" and TRF laser "T" models


See "Report on laser safety" on page xi.

"A"-model instruments are equipped with an alpha laser to support Laser Alpha LISA, AlphaScreen® SureFire and similar applications.

Light source	100 mW, 680 nm laser
Detector	PMT
Wavelength selection	Filters (top)
Read speed	96 well: 30 seconds 384 well: 1 minute 50 seconds 1536 well: 7 minutes 20 seconds

TRF Laser

WARNING

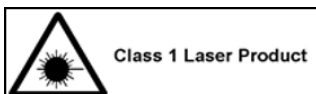

Class 1 Laser Product. Alpha laser "A" and TRF laser "T" models

See "Report on laser safety" on page xi.

"T" models provide a 337 nm nitrogen laser for TR-FRET/HTRF/TRF applications. Gen5 version 3.06 or higher is required. The TRF laser in enhanced sweep mode is more sensitive than the Xenon flash lamp.

384-Well Plates	Xenon Flash	TRF Laser
Detection limit	40 fM	5 fM



Dedicated barcoded filter cubes are required:

PN	Description	Barcode #	Label	Notes
1035018	Cube 18 EX 337 dual PMT	018	TRF Laser: EX 377	For running dual PMT TRF laser assay, e.g. HTRF
1035126	Cube 126 337/620 LUM ASBY	126	337/620: LUM	For running Eu testing during instrument qualification

Gen5 Software

BioTek Gen5 software supports all Synergy Neo 2 reader models.

Note: “T”-models with TRF Laser and “M”-models with dedicated broadband luminescence fiber require Gen5 version 3.06 or higher.

Use Gen5 to control the reader, the dispense module (if equipped), and the stacker (if equipped); perform data reduction and analysis on the measurement values; print or export results; and more. This section provides brief instructions for working with Gen5 to create protocols and experiments and read plates in **Standard Mode** for versions 3.06 and higher. Refer to the Gen5 Help system for more information.



Simple mode can be used after instrument installation is completed. Simple mode is a basic workflow designed for faster set up and execution of experiments; it limits user choices and guides users toward success using interactive dialogs.

Protocols and Experiments

In Gen5, a protocol contains instructions for controlling the reader and (optionally) instructions for analyzing and reporting or exporting the data retrieved from the reader. At a minimum, a protocol must specify the procedure for the assay you wish to run. After creating a protocol, create an experiment that references the protocol. You'll run the experiment to read plates and analyze the data.

These instructions briefly describe how to create a protocol in Gen5. See the Gen5 Help system for complete instructions.

1. Create a new protocol.
2. Open the Procedure dialog. If prompted to select a reader, select the **Synergy Neo 2** and click **OK**.
3. Select a Plate Type.

The Plate Type selected in Gen5 must match the assay plate. Otherwise, the results of the read may be invalid. Also, a collision may result with the transport system because of differing plate dimensions. Be sure to select the **Use lid** check box, if applicable.

4. Add steps to the procedure for shaking or heating the plate, dispensing fluid, reading the plate, and more. Click **Validate** to verify that the reader supports the defined steps, and then click **OK**.

Optionally, perform the next steps to analyze and report the results:

5. Open the Plate Layout dialog and assign blanks, samples, controls, and/or standards to the plate.
6. Open the Data Reduction dialog to add data reduction steps. Categories include Transformation, Well Analysis, Curve Analysis, and Qualitative Analysis/QC.
7. Create a report or export template via the Report/Export Builders.
8. Select **File > Save** and give the file an identifying name.

These instructions briefly describe how to create an experiment and then read a plate in Gen5. See the Gen5 Help system for complete instructions.

1. Create a new experiment using an existing protocol.
2. Select the desired protocol and click **OK**.
3. Select a plate in the menu tree and read it.
4. When the read is complete, measurement values appear in Gen5. Select the desired data set from the Data list.
5. Select **File > Save** and give the file an identifying name.

Dispense Module Control

This section applies only to models with injectors.

Gen5 is used to perform several dispense functions, such as initialize, dispense, prime, and purge. The Prime and Purge functions are introduced here. See the Gen5 Help system for more information.

Priming routines are used to fill the dispense tubing with fluid so that air is not dispensed first, resulting in an incorrect dispense volume. The purging routines are used to clean the fluid paths. See also **Flushing/Purging the Fluid Path** in the **Periodic Maintenance** chapter.

Prime

Before running an experiment with a dispense step, prime the system with the fluid to be used.

1. Place the priming plate on the carrier.
2. Fill the supply bottle with a sufficient volume of the fluid to be used for the prime and the assay. Insert the appropriate inlet tube into the bottle.
3. In Gen5, select **System > Instrument Control > Synergy Neo 2** and click the **Prime** tab.
4. Select the Dispenser number (1 or 2) associated with the supply bottle.
5. Enter the Volume to be used for the prime. The minimum recommended prime volume is 2000 μL .
6. Select a prime Rate, in $\mu\text{L}/\text{second}$.
7. Click **Prime** to start the process.
8. When finished, carefully remove the priming plate from the carrier and empty it.

If the priming plate is empty, the prime volume was too low or there is a problem with the dispense system.

Purge

To save reagent, Gen5 provides the option to purge fluid from the dispense tubing back into the supply bottle.

1. In Gen5, select **System > Instrument Control > Synergy Neo 2** and click the **Prime** tab.
2. Select the Dispenser number (1 or 2) associated with the supply bottle.
3. Enter the desired purge Volume in μL (e.g., 2000).
4. Select a prime Rate in $\mu\text{L}/\text{second}$.
5. Click **Purge** to start the process.

Recommendations for Optimum Performance

General

- Microplates should be clean and free from dust or bottom scratches. Use new microplates from sealed packages. Do not allow dust to settle on the surface of the solution; use microplate covers or seals when not reading the plate. If reading plates with covers still installed, do not forget to compensate for this. See the Gen5 Help.
- Filter solutions to remove particulates that could cause erroneous readings.

- Although the Synergy Neo 2 supports standard flat, U-bottom, and V-bottom microplates, the reader achieves optimum performance with flat-bottomed wells when running in Absorbance mode. See **Appendix A, Specifications** for more information on the supported plates.
- Non-uniformity in the optical density of the well bottoms can cause loss of accuracy, especially with U- and V-bottom polyvinyl microplates. Check for this by reading an empty microplate. Dual wavelength readings can eliminate this problem, or bring the variation in density readings to within acceptable limits for most measurements.
- Inaccuracy in pipetting has a large effect on measurements, especially if smaller volumes of liquid are used. For best results in most cases, use at least 100 μL per well in a 96-well plate and 25 μL in a 384-well plate.
- Pipetting solution into 384-well plates often traps air bubbles in the wells, which may result in inaccurate readings. A dual-wavelength reading method usually eliminates these inaccuracies. For best results, however, remove the air bubbles by degassing the plate in a vacuum chamber or spinning the plate in a centrifuge before reading.
- The inclination of the meniscus can cause loss of accuracy in some solutions, especially with small volumes. Shake the microplate before reading to help bring it within acceptable limits. Use Tween 20, if possible (or some other wetting agent) to normalize the meniscus for absorbance measurements. Some solutions develop menisci over a period of several minutes. This effect varies with the brand of microplate and the solution composition. As the center of the meniscus drops and shortens the light path, the density readings change. The meniscus shape will stabilize over time.
- It is the user's responsibility to understand the volumetric limits of the plate type in use as it applies to the assay being run.
- Use UV transparent microplates for UV wavelength reads.
- Use of liquids with concentrations of acids, corrosives, or solvents of 3 percent and greater can begin attacking the materials inside the instrument's chamber. Running multiple plates with concentrations < 3% in long kinetics may also have a destructive effect. If the experiment is incubated, it will accelerate the deterioration of chamber components. When in doubt about the use of acids, corrosives, or solvents, please contact Technical Support.

Dimethyl sulfoxide (DMSO) vapor can create a removable deposit on optical surfaces, which can trigger instrument self-test errors. Using **DMSO assay concentrations of 2% or below** is recommended. Limit long exposure in kinetic assays or incubated assays when possible. For questions about the use of DMSO, contact Technical Support.

Read Direction

The Synergy Neo 2 performs most reads in a column-wise direction, that is, moving from well A1, to B1, then C1, and so on. The exception to this is for reads that use sweep speed and random well reads, which are read in a row-wise fashion.

Luminescence Measurements

For highly sensitive luminescence assays using white plates, add a Delay step to your procedure to "dark adapt" the plates in the Synergy Neo 2's reading chamber before taking measurements.

Monochromator-Based Fluorescence Systems

Although Time-Resolved Fluorescence can be performed with the monochromator, the filter-based fluorescence system is more sensitive for TRF and is the better choice.

Models with Injectors

- To keep the dispense system in top condition, flush and purge the fluid lines with deionized (DI) water every day or upon completion of an assay run, whichever is more frequent. Some reagents may crystallize or harden after use, clogging the fluid passageways. Flushing the tubing at the end of each day, letting the DI water soak, and then purging the lines at the beginning of each day ensures optimal performance of the dispense system. See the **Periodic Maintenance** chapter for more information.
- When dispensing volumes less than or equal to 20 μL /well, we recommend specifying a tip prime volume that is equal to the dispense volume. For dispense volumes greater than 20 μL /well, we recommend a tip prime volume of 20 μL .
- To avoid spillage and possible contamination of the instrument, empty the tip prime trough frequently and do not exceed the total fluid volume of the plate well when dispensing.

Using 384-Well Microplates

When using a 384-well microplate, you can use the Gen5 Auto Map feature to ensure you are using an accurate plate map for your reads. See the Gen5 Help for more information.

Incubation and Partial Plates

When performing a partial plate read that includes an incubation step, the following recommendations can reduce the effects of evaporation of your samples:

- Use microplate lids.
- Fill unused wells with fluid.
- Cluster your sample wells rather than spacing them throughout the plate.
- Place your sample wells in the center of the plate. This placement may lead to less evaporation than if you place the samples in wells on the edge of the plate.

Alpha Laser Detection

When using alpha detection, ensure that the relative humidity remains below 85% and the ambient temperature is less than 30°C.

Periodic Maintenance

This chapter provides instructions for maintaining the Synergy Neo 2 and external dispense module (if used) in top condition, to ensure that they continue to perform to specification.

Periodic Maintenance Overview	62
Recommended Maintenance Schedule	62
Warnings and Precautions	63
Clean Exposed Surfaces	64
Inspect/Clean Filter Cubes	65
Flush/Purge the Fluid Path	65
Run a Dispense Protocol (Optional)	67
Empty/Clean the Tip Priming Trough	68
Clean the Priming Plate	68
Clean the Dispense Tubes and Injectors	68

Periodic Maintenance Overview

A general maintenance regimen for all Synergy Neo 2 models includes periodically cleaning all exposed surfaces and, except for NEO2MON models, inspecting, and if necessary, cleaning the filter cubes.

For models with the external dispense module, additional tasks include flushing/purging the fluid path and cleaning the tip prime trough, priming plate, supply bottles, dispense tubing, and injectors.

Daily Cleaning for the Dispense Module

To ensure accurate performance and a long life for the dispense module and injectors, flush and purge the fluid lines with deionized (DI) water every day or after completing an assay run, whichever is more frequent. Some reagents may crystallize or harden after use and clog the fluid passageways. Take special care when using molecules that are active at very low concentrations (e.g., enzymes, inhibitors). Remove any residual reagent in the dispense lines using a suitable cleaning solution (review the reagent's package insert for specific recommendations).

Flushing the tubing at the end of each day, letting the DI water soak overnight, and then purging the lines at the beginning of each day ensures optimal performance of the dispense system. BioTek recommends performing a visual inspection of the dispense accuracy before running an assay protocol that includes dispense steps.

BioTek also recommends flushing the module with DI water before conducting the decontamination procedure described in the **As Needed Maintenance** chapter.

Models with injectors: Accumulated algae, fungi, or mold may require decontamination. See the **As-Needed Maintenance** chapter for complete decontamination instructions.

Recommended Maintenance Schedule

This table recommends maintenance tasks and the frequency with which each task should be performed.

The risk and performance factors associated with your assays may require that some of all of the maintenance procedures be performed more frequently than shown here.

Task	Daily	Quarterly	As Needed
All models:			
PM-1 Clean exposed surfaces			✓
PM-2 Inspect/clean filter cubes*		✓	
Decontamination	<i>before shipment or storage</i>		
Models with injectors:			
PM-3 Flush/purge the fluid path	✓		
PM-4 (Optional) Run Dispense protocol			✓
PM-5 Empty/clean tip prime trough	✓		
PM-6 Clean priming plate			✓

* except NEO2MON models

Warnings and Precautions

Read the following before performing any maintenance procedures:

WARNING **Internal Voltage.** Always turn off the power switch and unplug the power supply before cleaning the outer surface of the instrument.

WARNING **Liquids.** Avoid spilling liquids on the instrument; fluid seepage into internal components creates a potential for shock hazard or instrument damage. If a spill occurs while a program is running, stop the program and turn off the instrument. Wipe up all spills immediately. Do not operate the instrument if internal components have been exposed to fluid.

CAUTION **Liquids.** Do not immerse the instrument, spray it with liquid, or use a dripping-wet cloth on it. Do not allow water or other cleaning solution to run into the interior of the instrument. If this happens, contact Technical Support.

CAUTION **Lubricants.** Do not apply lubricants to moving parts. Lubricant on components in the carrier compartment will attract dust and other particles, which may cause the instrument to produce an error.

WARNING

Potential Biohazards. Wear protective gloves when handling contaminated instruments. Gloved hands should be considered contaminated at all times; keep gloved hands away from eyes, mouth, nose, and ears.

Mucous membranes are considered prime entry routes for infectious agents. Wear eye protection and a surgical mask when there is a possibility of aerosol contamination. Intact skin is generally considered an effective barrier against infectious organisms; however, small abrasions and cuts may not always be visible. Wear protective gloves when handling contaminated instruments.

WARNING

The instrument with all available modules weighs approximately **100 lbs. (45.36 kg)**. Use two people when lifting and carrying the instrument.

Clean Exposed Surfaces

Exposed surfaces may be cleaned (not decontaminated) with a cloth moistened (not soaked) with water or water and a mild detergent. You'll need:

- Deionized or distilled water
 - Clean, lint-free cotton cloths
 - Mild detergent (optional)
1. Turn off and unplug the instrument.
 2. Moisten a clean cotton cloth with water, or with water and mild detergent, then thoroughly wring it out so that liquid does not drip from it. **Do not soak the cloth.**
 3. Wipe the plate carrier and all exposed surfaces of the instrument.
 4. Wipe all exposed surfaces of the dispense module (if used).
 5. If detergent was used, wipe all surfaces with a cloth moistened with water.
 6. Use a clean, dry cloth to dry all wet surfaces.

Models with injectors: If the Tip Priming Trough overflows or other spills occur inside the instrument, wipe the carrier and the surface beneath the carrier with a dry cotton cloth. The internal chamber and probes are not customer-accessible. If overflow is significant, contact Technical Support with any questions about your particular model.

Inspect/Clean Filter Cubes

* except NEO2MON-N models

Ambient laboratory air is used to cool the flash bulb, and the filter cubes can become dusty as a result. Filters should be inspected and cleaned at least every three months. You'll need:

- Isopropyl, ethyl, or methyl alcohol
- 100% pure cotton balls or high-quality lens-cleaning tissue
- Cloth gloves
- Magnifying glass

Do not touch the filters with your bare fingers!

1. Open the access door on the front of the instrument, and the side door, if equipped with a bottom filter cube. Slide the filter cubes out of their compartments.
2. Inspect the glass filters for speckled surfaces or a "halo" effect. This may indicate deterioration due to moisture exposure over a long period of time.

If you have any concerns about the quality of the filters, contact Technical Support.

3. Using cotton balls or lens-cleaning tissue moistened with a small amount of high-quality alcohol, clean each filter by lightly stroking its surface in one direction.
4. Use a magnifying glass to inspect the surface; remove any loose threads left from the cotton ball.
5. Replace the filter cubes and close the door.

Flush/Purge the Fluid Path

Applies only to Synergy Neo 2 models with injectors.

At the end of each day that the dispense module is in use, flush the fluid path using the Gen5 priming utility. Leave the fluid to soak overnight or over a weekend, and then purge the fluid before using the instrument again.

This flushing and purging routine is also recommended before disconnecting the outlet tubes from the reader to prevent a spill. It is required before decontamination to remove any assay residue prior to applying isopropyl alcohol or sodium hypochlorite.

To flush the fluid path:

1. Fill two supply bottles with deionized or distilled water. Insert the supply (inlet) tubes into the bottles.
2. Place the priming plate on the carrier.
3. Select **Prime/Dispense** from the Instrument Control panel or from the Gen5 main screen, select **System > Instrument Control > Synergy Neo 2**.
4. Click the **Prime** tab and select **Dispenser 1**.
5. Set the Volume to **5000 µL**. Keep the default prime rate.
6. Click **Prime** to start the process. When the process is complete, carefully remove the priming plate from the carrier and empty it.
7. Repeat the process for Dispenser 2.

Leave the water in the system overnight or until the instrument will be used again. Purge the fluid from the system (see below) and then prime with the dispense reagent before running an assay.

To purge the fluid from the system:

1. Place the inlet tubes in empty supply bottles or a beaker.
2. Select **Prime/Dispense** from the Instrument Control panel or from the main menu **System > Instrument Control > Synergy Neo 2**.
3. Click the **Prime** tab and select **Dispenser 1**.
4. Set the Volume to **2000 µL**.
5. Click **Purge** to start the process.
6. When the purge is complete, repeat the process for Dispenser 2.

After purging the system, you may wish to run a quick Dispense protocol to visually verify the dispense accuracy (see below) or the more thorough Dispense Accuracy and Precision Tests (see **Instrument Qualification**).

Run a Dispense Protocol (Optional)

Applies only to Synergy Neo 2 models with injectors.

After flushing/purging the system and before running an assay that requires dispense, visually verify the dispense function.

1. Create a Dispense protocol in Gen5:
 - a. Create a new protocol with the plate type set to match the plate you will use.
 - b. Add a Dispense step with the following parameters:
 - Select Dispenser **1**.
 - Set Tip Priming to **Before this dispense step** and Volume to **10 µL**.
 - Set the Dispense Volume to **100 µL** (or an amount to match your assay protocol).
 - Adjust the Rate to support the dispensing volume.
 - Click **OK** to close the dialog and add the Dispense step to the procedure.
 - c. Add another Dispense step with the same parameters; select Dispenser **2**.
 - d. Add a Read step with the following parameters (Gen5 requires a Read step in a Dispense protocol):
 - Select any Detection Method.
 - Set the Read Type to **Endpoint**.
 - Click **Full Plate**, click **Clear All**, then select well **A1**. Click **OK**.
 - Select any wavelength or define one Filter Set.
 - Click **OK** to close the dialog and add the Read step to the procedure.
 - e. Click **OK** to close the procedure.
 - f. Select **File > Save** and give the protocol an identifying name, such as "Dispense Observation."
2. Fill the reagent bottles with a DI water–Tween solution (e.g., add 1 mL Tween 20 to 1000 mL of deionized water).
3. Create a new experiment using the "Dispense Observation" protocol.
4. Click **Read** and follow the prompts.

5. When the procedure is complete, visually assess the fluid level in the wells for accuracy. If the well volume appears to be unevenly distributed, clean the internal dispense tubes and injectors.

Empty/Clean the Tip Priming Trough

Applies only to Synergy Neo 2 models with injectors.

The tip priming trough is a removable cup located in the rear pocket of the microplate carrier, used for performing the Tip Prime. The trough holds about 1.5 mL of liquid and must be periodically emptied and cleaned by the user. Gen5 will instruct you to do this at the start of an experiment that requires dispensing.

1. Extend the microplate carrier and carefully remove the tip priming trough from the carrier.
2. Wash the trough in hot, soapy water. Use a small brush to clean in the corners.
3. Rinse the trough thoroughly and allow it to dry completely.
4. Replace the trough in the microplate carrier.

Clean the Priming Plate

Applies only to Synergy Neo 2 models with injectors.

Clean the priming plate regularly to prevent bacteria growth and residue buildup. Wash the plate in hot, soapy water, using a small brush to clean in the corners. Rinse thoroughly and allow it to dry completely.

Clean the Dispense Tubes and Injectors

Applies only to Synergy Neo 2 models with injectors.

The Synergy Neo 2's dispense tubes and injectors require routine cleaning, at least quarterly and possibly more frequently depending on the type of fluids dispensed.

Required Materials

- Protective gloves
- Safety glasses
- Mild detergent

- Clean, lint-free cotton cloths
- Deionized or distilled water
- Stylus (stored in a plastic cylinder affixed to the rear of the dispense module or reader) (PN 2872304)

Remove the Dispense Tubes and Injector Holders

1. Open the door on the front of the reader.
2. Grasp the injector tip holder by the tab and pull it up out of its socket.
3. Using your fingers, remove the thumbscrews securing the light shield to the top of the reader and slide the shield up the outlets tubes.
4. Slide the injector tip holder through the hole in the top of the reader.
5. Turn each tube's thumbscrew counterclockwise and gently pull each tube from its injector tip.
6. On the dispense module, turn each outlet tube's thumbscrew counterclockwise to disconnect it from the injector.

Clean the Dispense Tubes and Injectors

Some reagents can crystallize and clog the tubing and injectors. Daily flushing and purging can help to prevent this, but more rigorous cleaning may be necessary if reagent has dried in the tubing or injectors.

To clean the dispense tubes, soak them in hot, soapy water to soften and dissolve any hardened particles. Flush each tube by holding it vertically under a stream of water.

To clean the injectors:

1. Gently insert the stylus into each injector tip to clear any blockages. (The stylus is stored in a plastic cylinder affixed to the rear of the dispense module.)
2. Stream water through the pipe to be sure it is clean. If the water does not stream out, try soaking in hot, soapy water and then reinserting the stylus.

Be careful not to bend the injector tips. A bent tip might not dispense accurately.

As-Needed Maintenance

This chapter contains maintenance and component-replacement procedures that need to be performed only occasionally.

Decontamination	71
Procedure for Models without a Dispenser	72
Procedure for Models with a Dispenser	73
Dispenser Syringe Replacement	76

Decontamination

- The Synergy Neo 2 requires decontamination prior to shipping, storage, and disposal.
- Decontamination is required by the U.S. Department of Transportation regulations.
- Persons performing the decontamination process must be familiar with the basic setup and operation of the instrument.
- BioTek recommends the use of the following decontamination solutions and methods based on our knowledge of the instrument and recommendations of the Centers for Disease Control and Prevention (CDC). Neither BioTek nor the CDC assumes any liability for the adequacy of these solutions and methods. Each laboratory must ensure that decontamination procedures are adequate for the biohazard(s) they handle.

WARNING **Internal Voltage.** Always turn off the power switch and unplug the power supply before cleaning the outer surface of the instrument.

WARNING The sodium hypochlorite (NaClO) solution is caustic; wear gloves and eye protection when handling the solution.

CAUTION Do not immerse the instrument, spray it with liquid, or use a dripping-wet cloth. Do not allow the cleaning solution to run into the interior of the instrument. If this happens, contact Technical Support.

WARNING Wear prophylactic gloves when handling contaminated instruments. Gloved hands should be considered contaminated at all times; keep gloved hands away from eyes, mouth, and nose. Eating and drinking while decontaminating instruments is not advised.

WARNING Mucous membranes are considered prime entry routes for infectious agents. Wear eye protection and a surgical mask when there is a possibility of aerosol contamination. Intact skin is generally considered an effective barrier against infectious organisms; however, small abrasions and cuts may not always be visible. Wear protective gloves when performing the decontamination procedure.

Required Materials

For all Synergy Neo 2 models:

- Sodium hypochlorite (NaClO, or bleach)
- 70% isopropyl alcohol (as an alternative to bleach)
- Deionized or distilled water

- Safety glasses
- Surgical mask
- Protective gloves
- Lab coat
- Biohazard trash bags
- 125-mL beakers
- Clean, lint-free cotton cloths

Additional materials for models with the dispense module:

- Phillips screwdriver
- Small brush for cleaning the tip priming trough and priming plate
- (Optional) Mild detergent

Procedure for Models without a Dispenser

1. Turn off and unplug the instrument.
2. Prepare an aqueous solution of 0.50% sodium hypochlorite (bleach). If the effects of bleach are a concern, 70% isopropyl alcohol may be used.

Check the percent NaClO of the solution you are using. Commercial products are typically 10.0% NaClO; prepare a 1:20 dilution. Household products are typically 5.0% NaClO; prepare a 1:10 dilution.

3. Wet a cloth with the NaClO solution or alcohol, then thoroughly wring it out so that liquid does not drip from it. Do not soak the cloth.
4. Open the plate carrier door and slide out the plate carrier.
5. Wipe the plate carrier and all exposed surfaces of the instrument.
6. Wait 20 minutes. Moisten a cloth with deionized (DI) or distilled water and wipe all surfaces of the instrument that have been cleaned with the bleach solution or alcohol.
7. Use a clean, dry cloth to dry all wet surfaces.
8. Reassemble the instrument as necessary.

9. Discard the used gloves and cloths using a biohazard trash bag and an approved biohazard container.

Procedure for Models with a Dispenser

Perform the Routine Procedure when the Synergy Neo 2 is functioning normally. If you are unable to perform a prime due to a system failure, perform the Alternate Procedure described on page 76.

Routine Procedure

If disinfecting with sodium hypochlorite (NaClO), be sure to flush repeatedly with deionized water to remove the NaClO. After disinfecting with sodium hypochlorite, perform the rinse procedure provided on page 74.

If disinfecting with alcohol, do not immediately prime with deionized water, because the drying effect of the alcohol is an important aspect of its disinfectant properties.

Clean Exposed Surfaces

1. Turn off and unplug the instrument.
2. Prepare an aqueous solution of 0.50% sodium hypochlorite (bleach). If the effects of bleach are a concern, 70% isopropyl alcohol may be used.

Check the percent NaClO of the solution you are using. Commercial products are typically 10.0% NaClO; prepare a 1:20 dilution. Household products are typically 5.0% NaClO; prepare a 1:10 dilution.

3. Open the plate carrier door and slide out the plate carrier.
4. Wet a cloth with the bleach solution or alcohol, then thoroughly wring it out so that liquid does not drip from it. Do not soak the cloth.
5. Wipe the plate carrier and the exposed surfaces of the dispenser.
6. Wait 20 minutes. Moisten a cloth with deionized (DI) or distilled water and wipe all surfaces that have been cleaned with the bleach solution or alcohol.
7. Use a clean, dry cloth to dry all wet surfaces.

8. If the dispenser is installed, purge any fluid (see **Flush/Purge the Fluid Path** on page 65) and detach the outlet tubes from the instrument. If it is not installed, attach only the dispenser's communication cable to the instrument. Remove the supply bottles and their holders.
9. Perform the decontamination procedures described below through page 75.

Decontaminate the Fluid Lines

1. Place a beaker with 20 mL of 0.5% sodium hypochlorite solution or 70% isopropyl alcohol near SYRINGE 1 on the dispenser.
2. Place the SYRINGE 1 inlet tube in the beaker.
3. If you have not already done so, detach the dispenser's outlet tubes from the instrument. Place the ends of the outlet tubes in an empty beaker and set the beaker next to the dispenser.
4. Launch Gen5 and from the main screen select **System > Instrument Control**, and click the **Prime** tab.
5. Select Dispenser **1**, enter a Volume of **5000 μ L**, and keep the default dispense Rate.
6. Place the priming plate on the carrier.
7. Run two prime cycles, for a total of 10,000 μ L.
8. Wait at least 20 minutes to allow the solution to disinfect the tubing.
9. Remove the inlet tube from the beaker of disinfectant solution.
10. From the Instrument Control dialog, change the Volume to 1000 μ L.
11. Run one prime cycle, to flush the disinfectant out of the fluid lines.
12. Empty the beaker containing the outlet tubes. Put the tubes back in the empty beaker.
13. If sodium hypochlorite (bleach) was used, perform Rinse the Fluid Lines.

Otherwise (or after performing the Rinse procedure), repeat steps 1–13 for SYRINGE 2/Dispenser 2.

Rinse the Fluid Lines

Perform this procedure only if decontamination was performed using sodium hypochlorite.

1. Place a beaker containing at least 30 mL of deionized water on the dispenser.
2. Place the SYRINGE 1 or 2 inlet tube in the beaker.
3. If you have not already done so, place the outlet tubes in an empty beaker.
4. From the Instrument Control dialog, select Dispenser **1** or **2**, set the Volume to **5000 μ L**, and keep the default dispense Rate.
5. Run five prime cycles, for a total of 25,000 μ L.
6. Pause for 10 minutes and then run one prime cycle with 5000 μ L. This delay will allow any residual sodium hypochlorite to diffuse into the solution and be flushed out with the next prime.
7. Empty the beaker containing the outlet tubes.
8. Wipe all surfaces with deionized water.
9. Discard the used gloves and cloths using a biohazard trash bag and an approved biohazard container.

Clean the Tubing and Injectors

Perform the procedures under **Clean the Dispense Tubes and Injectors** in **Periodic Maintenance**.

Decontaminate the Tip Priming Trough and Priming Plate

1. Remove the tip priming trough from the instrument's microplate carrier.
2. Wash the tip priming trough and priming plate in hot, soapy water. Use a small brush or cloth to clean the corners of the trough and plate.
3. To decontaminate, soak the trough and plate in a container of 0.5% sodium hypochlorite or 70% isopropyl alcohol for at least 20 minutes.
 - If decontaminating in a bleach solution, thoroughly rinse the trough and plate with DI water.
 - If decontaminating with alcohol, let the trough and plate air dry.
4. Discard the used gloves and cloths using a biohazard trash bag and an approved biohazard container.

Alternate Procedure

If you are unable to prime the Synergy Neo 2 due to a system failure, decontaminate the instrument and the dispenser as follows:

1. Perform the procedures under **Clean the Dispense Tubes and Injectors** in **Periodic Maintenance**.
2. Prepare an aqueous solution of 0.50% sodium hypochlorite (bleach). If the effects of bleach are a concern, 70% isopropyl alcohol may be used.

Check the percent NaClO of the solution you are using. Commercial products are typically 10.0% NaClO; prepare a 1:20 dilution. Household products are typically 5.0% NaClO; prepare a 1:10 dilution.

3. Slide the microplate carrier out of the instrument.
4. Wet a cloth with the bleach solution or alcohol, then thoroughly wring it out so that liquid does not drip from it. Do not soak the cloth.
5. Use the cloth to wipe:
 - All exterior surfaces of the instrument
 - All surfaces of the plate carrier
 - The exposed surfaces of the dispenser, including the syringe valves
6. Remove the tubing and the syringes from the dispenser and soak them in the bleach or alcohol solution. Wait for 20 minutes.
7. Moisten a cloth with DI or distilled water and wipe all surfaces that have been cleaned with the bleach solution or alcohol.
8. Rinse all tubing and the syringes with DI water.
9. Use a clean, dry cloth to dry all surfaces on the instrument and the dispenser.
10. Reassemble the dispenser as necessary.
11. Discard the used gloves and cloths using a biohazard trash bag and an approved biohazard container.

Dispenser Syringe Replacement

Refer to the **Periodic Maintenance** chapter for cleaning procedures you must perform regularly and also in the case of poor performance (for example, when Dispense Accuracy

and Precision tests fail). If cleaning the dispenser does not eliminate performance problems, or if a syringe is obviously leaking, perform these instructions to replace a faulty syringe. Contact Technical Support to order replacement syringes.

To change a syringe, first use Gen5 to put the syringe in its maintenance position.

Syringe Maintenance Position

Do not change the syringe position or calibrate the dispensers unless instructed to do so as part of installation, upgrade, or maintenance.

Gen5 provides access to syringe setup functions for maintenance and calibration purposes. When a syringe needs to be installed or replaced, it must first be moved to its “maintenance position.”

1. From the Gen5 main screen, select **System > Instrument Control > Synergy Neo 2** and click the **Prime** tab.
2. Select the appropriate Dispenser number (**1** or **2**) associated with the syringe.
3. Click **Maintenance**. The syringe plunger will move to its furthest-from-home position. The syringe can then be disconnected from the drive bracket and unscrewed from the valve.

Replace the Syringe

After using Gen5 to move the syringe into its maintenance position:

1. Using your fingers, unscrew the bottom thumbscrew that secures the syringe, underneath the bracket. Retain this bottom thumbscrew; it is needed for the replacement syringe.
2. Unscrew the top thumbscrew to disengage the syringe from the valve.
3. Remove the new syringe from its protective box. (The syringe should already be assembled in one piece; if it is not, see “Install the Dispenser” in the **Installation** chapter.
4. Hold the syringe vertically with the threaded end at the top. Screw the top of the syringe into the bottom of the syringe valve. Finger-tighten only.
5. Carefully pull down the bottom of the syringe until it rests inside the hole in the bracket.
6. Pass the thumbscrew (used to hold the old syringe) up through this hole and thread it into the bottom of the syringe. Hold the syringe from rotating while tightening the thumbscrew. Finger-tighten only.

7. From the Gen5 main screen, select **System > Instrument Control > Synergy Neo 2**. Click the **Prime** tab and click **Initialize**.

Instrument Qualification

This chapter describes the tests that BioTek has developed for complete qualification of all models of the Synergy Neo 2. This chapter introduces the various test methods, describes the materials and relevant Gen5 protocols used to execute the tests, explains how to analyze test results, and provides troubleshooting tips in the event of a failure.

Instrument Qualification Overview	80
IQ/OQ/PQ Description	80
Recommended Qualification Schedule	81
System Test	83
Plate Shaker Test	84
Absorbance Testing Overview	85
Absorbance Plate Test	86
Absorbance Liquid Tests	91
Fluorescence Testing	98
Luminescence Test	137
Dispense Module Tests	143
Alpha Detection Test	151

Instrument Qualification Overview

This chapter contains the recommended Installation Qualification (IQ), Operational Qualification (OQ), and Performance Qualification (PQ) procedures for all models of the Synergy Neo 2 Microplate Reader.

Every Synergy Neo 2 reader and external dispense module is fully tested prior to shipment and should operate properly upon initial setup. If you suspect that a problem occurred during shipment, if you have received the equipment after returning it to the factory for service, and/or if regulatory requirements dictate that you qualify the equipment on a routine basis, perform the procedures outlined in this chapter.

A Product Qualification Package (PN 1350526N) for the Synergy Neo 2 is available for purchase. The package contains complete procedures, Gen5 protocols, checklists, and logbooks for performing Installation Qualification, Operational Qualification, Performance Qualification, and Periodic Maintenance.

IQ/OQ/PQ Description

Installation Qualification confirms that the reader and its components have been supplied as ordered and ensures that they are assembled and configured properly for your lab environment.

- The recommended IQ procedure consists of setting up the instrument and its components as described in **Chapter 2, Installation**, and performing the System Test. For models with injectors, a quick “Injector Test” is also performed, to ensure that the dispense module is properly installed and there are no leaks.
- The IQ procedure should be performed initially (before the reader is used for the first time).
- The successful completion of the IQ procedure verifies that the instrument is installed correctly. The Operational Qualification procedure should be performed immediately following the successful IQ.

Operational Qualification confirms that the equipment operates according to specification initially and over time.

- The recommended OQ procedure consists of performing the system test, Absorbance Plate Test, luminescence test, a series of liquid tests, and, if the external dispense module is used, the Dispense Accuracy and Precision Tests.
- The OQ procedure should be performed initially (before first use) and then routinely; the recommended interval is annually. It should also be performed after any major repair or upgrade to the hardware or software.

- Although out-of-tolerance failures will be detected by the OQ tests, results should be compared with those from the routine Performance Qualification tests and previous OQ tests to monitor for trends.
- The successful completion of the OQ procedure, in combination with results that are comparable to previous PQ and OQ tests, confirms that the equipment is operating according to specification initially and over time.

Performance Qualification confirms that the reader consistently meets the requirements of the tests performed at your laboratory.

- The recommended PQ procedure consists of performing the System Test, Absorbance Plate Test, Luminescence test, a series of Liquid Tests, and, if the external dispense module is used, the Dispense Accuracy and Precision Tests.
- Your facility's operating policies may also require that you routinely perform an actual assay, to confirm that the reader will consistently give adequate results for the assays to be run on it.
- These tests should be performed routinely; the recommended interval is monthly or quarterly, depending on the test. This frequency may be adjusted depending on the trends observed over time.
- The successful completion of the PQ procedure confirms that the equipment is performing consistently under normal operating conditions.

Recommended Qualification Schedule

This table defines the recommended intervals for qualification for an instrument used two to five days a week. The schedule assumes that the Synergy Neo 2 is properly maintained as outlined in the **Periodic Maintenance** section.

The risk and performance factors associated with your assays may require that the Operational and Performance Qualification procedures be performed more or less frequently than shown here.

Tasks/Tests	IQ	OQ	PQ	
	Initially	Initially/ Annually	Monthly	Quarterly
All models:				
Unpacking, installation, setup, and verification	✓			

Tasks/Tests	IQ	OQ	PQ	
	Initially	Initially/ Annually	Monthly	Quarterly
System Test	✓	✓	✓	
Absorbance capability:				
Absorbance Plate Test		✓	✓	
Absorbance Liquid Test 1 or Liquid Test 2**		✓		✓
(Optional) Absorbance Liquid Test 3 or Absorbance Plate Test at 340 nm with PN 7260551		✓		✓
Fluorescence capability:				
Corners, Sensitivity, Linearity Tests		✓	✓	
Fluorescence Polarization Test*		✓		✓
Time-Resolved Fluorescence Test -Flash*		✓		✓
Time-Resolved Fluorescence Test -Laser*		✓		✓
Luminescence capability:				
Luminescence Test		✓	✓	
Dispenser system*:				
Injection System Test	✓			
Dispense Accuracy/Precision Test		✓		✓
Alpha Laser capability*:				
Alpha Detection Test		✓		

* If applicable to your reader model

** If you have Absorbance Test Plate PN #7260522, perform Liquid Test 1. Otherwise, perform Liquid Test 2.

System Test

Each time the Synergy Neo 2 is turned on, it automatically performs a series of tests on the reader's motors, lamp, the PMTs, and various subsystems. The duration of this system test depends on the reader model, and can take a few minutes to complete. If all tests pass, the microplate carrier is ejected and the green LED on the carrier switch remains on.

If any test results do not meet the internally coded Failure Mode Effects Analysis (FMEA) criteria, the reader beeps repeatedly and the red LED on the carrier switch flashes. If this occurs, press the carrier eject button to stop the beeping. If necessary, initiate another system test using Gen5 to try to retrieve an error code from the reader. Refer to **Appendix B, Error Codes** for information on error codes and for troubleshooting tips.

If the power-up system test fails, when you initiate a system test using Gen5, Gen5 displays a message stating that the reader has a pending system test report. Click **OK** in the message box to review the report; it contains information obtained up to the point of the failure.

1. Turn on the reader and launch Gen5.
2. If your assays use incubation, we recommend enabling Temperature Control and allowing the incubator to reach its set point before running the system test. To access this feature, from the Gen5 main screen, select **System > Instrument Control** and click the **Pre-Heating** tab. Wait until the display indicates that the temperature set point has been reached.
3. Select **System > Diagnostics > Run System Test**.

If the test fails during execution, a message box appears in the software. Close the box; the test report contains the error code that was generated by the failure.

4. When the test is complete, a dialog appears, requesting additional information. Enter your user name and other information (if desired) and then click **OK**.
5. The results report appears. Scroll down toward the bottom of the report; it shows either "SYSTEM TEST PASS" or "SYSTEM TEST FAIL *** ERROR (error code) DETECTED."
6. Print the report if desired.
 - Gen5 stores the results in a database, so the results can be retrieved at any time. We recommend that you print and save the reports to document that the test was performed. You can also save the report as an HTML file.

7. If the test failed, look up the error code in **Appendix B, Error Codes** to determine its cause. If the cause is something you can fix, turn off the reader, fix the problem, and then turn the reader back on and retry the test.

If the test continues to fail, or if the cause is not something you can fix, contact Technical Support.

Plate Shaker Test

This test verifies that the multispeed plate shaker is operating properly. The test involves creating and running a protocol with shaking enabled for a duration of 30 seconds. The sound of the carrier shaking is all that needs to be confirmed to verify that the plate shaker is operating properly.

Absorbance Testing Overview

Tests for the absorbance system use a combination of solid-state Absorbance Test Plates and liquid plates. The test plates and the materials used for creating the liquid plates are available for purchase.

To qualify the absorbance system for the Synergy Neo 2, you should perform:

- Absorbance Liquid Test 1 *and* Absorbance Plate Test (using 7260522) *or*
- Absorbance Liquid Test 2

Optionally, to qualify operation in the UV range, you should also perform:

- Absorbance Liquid Test 3 *or* Absorbance Plate Test at 340 nm (using 7260551)

Absorbance Liquid Tests

Test Methods

Absorbance Liquid Test 1 confirms repeatability and alignment of the reader when a solution is used in the microplate. If these tests pass, then the lens placement and optical system cleanliness are proven. For the Repeatability portion of this test, two columns containing a color-absorbing solution are read five times at 405 nm. For each well, an "allowed deviation" is determined based on its Mean OD and the reader's repeatability specification. Each well's Standard Deviation must be less than its Allowed Deviation to pass. To confirm the reader's mechanical alignment, the plate is rotated 180 degrees in the carrier (e.g., A1 is now in the H12 position) and the same two columns are read. The initial and new OD readings are compared, using the reader's accuracy specification. If the two readings in the same well do not meet specification, the reader may be out of alignment.

If an Absorbance Test Plate is not available, **Absorbance Liquid Test 2** may be conducted to test the instrument's alignment, repeatability, and accuracy by preparing a series of solutions of varying OD values. "Absorbance Liquid Test 2" on page 93.

Absorbance Liquid Test 3 is an optional test offered for those sites that must have proof of linearity at 340 nm. (Alternatively, the BioTek 340 nm Absorbance Test Plate may be used.) This test is optional since the instrument has good "front-end" linearity throughout the specified wavelength range. While the absolute values of the OD cannot be determined by this test, the results will indicate if there is adequate repeatable absorbance and a linear slope. This method is dependent upon proper dye dilution and a skilled pipetting technique. It is expected that the first dilution (mid-level solution) will have an absorbance value near 75% of that of the stock (high-level) solution, and that the second dilution (low-level solution) will have an absorbance value near 50% of that of the stock solution.

Absorbance Plate Test

This test uses Absorbance Test Plate (PN 7260522) to confirm the mechanical alignment; optical density accuracy, linearity, and repeatability; and wavelength accuracy of the Synergy Neo 2. The Absorbance Test Plate compares the reader's optical density and wavelength measurements to NIST-traceable values.

An alternate method for confirming accuracy, linearity, and repeatability is Liquid Test 2, described on page 93.

To run this test, you need Absorbance Test Plate (PN 7260522), with its accompanying data sheet.

- The Absorbance OD Standards section contains NIST-traceable standard OD values for the filters at several different wavelengths. We recommend testing at six wavelengths—those at or close to the wavelengths used in your assays.
- The Wavelength Accuracy Standards section contains Expected Peak wavelength values for the filter in position C6 on the plate. Each value has a valid test range and expected peak range recorded on the accompanying data sheet. For example, an Expected Peak value may be 586 nm with tolerance values of $-6/+4$ (or a test range of 580 to 590 nm).

Absorbance Test Plate PN 7260551 can be used to confirm optical density accuracy, linearity, and repeatability at 340 nm and is offered as an alternative to conducting Absorbance Liquid Test 3.

The instructions provided below are guidelines. Refer to the Gen5 Help system for more information.

Define Absorbance Test Plate Parameters

1. Obtain the certificates that came with the Test Plate.
2. Start Gen5 and from the main screen select **System > Diagnostics > Test Plates > Add/Modify Plates**.
3. Click **Add**. The Absorbance Test Plate dialog appears.
4. Select the appropriate Plate Type and enter the plate's serial number.

5. Enter the Last Certification and Next Certification dates from the calibration sticker on the Test Plate.
6. If the wavelength values in the top row of the grid in Gen5 are appropriate for your tests, enter the OD values from the data sheet into the grid. Make sure you enter the correct value for each well/wavelength combination.

If you need to add, change, or delete any wavelength values, click **Wavelength List**. Click the Gen5 **Help** button for assistance.

7. Select the number of Peak Wavelength tests to run (up to 4), based on the number of expected peak wavelength values provided on the certificate.
8. Enter the Expected Peak value(s) from the certificate and set the Test Range – and + values.

The glass in each test plate is unique. Therefore, the expected peaks may differ slightly from plate to plate.

- If the C6 filter is Holmium or Erbium glass, the certificate contains two Spectral Bandpass tables. For wavelengths greater than 285 nm, we recommend performing the test using the 5.0 nm table. For wavelengths in the 230–285 nm range, we recommend using the 2.4 nm table.
 - For the Erbium glass, any peak can be used.
 - For the Holmium glass, use the expected peak values closest to 242, 279, 362, 417, and 538 nm. For example, if your certificate looks like the one below, you might choose to run the test at four of the five highlighted Expected Peak/Test Range combinations:

2.4 nm Spectral Bandpass	
Expected Peak	Test Range
241	–5+5
278	–6+4
287	–4+6
334	–5+5

5.0 nm Spectral Bandpass	
Expected Peak	Test Range
242	–5+5
279	–6+4
288	–4+6
334	–5+5

360	-5+5	362	-5+5
417	-5+5	417	-5+5
484	-5+5	485	-5+5
537	-5+5	538	-5+5
642	-5+5	643	-5+5

- If your C6 filter is Didymium glass, a single peak wavelength value is provided. Enter this value and set the Test Range – and + values so the range displayed in parentheses is 580 to 590, as shown here:

Peak wavelength tests: ▼

Expected Peak - Test Range +

 ▼ ▼ (580 to 590)

9. Review all the values you entered. Click **OK** to save the data.

The information you just entered is available on Gen5 each time the Absorbance Plate Test is performed. It may need to be modified after the annual recertification of your test plate.

Run the Absorbance Plate Test

1. From the Gen5 main screen, select **System > Diagnostics > Test Plates > Run**. If prompted, select the desired Test Plate and click **OK**.
2. When the Absorbance Test Plate Options dialog appears, select **Perform Peak Wavelength Test** if it is not already selected.
3. Highlight the wavelength(s) to be included in this test.

Select only those wavelengths most appropriate for your use of the reader.

4. (Optional) Enter any Comments.
5. Click **Start Test**.
6. Place the Test Plate in the microplate carrier so that well A1 is in the right-rear corner of the carrier.

7. Click **OK** to run the test.
8. When the test is completed, the results report appears. Scroll through the report; every result should show "PASS". Because Synergy Neo 2 absorbance specifications extend only to 2.5 OD, any test results above this number do not state "PASS" or "FAIL," but display only the numerical results. This information is for reference only and does not imply any guarantee of accuracy above 2.5 OD. See **Results and Troubleshooting Tips** for information on results and troubleshooting tips in the event of failures.

Gen5 stores the results in a database; they can be retrieved any time. We recommend you print and save the report to document that the test was performed.

Results and Troubleshooting Tips

The Absorbance Test Plate Report contains results for the following:

- **Peak Absorbance:** When the test is performed, the C6 filter is scanned at the test range(s) defined by the user in the Absorbance Test Plate dialog. To verify wavelength accuracy, the wavelength of the maximum absorbance is compared with the peak wavelength value entered in the software, which comes from the Peak Wavelength Certificate supplied with the Test Plate. The accuracy of the wavelength should be ± 3 nm (± 2 nm instrument, ± 1 nm filter allowance). If the reader fails this test:
 - Make sure the information entered into Gen5 matches the Test Plate's Peak Wavelength Certificate.
 - Verify that the Test Plate has a filter in location C6. (Test Plates with the part number 9000547 and 7260551 do not have this filter.) This test plate was designed for use with readers that use interference filters to generate the specific wavelength required. The Peak Absorbance filter in C6 is used for testing the monochromator.
 - Check the C6 filter to make sure it is clean. If needed, clean it with lens paper. Do not remove the filter from the Test Plate, and do not use alcohol or other cleaning agents.
 - Verify that the Test Plate is within its calibration certification period. If it is out of date, contact BioTek to schedule a recertification.
 - Check the microplate carrier to ensure it is clear of debris and that the plate is held securely in place by the retainer when the carrier moves into the reader.
- **Alignment:** This test measures the alignment of the microplate carrier with the optical path. A reading greater than 0.015 OD represents an out-of-alignment condition.

Wells A1, A12, H1, and H12 are the only valid alignment holes for the reader on the PN 7260522 Test Plate.

If the reader fails this test:

- Ensure that the Test Plate is correctly seated in the microplate carrier.
- Check the four alignment holes (A1, A12, H1, H12) to ensure they are clear of debris.
- Check the microplate carrier to ensure it is clear of debris and, again, that the plate is held securely in place by the retainer when the carrier moves into the reader.

- **Accuracy:** Accuracy is a measure of the optical density of Test Plate wells C1, D4, E2, F5, G3, and H6 as compared with known standard values contained in the Standards Certificate that accompanies each Test Plate.

If the reader fails this test:

- Verify that the filter calibration values entered in Gen5 are the same as those on the Test Plate's Standards Certificate.
- Check the filters on the Test Plate to ensure they are clean. If necessary, clean them with lens paper. Do not remove the filters from the test plate, and do not use alcohol or other cleaning agents.
- Verify that the Test Plate is within its calibration certification period. If it is out of date, contact Technical Support to schedule a recertification.

- **Repeatability:** Repeatability is a measure of the instrument's ability to read the same well with minimum variation between two reads with the well in the same location. If the reader fails this test:

- Check the filters on the Test Plate to ensure there is no debris that may have shifted between readings and caused changes.
- Check the microplate carrier to ensure it is clear of debris.

Linearity of the optical density readings is confirmed by default if the optical density readings are accurate. To further verify this, you can perform a regression analysis on the Test Plate OD values in a spreadsheet program such as Microsoft Excel. An R Squared value of at least 0.9900 is expected.

Absorbance Liquid Tests

Conducting liquid tests confirms the reader's ability to perform to specification with liquid samples. Liquid testing differs from testing with the Absorbance Test Plate in that liquid in the wells has a meniscus, whereas the Test Plate's neutral density glass filters do not. The optics characteristics may differ in these two cases, thus alerting the operator to different types of problems.

Absorbance Liquid Test 1

Absorbance Liquid Test 1 confirms repeatability and alignment of the reader when a solution is used in the microplate. If these tests pass, then the lens placement and optical system cleanliness are proven.

Materials

Manufacturer part numbers are subject to change.

- New 96-well, clear, flat-bottom microplate (Corning Costar #3590 recommended)
- Stock Solution A or B, which may be formulated by diluting a dye solution available from BioTek (A) or from the materials listed below (B).

Solution A

- BioTek QC Check Solution No. 1 (PN 7120779, 25 mL; PN 7120782, 125 mL)
 - Deionized water
 - 5-mL Class A volumetric pipette
 - 100-mL volumetric flask
1. Pipette a 5-mL aliquot of BioTek QC Check Solution No. 1 into a 100-mL volumetric flask.
 2. Add 95 mL of DI water; cap and shake well. The solution should measure approximately 2.000 OD when using 200 μ L in a flat-bottom microwell.

Solution B

- Deionized water
- FD&C Yellow No. 5 dye powder (typically 90% pure)
- Tween 20 (polyoxyethylene (20) sorbitan monolaurate) or BioTek wetting agent (PN 7773002) (a 10% Tween solution)

- Precision balance with capacity of 100 g minimum and readability of 0.001 g
 - Weigh boat
 - 1-liter volumetric flask
1. Weigh out 0.092 g of FD&C Yellow No. 5 dye powder into a weigh boat.
 2. Rinse the contents into a 1-liter volumetric flask.
 3. Add 0.5 mL of Tween 20, or 5 mL of BioTek's wetting agent.
 4. Fill to 1 liter with DI water; cap and shake well. The solution should measure approximately 2.000 OD when using 200 μ L in a flat-bottom microwell.

Prepare the Plate

Be sure to use a new microplate, because fingerprints or scratches may cause variations in readings.

1. Using freshly prepared stock solution (Solution A or B), prepare a 1:2 dilution using deionized water (one part stock, one part deionized water; the resulting solution is a 1:2 dilution).
2. Pipette 200 μ L of the concentrated solution (A or B) into the first column of wells in the microplate.
3. Pipette 200 μ L of the diluted solution into the second column of wells.

After pipetting the test solution into the microplate and before reading the plate, we strongly recommend shaking the plate for four minutes. This will allow any air bubbles in the solution to settle and the meniscus to stabilize. Alternatively, wait 20 minutes after pipetting the test solution before reading the plate.

Read the Plate

1. Using Gen5, read the microplate **five times** at 405 nm using the Normal read mode, single wavelength, no blanking. Save the data after each read ("Normal" plate position).
2. Without delay, rotate the microplate 180 degrees so that well A1 is in the "H12" position. Read the plate **five more times**, saving the data after each read ("Turnaround" plate position).

3. Print out the ten sets of raw data, or export them to an Excel spreadsheet.
See the calculations described below to obtain results for the tests performed.

Analyze the Results

1. The plate is read five times in the “Normal” position at 405 nm. Calculate the Mean OD and Standard Deviation of those five reads for each well in columns 1 and 2.
2. For each well in columns 1 and 2, calculate the Allowed Deviation using the repeatability specification for a 96-well plate: $\pm 1\% \pm 0.005$ OD from 0.000 to 2.000 OD (Mean * 0.010 + 0.005). For each well, its standard deviation should be less than its allowed deviation.

Example: Five readings in well A1 of 0.802, 0.802, 0.799, 0.798, and 0.801 result in a mean of 0.8004 and a standard deviation of 0.0018. The mean multiplied by 1.0% (0.8004 * 0.010) equals 0.008, and when added to 0.005 equals 0.013; this is the allowed deviation for well A1. Since the standard deviation for well A1 is less than 0.013, the well meets the test criteria.

3. The plate is read five times in the “Turnaround” position at 405 nm. Calculate the Mean OD of those reads for each well in columns 11 and 12.
4. Perform a mathematical comparison of the Mean values for each microwell in its Normal and Turnaround positions (that is, compare A1 to H12, A2 to H11, B1 to G12,... H2 to A11). To pass the test, the differences in the compared mean values must be within the accuracy specification for a 96-well microplate: $\pm 1.0\% \pm 0.010$ OD from 0.000 to 2.000 OD.

Example: If the mean value for well A1 in the Normal position is 1.902 with a specified accuracy of $\pm 1.0\% \pm 0.010$ OD, then the expected range for the mean of the well in its Turnaround (H12) position is 1.873 to 1.931 OD. $1.902 \times 0.010 + 0.010 = 0.029$; $1.902 - 0.029 = 1.873$; $1.902 + 0.029 = 1.931$.

Repeatability Specification:

- $\pm 1.0\% \pm 0.005$ OD from 0.000 to 2.000 OD
- $\pm 3.0\% \pm 0.005$ OD from 2.000 OD to 2.500 OD

Accuracy Specification:

- $\pm 1.0\% \pm 0.010$ OD from 0.000 to 2.000 OD
- $\pm 3.0\% \pm 0.010$ OD from 2.000 OD to 2.500 OD

Absorbance Liquid Test 2

The recommended method for testing the instrument’s alignment, repeatability, and accuracy is to use the Absorbance Test Plate. If the Test Plate is not available, however, Liquid Test 2 can be used for these tests.

Materials

- A new 96-well, clear, flat-bottom microplate (Corning Costar #3590 is recommended)
- Ten test tubes, numbered consecutively, set up in a rack
- Calibrated hand pipette (Class A volumetric pipette recommended)
- Solution A or B (see the instructions for Liquid Test 1)
- A 0.05% solution of deionized water and Tween 20

Prepare the Dilutions

Create a percentage dilution series, beginning with 100% of the original concentrated stock solution (A or B) in the first tube, 90% of the original solution in the second tube, 80% in the third tube, all the way to 10% in the tenth tube. Dilute using the 0.05% solution of deionized water and Tween 20. This solution can also be made by diluting the wetting agent 200:1.

Test Tube Dilutions for Liquid Test 2

Tube Number	1	2	3	4	5	6	7	8	9	10
Volume of Original Concentrated Solution (mL)	20	18	16	14	12	10	8	6	4	2
Volume of 0.05% Tween Solution (mL)	0	2	4	6	8	10	12	14	16	18
Absorbance expected if original solution is 2.0 at 200 μ L	2.0	1.8	1.6	1.4	1.2	1.0	0.8	0.6	0.4	0.2

The choice of dilutions and the absorbance of the original solution can be varied. Use this table as a model for calculating the expected absorbances of a series of dilutions, given a different absorbance of the original solution.

Prepare the Plate

- Pipette 200 μ L of the concentrated solution from Tube 1 into each well of the first column, A1 to H1, of a new flat-bottom microplate.
- Pipette 200 μ L from each of the remaining tubes into the wells of the corresponding column of the microplate (Tube 2 into wells A2 to H2, Tube 3 into wells A3 to H3, and so on).

Linearity and Repeatability Tests

1. Using Gen5, read the microplate prepared above five times using Normal mode, dual wavelength at 450/630 nm. Save the data after each read.

Do not discard the plate; you will use it for the Alignment test.

2. Print out the five sets of Delta OD data, or export them to an Excel spreadsheet.
3. Calculate the results for Linearity:
 - Calculate the mean absorbance for each well, and average the means for each concentration.
 - Perform a regression analysis on the data to determine if there is adequate linearity.

Since it is somewhat difficult to achieve high pipetting accuracy when conducting linear dilutions, an R Square value of at least 0.9900 is considered adequate.

4. Calculate the results for Repeatability:
 - Calculate the mean and standard deviation for the five readings taken in Step 1 at each concentration. Only one row of data needs to be analyzed.
 - For each mean below 2.000 OD, calculate the allowed deviation using the repeatability specification for a 96-well plate of $\pm 1.0\% \pm 0.005$ OD. If above 2.000 OD, apply the $\pm 3.0\% \pm 0.005$ specification.
 - The standard deviation for each set of readings should be less than the allowed deviation.

Example: Absorbance readings of 1.950, 1.948, 1.955, 1.952, and 1.950 will result in a mean of 1.951, and a standard deviation of 0.0026. The mean (1.951) multiplied by 1.0% (1.951×0.010) = 0.0195, which, when added to the 0.005 ($0.0195 + 0.005$) = 0.0245 OD, which is the allowed deviation. Since the standard deviation is less than this value, the reader meets the test criteria.

Repeatability Specification:

$\pm 1.0\% \pm 0.005$ OD from 0.000 to 2.000 OD

$\pm 3.0\% \pm 0.005$ OD from 2.000 OD to 2.500 OD

Alignment Test

1. Using the plate prepared for the Linearity Test on the previous page, conduct a Turnaround test by reading the plate five times with the A1 well in the H12 position. Save the data after each read.

This test results in values for the four corner wells that can be used to determine alignment.

2. Calculate the means of the wells A1 and H1 in the Normal plate position (data from Linearity Test) and in the Turnaround position (from Step 1).
3. Compare the mean reading for well A1 to its mean reading when in the H12 position. Next, compare the mean values for the H1 well to the same well in the A12 position. The difference in the values for any two corresponding wells should be within the accuracy specification for the instrument.

Example: If the mean of well A1 in the normal position is 1.902, where the specified accuracy is $\pm 1.0\% \pm 0.010$ OD, then the expected range for the mean of the same well in the H12 position is 1.873 to 1.931 OD. ($1.902 \times 1.0\% = 0.019 + 0.010 = 0.029$, which is added to and subtracted from 1.902 for the range.)

If the four corner wells are within the accuracy range, the reader is in alignment.

Accuracy Specification:

$\pm 1.0\% \pm 0.010$ OD from 0.000 to 2.000 OD

$\pm 3.0\% \pm 0.010$ OD from 2.000 OD to 2.500 OD

Absorbance Liquid Test 3

Absorbance Liquid Test 3 is provided for sites requiring proof of linearity at 340 nm. This test is optional because the Synergy Neo 2 has good "front end" linearity throughout its wavelength range. As an alternative, the 340 nm Absorbance Test Plate (PN 7260551) may be used for this test.

Materials

Manufacturer part numbers are subject to change.

- New 96-well, clear, flat-bottom microplate (Corning Costar #3590 recommended); alternatively, a UV transparent microplate may be used
- Calibrated hand pipette(s)
- Beakers and graduated cylinder

- Precision balance with readability to 0.01 g
- Buffer solution described below

Buffer Solution

- Deionized water
 - Phosphate-Buffered Saline (PBS), pH 7.2–7.6, Sigma tablets, #P4417 (or equivalent)
 - β -NADH Powder (β -Nicotinamide Adenine Dinucleotide, Reduced Form) Sigma bulk catalog number N 8129, or preweighed 10-mg vials, Sigma number N6785-10VL (or BioTek PN 98233). Store the powder according to the guidelines on its packaging.
1. Prepare a PBS solution from the Sigma tablets.
 2. In a beaker, mix 50 mL of the PBS solution with 10 mg of the β -NADH powder and mix thoroughly. This is the **100% Test Solution**.
 3. (Optional) Read a 150- μ L sample of the solution at 340 nm; it should be within 0.700 to 1.000 OD. If low, adjust up by adding more powder. Do not adjust if slightly high.

Prepare the Plate

1. Prepare the **75% Test Solution** by mixing 15 mL of the 100% Test Solution with 5 mL of the PBS Solution.
2. Prepare the **50% Test Solution** by mixing 10 mL of the 100% Test Solution with 10 mL of the PBS Solution.
3. Carefully pipette the three solutions into a **new** 96-well microplate:
 - 150 μ L of the 100% Test Solution into all wells of columns 1 and 2
 - 150 μ L of the 75% Test Solution into all wells of columns 3 and 4
 - 150 μ L of the 50% Test Solution into all wells of column 5 and 6

Read the Plate

1. Using Gen5, read the microplate **five times** using Normal mode, single wavelength at 340 nm, no blanking. Save the data after each read.
2. Print out the five sets of raw data, or export them to an Excel spreadsheet.

Analyze the Results

The plate is read five times at 340 nm.

1. For each well, calculate the Mean OD and Standard Deviation of the five readings.
2. For each mean calculated in step 1, calculate the allowed deviation using the repeatability specification for a 96-well plate: $\pm 1\% \pm 0.005$ OD from 0.000 to 2.000 OD ($\text{Mean} \times 0.010 + 0.005$). For each well, its standard deviation should be less than its allowed deviation.

Example: Five readings in well A1 of 0.802, 0.802, 0.799, 0.798, and 0.801 result in a mean of 0.8004 and a standard deviation of 0.0018. The mean multiplied by 1.0% (0.8004×0.010) equals 0.008, and when added to 0.005 equals 0.013; this is the allowed deviation for well A1. Since the standard deviation for well A1 is less than 0.013, the well meets the test criteria.

3. Calculate the results for Linearity:
 - For each of the three Test Solutions, calculate the average Mean OD for the wells containing that solution (mean of wells A1 to H2, A3 to H4, and A5 to H6).
 - Perform a regression analysis on the data to determine if there is adequate linearity. The three average Mean OD values are the “Y” values. The solution concentrations are the “X” values (1.00, 0.75, 0.50).

Since it is somewhat difficult to achieve high pipetting accuracy when conducting linear dilutions, an R Square value of at least 0.9900 is considered adequate.

Repeatability Specification:

$\pm 1.0\% \pm 0.005$ OD from 0.000 to 2.000 OD

$\pm 3.0\% \pm 0.005$ OD from 2.000 OD to 2.500 OD

Fluorescence Testing

This section applies to all models except NEO2MON.

Fluorescence Testing Overview

Two options are provided for testing the Synergy Neo2 Multi-Mode Microplate Reader

fluorescence system. One uses a solid state Fluorescence Test Plate (package PN 1400501¹). The other uses liquid plates.

Kits containing the microplates and solutions required for the Liquid Tests are available for purchase; see Materials for Conducting Liquid Tests on page 5.

Fluorescence Test Plate

The Fluorescence Test Plate simplifies the process for conducting fluorescence intensity, fluorescence polarization, and time-resolved fluorescence qualification tests on the Synergy Neo2 Multi-Mode Microplate Reader. The test plate is solid and, therefore, immune to pipetting errors, evaporation issues, and costs experienced with conventional liquid tests.

The test plate package includes Gen5 protocols designed specifically for use with the test plate. The protocols include embedded Microsoft Excel spreadsheets to automatically calculate results and determine pass/fail. The protocols and their spreadsheets were fully validated in accordance with BioTek Product Validation policies and procedures.

The package also contains a user guide that describes the test methods, helps you get started using the plate, and provides important information for cleaning and maintaining the test plate. The guide also provides troubleshooting tips and information on the annual recalibration program.

Results Analysis

Refer to the *Fluorescence Test Plate User Manual* for descriptions of the data reduction calculations for each test. The tests must meet the following criteria to pass:

Fluorescence Intensity (FI) Tests	
Corners	%CV < 3.0
Linearity	R ² >= 0.9500
Sensitivity, filter-based system:	
Top optics, Sodium Fluorescein analogue	Detection Limit <= 5.0 pM
Bottom optics, Sodium Fluorescein analogue	Detection Limit <= 5.0 pM
Top optics, Methylumbelliferone analogue	Detection Limit <= 160.0 pg/mL
Bottom optics, Methylumbelliferone analogue	Detection Limit <= 160.0 pg/mL
Sensitivity, monochromator-based system:	

¹Fluorescence Test Plate PN 7092092 cannot be used for these tests.

Fluorescence Intensity (FI) Tests	
Top optics, Sodium Fluorescein analogue	Detection Limit \leq 20.0 pM
Bottom optics, Sodium Fluorescein analogue	Detection Limit \leq 20.0 pM
Top optics, Methylumbelliferone analogue	Detection Limit \leq 160.0 μ g/mL
Bottom optics, Methylumbelliferone analogue	Detection Limit \leq 160.0 μ g/mL
TRF - FP tests:	
Time-Resolved Fluorescence - Flash Test	Detection Limit \leq 250.0 fM
Time-Resolved Fluorescence - TRF Laser Test	Detection Limit \leq 40 fM
Fluorescence Polarization (FP) Test (filter-based system, dual PMT)	HPR Polarization $>$ 340 mP LPR Standard Deviation $<$ 5
Fluorescence Polarization (FP) Test (filter-based system, single PMT)	HPR Polarization $>$ 340 mP LPR Standard Deviation $<$ 5

Fluorescence Liquid Tests

BioTek has developed a series of liquid tests for verifying fluorescence performance.

Fluorescence Capability	Applicable Liquid Test(s)
Filter-Based Fluorescence Intensity	Corners, Sensitivity, Linearity
Monochromator-Based Fluorescence Intensity	Corners, Sensitivity, Linearity
Fluorescence Polarization (for both single and dual PMT)	"FP"
Time-Resolved Fluorescence	"TRF Flash" or "TRF Laser"

- The **Corners Test** uses fluorescence compounds to verify that the plate carrier is properly aligned in relation to the fluorescence probe(s).
- The **Sensitivity Test** uses a fluorescence compound and buffer solution to test the fluorescence reading capability of the instrument. The ability to detect specific compounds at the required limit of detection ensures that the filters, optical path, and PMT(s) are all in working order. This test verifies that the difference between the concentration well under investigation and the mean of the median buffer well is statistically distinguishable.

- The **Linearity Test** verifies that the system is linear, that is, the signal changes proportionally with changes in concentration (R2 value). Proving that the system is linear allows the Sensitivity Test to be run on two points instead of using serial dilutions.
- The **FP Test** verifies the ability of the instrument to measure polarization of the solution. It verifies the polarizers are installed in the proper orientation and the mechanism is in proper order.
- The **TRF Test** verifies the performance of the xenon flash or TRF laser, and that the filters, optical path, and PMTs are all in working order.

The tests presented in this section require specific microplates, solutions, wavelengths, mirrors, and filters. Your laboratory may require a deviation from some of these tests. For example, you may wish to use a different fluorescing solution or microplate.

1. If deviation from the tests as presented in this section is required, the following steps should be taken the first time each test is run:
2. Perform the tests exactly as described on the following pages.
3. Rerun the tests using your particular solutions, filter cubes, microplates, and so on. If results are comparable, then the results from these tests will be your baseline for future tests.
4. Document your new test procedure(s), and save all test results.

Required Materials

Kits containing the microplates and solutions required for the Liquid Tests are available for purchase; see Materials for Conducting Liquid Tests on page 5.

Microplates should be perfectly clean and free from dust or bottom scratches. Use new microplates from sealed packages.

Manufacturer part numbers are subject to change.

All Tests:

- Deionized or distilled water
- Various beakers, graduated cylinders, and pipettes

- 95% ethanol (for cleaning clear-bottom plates)
- Aluminum foil
- (Optional, but recommended) 0.45-micron filter
- (Optional) Black polyethylene bag(s) to temporarily store plate(s)
- Gen5 protocols listed in the next table and described starting on page 103

For the Filter-Based Fluorescence System	
Synergy Neo2_FI_T_SF.prt	Corners, Sensitivity, and Linearity tests for top optics, using sodium fluorescein
Synergy Neo2_FI_B_SF.prt	Corners, Sensitivity, and Linearity tests for bottom optics, using sodium fluorescein
Synergy Neo2_FI_T_MUB.prt	Alternative top optics test, using methylumbelliferone
Synergy Neo2_FI_B_MUB.prt	Alternative bottom optics test, using methylumbelliferone
Synergy Neo2_Dual PMT_FP.prt	Fluorescence Polarization test using dual PMT
Synergy Neo2_Single PMT_FP.prt	Fluorescence Polarization test using single PMT
Synergy Neo2_TRF.prt	Time-Resolved Fluorescence -flash test
Synergy Neo 2_TRF-N2.prt	Time-Resolved Fluorescence -laser test

For the Monochromator-Based Fluorescence System	
Synergy Neo2_M_FI_T_SF.prt	Corners, Sensitivity, and Linearity tests for top optics, using sodium fluorescein
Synergy Neo2_M_FI_B_SF.prt	Corners, Sensitivity, and Linearity tests for bottom optics, using sodium fluorescein
Synergy Neo2_M_FI_T_MUB.prt	Alternative top optics test, using methylumbelliferone

Corners/Sensitivity/Linearity Tests

Manufacturer part numbers are subject to change.

The materials listed here are for use with Sodium Fluorescein. Methylumbelliferone can be used as an alternate or supplemental method for performing these tests. See page 126.

Buffer:

- NIST-traceable Sodium Borate Reference Standard (pH 9.18) (e.g., Fisher-Scientific 1 L Sodium Borate Mfr. #159532, or equivalent), **or**
- Phosphate-Buffered Saline (PBS), pH 7.2–7.6 (e.g., Sigma tablets, Mfr. #P4417, or equivalent) and pH meter or pH indicator strips with pH range 4 to 10
- Sodium Fluorescein Powder (1 mg vial, BioTek PN 98155)
- **If testing both Top and Bottom optics:** A new, clean 96-well glass-bottom Greiner SensoPlate (Mfr. #655892), a clean Hellma Quartz 96-well titration plate (Mfr. #730.009.QG), or equivalent
- **Top optics:** A new, clean 96-well solid black microplate, such as Corning Costar #3915, or equivalent
- Filter cubes:
 - Upper Top: 3 or equivalent (dual PMT instruments)
 - Lower Top: 107 or equivalent (485/528)

Fluorescence Polarization (FP) Test

- A new, clean, 96-well solid black microplate
- The recommended test solutions are available from Invitrogen Corporation in their “FP One-Step Reference Kit” (PN P3088) or from BioTek (PN 7160014). This kit includes:
 - (Green) Polarization Reference Buffer, 15 mL
 - Green Low Polarization Reference, 4 mL
 - Green High Polarization Reference, 4 mL

The Invitrogen kit also includes two red polarization solutions; these are not used.

- Filter cubes:
 - Upper Top: 4 or equivalent (dual FP)
 - Lower Top: 61 or equivalent (FP 485/528) (dual FP)
 - 108 or equivalent (single PMT)

Time-Resolved Fluorescence (TRF) Test

- 15-mL conical-bottom, polypropylene sample tube
- Filter cubes:

Flash test	TRF laser test
<ul style="list-style-type: none"> ■ Cube 3 or equivalent (dual PMT instruments) ■ Cube 112 or equivalent (360/620) 	<ul style="list-style-type: none"> ■ Cube 126: 337/620; LUM ■ Cube 3 or equivalent

- A new, clean 96-well solid white microplate
- The recommended test solution (FluoSpheres carboxylate-modified microspheres, 0.2 μm europium luminescent, 2 μL) is available from Invitrogen Corporation (PN F20881) or from BioTek (PN 7160011)

Fluorescence Test Solutions

Corners/Sensitivity/Linearity Tests Solution

If using the BioTek sodium fluorescein powder (PN 98155), be sure to hold the vial upright and open it carefully; the material may be concentrated at the top. If a centrifuge is available, spin down the tube before opening.

When diluting the sodium fluorescein powder in buffer, it takes time for the powder to completely dissolve. Allow the solution to dissolve for five minutes, with intermittent vortexing, before preparing the titration dyes.

Wrap the vial containing the stock solution in foil to prevent exposure to light. Discard any open, unused solution after seven days.

1. The Sodium Borate solution does not require further preparation; proceed to step 2. If you are using PBS, prepare the solution:
 - (Optional, but recommended) Using a 0.45-micron filter, filter 200 mL of deionized or distilled water.
 - Follow the manufacturer's instructions on the PBS packaging to create 200 mL, dissolving the necessary amount of PBS into the filtered water.
 - Stir the solution (preferably using a stir table) until the PBS is completely dissolved.
 - Check the pH; it should be between 7.2 and 7.6 at 25°C.
2. Prepare the sodium fluorescein stock solution:
 - Add 2.0 mL of the buffer solution to the 1 mg Sodium Fluorescein (SF) vial. This yields a 1.3288 mM stock solution.
 - Ensure that the dye has completely dissolved and is well mixed.
3. Carefully prepare the dilutions. Label each with "SF" and the concentration:

Mix This SF Solution:	With Buffer:	To Make:	
0.53 mL of 1.3288 mM stock solution	13.47 mL	50.2 μ M	
110 μ L of 50.2 μ M SF	13.89 mL	400 nM	
3.5 mL of 400 nM SF	10.5 mL	100 nM	
0.46 mL of 100 nM SF	13.54 mL	3.3 nM	<i>Corners Test</i>
4.24 mL of 3.3 nM SF	9.76 mL	1 nM	<i>Sensitivity/Linearity Tests</i>

Fluorescence Polarization (FP) Test Solution

As described in **Fluorescence Polarization (FP) Tests** on page 103, the recommended test solutions are available from Invitrogen Corporation or from BioTek. They do not require additional preparation.

Time-Resolved Fluorescence (TRF) Test Solution

As described in **Time-Resolved Fluorescence (TRF) Tests** on page 104, the recommended test solutions are available from Invitrogen Corporation or from BioTek.

- Shake the FluoSpheres container vigorously for 30 seconds prior to pipetting. Alternatively, sonicate or vortex the container.
- Mix 10 μL of FluoSpheres with 10 mL of deionized water, in a 15 mL conical-bottom, polypropylene sample tube. This yields a 20 nM equivalent suspension.
- Shake the vial vigorously for 30 seconds prior to pipetting. Alternatively, sonicate or vortex the container.
- Mix 10 μL of 20 nM suspension with 10 mL of deionized water, in a 15 mL conical-bottom, polypropylene sample tube. This yields a 20 pM equivalent suspension.
- Refrigerate any unused portions of the FluoSpheres. The temperature must be between $+2^{\circ}\text{C}$ to $+6^{\circ}\text{C}$.

The prepared TRF plate can be kept for a maximum of seven days, if covered and stored in the dark between $+2^{\circ}\text{C}$ to $+6^{\circ}\text{C}$.

Allow the plate to sit at room temperature for approximately 15 minutes prior to use.

Fluorescence Test Procedure

Top Optics

1. If you have not already done so, create the Gen5 protocols as described starting on page 116.
2. If you have not already done so, prepare the solutions for the tests you plan to perform. See **Fluorescence Test Solutions** on page 104.

Refer to the pipette maps starting on page 109 for the remaining steps.

3. Perform the Corners/Sensitivity/Linearity tests using the Top optics of the filter-based fluorescence system:
 - Pipette the solutions for the Corners and Sensitivity/Linearity Tests into a clean, new 96-well solid black plate.
 - Create an experiment based on the **Synergy Neo2_FI_T_SF.prt** protocol.
 - Perform the steps under **Determine the Optimal Read Height** to determine and set the optimal read height for the filter cube/fluid height combinations.

- Read the plate, and then save the experiment.
 - If you intend to run this test in the future with the same filter cube, select **File > Save Protocol As**, and save the updated protocol with the adjusted read height. You must, however, rerun these steps if you change the filter cube used with this test.
4. Perform the Corners/Sensitivity/Linearity tests using the Top optics of the monochromator-based fluorescence system:
 - Create an experiment based on the **Synergy Neo2_M_FI_T_SF.prt** protocol. Read the plate used in step 3, and then save the experiment.
 5. To test Fluorescence Polarization capability:
 - Pipette the solutions for the “FP” test into a new 96-well solid black plate.
 - Create an experiment based on the **Synergy Neo 2_Dual PMT_FP.prt** or **Synergy Neo 2_Single PMT_FP.prt** protocol. Read the plate and then save the experiment.
 6. To test the Time-Resolved Fluorescence - Flash capability:
 - Pipette the solutions for the “TRF” test into a new 96-well solid white plate.
 - Create an experiment based on the **Synergy Neo2_TRF.prt** protocol. Read the plate and then save the experiment.
 7. To test the Time-Resolved Fluorescence - Laser capability:
 - Pipette the solutions for the “TRF” test into a new 96-well solid white plate.
 - Create an experiment based on the **Synergy Neo2_TRF_N2.prt** protocol. Read the plate and then save the experiment.
 8. Calculate and evaluate results as described under **Results Analysis**, starting on page 112.

Bottom Optics

1. If you have not already done so, create the Gen5 protocols as described starting on page 116.
2. If you have not already done so, prepare the solutions for the tests you plan to perform. See "**Fluorescence Test Solutions**" on page 104 on page 104.

Refer to the pipette maps starting on page 109 for the remaining steps.

3. If applicable, perform the Corners/Sensitivity/Linearity tests using the Bottom optics of the filter-based fluorescence system:
 - Pipette the solutions for the Corners and Sensitivity/Linearity Tests into a clean, 96-well glass-bottom plate.
 - Create an experiment based on the **Synergy Neo2_FI_B_SF.prt** protocol.
 - Perform the steps under Determine the Optimal Read Height to determine and set the optimal read height for the filter cube/fluid height combinations.
 - Read the plate and then save the experiment.
 - If you intend to run this test in the future with the same filter cube, select **File > Save Protocol As**, and save the updated protocol with the adjusted read height. You must, however, rerun these steps if you change the filter cube used with this test.
4. Perform the Corners/Sensitivity/Linearity tests for the monochromator-based fluorescence system:
 - If you skipped step 3, pipette the solutions for the Corners and Sensitivity/Linearity Tests into a clean, new 96-well glass-bottom plate.
 - Create an experiment based on the **Synergy Neo2_M_FI_B_SF.prt** protocol. Read the plate from step 3 (or the newly created plate described above) and then save the experiment.
5. Calculate and evaluate results as described under **Results Analysis**, starting on page 112.

Determine the Optimal Read Height

1. Select **Protocol > Procedure**, and open the "Sensitivity Read" step.
2. Click **Options**, clear **Automatic Gain Adjustment**, and click **OK**.
3. Click **Show advanced options** to display the Read Height information.
4. Click **Auto-Adjust**. Set the Test Well to **D7**, and click **Start Calibration**. When prompted, place the plate on the carrier and click **OK**.
5. When calibration is complete, a graph appears. Click **Select** to use the Optimal Height.
6. Click **Options** again, and select **Automatic Gain Adjustment**.
7. Record the optimal read height, then click **OK** to close the read step and return to the Procedure dialog.

8. Open the "Sensitivity Read Buffer" step. Set the read height to match the recorded value from step 7 above, then return to the Procedure dialog.
9. Open the "Corners Read" step. Set the read height to match the recorded value from step 7 above, then return to the Procedure dialog.
10. Open the "Linearity Read" step. Repeat steps 2–6 above; set the Test Well to **C1** in step 4.
11. Return to the Procedure dialog, then click **OK** to close the Procedure and return to the experiment.

Pipette Maps

Seal the plates with foil or store them in black polyethylene bags until use. When using a clear-bottom plate, if the base of the plate is touched, clean the entire base with alcohol (95% ethanol) and then wipe with a lint-free cloth. Before placing the plate in the instrument, blow the bottom of the plate with an aerosol duster.

Perform these steps carefully, and refer to the grid on the next page.

For the **Corners** test:

- Pipette 200 μL of the **3.3 nM SF** solution into wells A1–A3, A10–A12, H1–H3, and H10–H12.
- *If using a Hellma plate:* Pipette 200 μL of buffer into the wells surrounding the 3.3 nM wells ("CBUF" in the grid).

For the **Sensitivity** test:

- Pipette 200 μL of the **1 nM SF** solution into well D7.
- Pipette 200 μL of the buffer solution into wells C9, D9, and E9.

For the **Linearity** test (wells C1–F5):

- Use a multichannel pipette with just four tips installed.
- Pipette 150 μL of buffer solution into wells C2–F5. Discard the tips.
- Pipette 150 μL of the **1 nM SF** solution into wells C1–F1.
- Pipette 150 μL of the 1 nM SF solution into wells C2–F2. Mix the wells using the pipette. Do not discard the tips.
- Aspirate 150 μL from wells C2–F2, and dispense into wells C3–F3. Mix the wells using the pipette. Do not discard the tips.

- Aspirate 150 μ L from wells C3–F3, and dispense into wells C4–F4. Mix the wells using the pipette. Do not discard the tips.
- Aspirate 150 μ L from wells C4–F4, and dispense into wells C5–F5. Mix the wells using the pipette. Do not discard the tips.
- Aspirate 150 μ L from wells C5–F5, and discard the tips.

	1	2	3	4	5	6	7	8	9	10	11	12
A	3300pM_200	3300pM_200	3300pM_200	CBUF					CBUF	3300pM_200	3300pM_200	3300pM_200
B	CBUF	CBUF	CBUF	CBUF					CBUF	CBUF	CBUF	CBUF
C	1000pM_150	500pM_150	250pM_150	125pM_150	62_5pM_150				BUF_200			
D	1000pM_150	500pM_150	250pM_150	125pM_150	62_5pM_150		1000pM_200		BUF_200			
E	1000pM_150	500pM_150	250pM_150	125pM_150	62_5pM_150				BUF_200			
F	1000pM_150	500pM_150	250pM_150	125pM_150	62_5pM_150							
G	CBUF	CBUF	CBUF	CBUF					CBUF	CBUF	CBUF	CBUF
H	3300pM_200	3300pM_200	3300pM_200	CBUF					CBUF	3300pM_200	3300pM_200	3300pM_200

Fluorescence Polarization (FP) Test

- Pipette 200 μ L of the (green) polarization buffer (BUF) into wells A6–H6.
- Pipette 200 μ L of the green high polarization reference (HPR) into wells A7–B7.
- Pipette 200 μ L of the green low polarization reference (LPR) into wells A8–H9.

	1	2	3	4	5	6	7	8	9	10	11	12
A						BUF_FP	HPR	LPR	LPR			
B						BUF_FP	HPR	LPR	LPR			
C						BUF_FP		LPR	LPR			
D						BUF_FP		LPR	LPR			
E						BUF_FP		LPR	LPR			
F						BUF_FP		LPR	LPR			
G						BUF_FP		LPR	LPR			
H						BUF_FP		LPR	LPR			

Time-Resolved Fluorescence (TRF) Test

- Pipette 200 μ L of deionized water into wells A6–C6.
- If you have not already done so, shake the vial of 20 pM europium suspension vigorously for 30 seconds prior to pipetting. Alternatively, sonicate or vortex the vial.
- Pipette 200 μ L of the 20 pM europium suspension (Eu) into well A8.

	1	2	3	4	5	6	7	8	9	10	11	12
A						BUF		Eu				
B						BUF						
C						BUF						
D												
E												
F												
G												
H												

Results Analysis

Corners Test

1. Calculate the Mean of the 12 wells containing the 3.3 nM SF test solution (A1–A3, A10–A12, H1–H3, and H10–H12).
2. Calculate the Standard Deviation for the same 12 wells.
3. Calculate the %CV: $(\text{Standard Deviation} / \text{Mean}) * 100$

The %CV must be < **3.0** to pass.

Sensitivity Test

1. Calculate the Mean and Standard Deviation of the 16 reads for each of the buffer wells (C9, D9, E9).
2. Among the three buffer wells, find the Median Standard Deviation and corresponding Mean.
3. Calculate the Mean for the 16 reads of the SF Concentration well (D7).
4. Calculate the Signal-to-Noise Ratio (SNR) using the Mean SF Concentration, Buffer Media STD with its corresponding Buffer Mean:
 $(\text{SF Mean} - \text{Buffer Mean}) / 3 * \text{Buffer STD}$
5. Calculate the Detection Limit, in pM, using the known concentration value of SF and the Calculated SNR: $1000 / \text{SNR}$

Filter-Based Fluorescence System		
Optics	Filter Cube	To pass, the Detection Limit must be:
Top	Cube 107 (or equivalent): 485/20, 528/20, 510 nm dichroic mirror	< 5 pM
Bottom	Cube 107 (or equivalent): 485/20, 528/20, 510 nm dichroic mirror	< 5 pM

Monochromator-Based Fluorescence System		
Optics	Wavelength	To pass, the Detection Limit must be:
Top/Bottom	EX 485 nm, EM 528 nm	≤ 20 pM

Linearity Test

1. Calculate the Mean of the four wells for each concentration in columns 1–5.
2. Perform linear regression using these values as inputs:

Filter- and Monochromator-Based Fluorescence System	
X	Y
1000	Mean of the 1000 pM wells
500	Mean of the 500 pM wells
250	Mean of the 250 pM wells
125	Mean of the 125 pM wells
62.5	Mean of the 62.5 pM wells

3. Calculate the R-Square value; it must be ≥ 0.9500 to pass.

Fluorescence Polarization (FP) Test

1. Using the raw data from the Parallel read:
 - Calculate the Mean Blank (wells A6–H6).
 - Calculate the Signal for each HPR well: Subtract the Mean Blank from its measurement value.
 - Calculate the Signal for each LPR well: Subtract the Mean Blank from its measurement value.
2. Using the raw data from the Perpendicular read:

- Calculate the Mean Blank (wells A6–H6)
 - Calculate the Signal for each HPR well: Subtract the Mean Blank from its measurement value.
 - Calculate the Signal for each LPR well: Subtract the Mean Blank from its measurement value.
3. Calculate the G-Factor for each LPR well:

$$\frac{\text{Parallel LPR Signal} * (1-0.02)}{\text{Perpendicular LPR Signal} * (1+0.02)}$$
 4. Calculate the Mean G-Factor.
 5. Calculate the Polarization value in mP for each HPR well (“PHPR”):

$$\frac{\text{Parallel HPR Signal} - \text{Mean G-Factor} * \text{Perpendicular HPR Signal} * 1000}{\text{Parallel HPR Signal} + \text{Mean G-Factor} * \text{Perpendicular HPR Signal}}$$

6. Calculate the Mean PHPR, in mP.

Filter Cubes	To pass, the Mean PHPR must be:
Dual PMT—Cube 61 (or equivalent); Single PMT—Cube 108 (or equivalent): 485/20, 528/20, 510 nm, dichroic mirror	> 340 mP

7. Calculate the Polarization value in mP for each LPR well (“PLPR”):

$$\frac{\text{Parallel LPR Signal} - \text{Mean G-Factor} * \text{Perpendicular LPR Signal} * 1000}{\text{Parallel LPR Signal} + \text{Mean G-Factor} * \text{Perpendicular LPR Signal}}$$

8. Calculate the Standard Deviation of the “PLPR,” in mP.

Filter Cubes	To pass, the Standard Deviation of the PLPR:
Dual PMT—Cube 61 (or equivalent); Single PMT—Cube 108 (or equivalent): 485/20, 528/20, 510 nm dichroic mirror	< 5

Time-Resolved Fluorescence (TRF) Test

1. Calculate the Mean and Standard Deviation of the 16 reads for each of the buffer wells (A6, B6, C6).
2. Among the three buffer wells, find the Media Standard Deviation and corresponding Mean.
3. Calculate the Mean for the 16 reads of the Eu Concentration well (A8).
4. Calculate the Signal-to-Noise Ratio (SNR) using the Mean Eu Concentration and Buffer Median STD with its corresponding Buffer Mean:

$$(\text{Eu Mean} - \text{Buffer Mean}) / (3 * \text{Buffer STD})$$
5. Calculate the Detection Limit, in fM:

$$20000 / (\text{Mean Eu} - \text{Mean DI water}) / (3 * \text{Standard Deviation DI water})$$

Light Source	Filter Cube	To pass, the Detection Limit must be:
Xenon flash	Cube 112 (or equivalent): 360/40, 620/40, 400 nm dichroic mirror	≤ 250 fM
TRF laser	Cube 126: 337/10, 620/40	≤ 40 fM

Troubleshooting Fluorescence Liquid Tests

If any tests fail, please try the following suggestions. If the test(s) continue to fail, print the results and contact Technical Support.

- Are the solutions fresh? Discard the plate and any opened, unused test solutions after seven days.
- Are the excitation/emission filters clean?
- Are you using the proper filter cube?
- If the Corners Test continues to fail, the hardware may be misaligned. Contact Technical Support.
- Are you using new/clean plates? If the base of a clear-bottom plate is touched, clean the entire base with alcohol (95% ethanol) and then wipe with a lint-free cloth. Before placing the plate in the instrument, blow the bottom of the plate with an aerosol

duster. If the test fails again, the optical probe(s) may need to be cleaned. Contact Technical Support.

- Review the pipetting instructions to verify the plate was correctly prepared.
- Does the Plate Type setting in the Gen5 protocol match the plate you used?
- For models with a dispenser, spilled fluid inside the reader may be fluorescing, which can corrupt your test results. If you suspect this is a problem, contact Technical Support.
- When testing Fluorescence Polarization capability using a solid black plastic microplate, if the standard deviation for the buffer wells is too high, try moving the buffer wells to another column. With some black plastic plates, the wells in the center of the plate may be slightly distorted due to the plate molding process, and this can affect the standard deviation.
- The Read steps in the protocols use the Gen5 Automatic Gain Adjustment feature to determine optimum sensitivity values for the plate. If an Auto Gain Result value is outside the range of 30–200, this may indicate a problem.

If the value is less than 30:

- The stock solution/dilution concentrations may be too high. Try creating fresh solutions/dilutions, and rerun the test using a new, clean plate.
- If all of the tests are passing but the Gain value is low, a PMT in your reader may just be very sensitive. Contact Technical Support to confirm that this may be the case.

If the value is greater than 200:

- The stock solution/dilution concentrations may be too low. Try creating fresh solutions/dilutions, and rerun the test using a new, clean plate.
- For injector models, spilled fluid inside the reader may be fluorescing, which can corrupt your test results. If you suspect this is a problem, contact Technical Support.
- The PMTs or optical path(s) may be deteriorating, or the optics or other hardware may be misaligned. Contact Technical Support.

Gen5 Protocol Reading Parameters

The information in the following tables represents the recommended reading parameters. It is possible that your tests will require modifications to some of these parameters, such as the Plate Type (see **Troubleshooting Tips** on page 115).

The Plate Type setting in each Gen5 protocol must match the plate you are actually using.

Synergy Neo 2_FI_T_SF.prt/Synergy Neo 2_FI_B_SF.prt

Parameter	Setting
Detection Method:	Fluorescence intensity
Read Type:	Endpoint
Plate Type:	Top Read: Costar 96 black opaque (#3915) Bottom Read: Greiner Sensoplate
Read Step 1:	
Kinetic:	Run Time: 0:00:45 Interval: 0:00:03 Reads: 16
Step Label:	Sensitivity Read
Read Well:	D7
Filter Sets:	Single PMT Filter set 1: 485/528 Optics position: Top/Bottom Gain: (top and bottom optics) Auto, Scale to High Wells, D7, 50000
Read Speed:	Normal Delay after plate movement: 350 msec Measurements per data point: 50
Read Height:	5.75 mm (top read); 9.25 mm (bottom read)
Dynamic Range:	Standard
Light source:	Xenon Flash
Lamp energy:	Low (faster)
Read Step 2:	
Kinetic:	Run Time: 0:01:35 Interval: 0:00:06 Reads: 16

Parameter	Setting
Step Label:	Sensitivity Read Buffer
Read Wells:	C9..E9
Filter Sets:	Single PMT Filter set 1: 485/528 Optics position: Top/Bottom Gain: (top and bottom optics) Auto, Use first filter set gain from FIRST Read Step
Read Speed:	Normal Delay after plate movement: 350 msec Measurements per data point: 50
Read Height:	5.75 mm (top read); 9.25 mm (bottom read)
Dynamic Range:	Standard
Light source:	Xenon Flash
Lamp energy:	Low (faster)
Read Step 3:	
Step Label:	Corners Read
Read Wells:	A1–A3, A10–A12, H1–H3, H10–H12
Filter Sets:	Single PMT Filter set 1: 485/528 Optics position: Top/Bottom Gain: (top and bottom optics) Auto, Scale to High Wells, A3, 50000
Read Speed:	Normal Delay after plate movement: 350 msec Measurements per data point: 50
Read Height:	5.75 mm (top read); 9.25 mm (bottom read)
Dynamic Range:	Standard
Light source:	Xenon Flash
Lamp energy:	Low (faster)
Read Step 4:	

Parameter	Setting
Step Label:	Linearity Read
Read Wells:	C1–F5
Filter Sets:	Single PMT Filter set 1: 485/528 Optics position: Top/Bottom Gain: (top and bottom optics) Auto, Scale to High Wells, C1, 50000
Read Speed:	Normal Delay after plate movement: 350 msec Measurements per data point: 50
Read Height:	5.75 mm (top read); 9.25 mm (bottom read)
Dynamic Range:	Standard
Light source:	Xenon Flash
Lamp energy:	Low (faster)

Synergy Neo 2_Single PMT_FP.prt and Synergy Neo 2_Dual PMT_FP.prt

This procedure contains one Read step using filters with Fluorescence Polarization enabled, inside a Plate Mode block.

Parameter	Default Setting
Detection Method:	Fluorescence polarization
Read Type:	Endpoint
Plate Type:	Costar 96 black opaque (#3915)
Synchronized Mode:	Plate Mode with Timing Control
Read Wells:	A6–H9

Parameter	Default Setting
Filter Sets:	Single/Dual PMT FP 485/528 Optics position: Top Gain: Single PMT—Auto, Scale to well: A9, Scale value: 20000 Gain: Dual PMT—Side PMT1/Top PMT2: Auto, Scale to well: A9, Parallel scale value: 20000, Requested polarization: 20
Read Speed:	Normal
Delay After Plate Movement:	0 msec
Measurements Per Data Point:	50
Read Height:	6.50 mm
Dynamic Range:	Standard
Light source:	Xenon Flash
Lamp energy:	Low (faster)

Synergy Neo 2_TRF.prt

Parameter	Default Setting
Detection Method:	Time-resolved fluorescence -flash
Read Type:	Endpoint
Plate Type:	Costar 96-well white opaque
Delay Step:	3 minutes
Read Step 1:	
Kinetic:	Run Time: 0:00:15 Intervale: 0:00:01 Reads: 16
Step Label:	Sensitivity Read
Read Well:	A8

Parameter	Default Setting
Filter Sets:	Single PMT 360/620 Optics position: Top Gain: Auto, Scale to High Wells, A8, 30000
Read Speed:	Normal
Delay after plate movement:	100 msec
Measurements per data point:	50
Read Height:	4.50 mm
Light source:	Xenon Flash
Lamp energy:	Low (faster)
Dynamic Range:	Standard
Delay:	300 μ sec
Data collection time:	1000 μ sec
Read Step 2:	
Kinetic:	Run Time: 0:00:45 Interval: 0:00:03 Reads: 16
Step Label:	Sensitivity Read Buffer
Read Wells:	A6–C6
Filter Sets:	Single PMT 360/620 Optics position: Top Gain: Auto, Use first filter set gain from FIRST Read Step
Read Speed:	Normal
Delay after plate movement:	100 msec
Measurements per data point:	50
Read Height:	4.50 mm
Light source:	Xenon Flash

Parameter	Default Setting
Lamp energy:	Low (faster)
Dynamic Range:	Standard
Delay:	300 μ sec
Data collection time:	1000 μ sec

Synergy Neo 2_M_FI_T_SF.prt and Synergy Neo 2_M_FI_B_SF.prt

Parameter	Default Setting
Detection Method:	Fluorescence
Read Type:	Endpoint
Plate Type:	Top: Costar 96 black opaque Bottom: Greiner SensoPlate
Read Step 1:	
Kinetic:	Run Time: 0:00:45 Interval: 0:00:03 Reads: 16
Step Label:	Sensitivity Read
Read Well:	D7
Wavelength:	Excitation: 485/14, Emission: 528/14
Optics Position:	Top/Bottom
Gain	Top: Auto, Scale to High Wells, D7, 50000
Read Speed:	Normal
Delay after plate movement:	100 msec
Measurements per data point:	50
Lamp Energy:	Low (faster)
Dynamic Range:	Standard
Read Height (for top optics):	5.00 mm
Read Step 2:	

Parameter	Default Setting
Kinetic:	Run Time: 0:01:35 Interval: 0:00:06 Reads: 16
Step Label:	Sensitivity Read Buffer
Read Wells:	C9–E9
Wavelengths:	Excitation: 485/14, Emission: 528/14
Optics Position:	Top/Bottom
Gain:	Top: Auto, Use first filter set gain from FIRST Read Step
Read Speed:	Normal
Delay after plate movement:	100 msec
Measurements per data point:	50
Lamp Energy:	Low (faster)
Dynamic Range:	Standard
Read Height (for top optics):	5.00 mm
Read Step 3:	
Step Label:	Corners Read
Read Wells:	A1–A3, A10–A12, H1–H3, H10–H12
Wavelengths:	Excitation: 485/14, Emission: 528/14
Optics Position:	Top/Bottom
Gain:	Top: Auto, Scale to High Wells, A3, 50000
Read Speed:	Normal
Delay after plate movement:	100 msec
Measurements per data point:	50
Lamp Energy:	Low (faster)
Dynamic Range:	Standard

Parameter	Default Setting
Read Height (for top optics):	5.00 mm
Read Step 4:	
Step Label:	Linearity Read
Read Wells:	C1–F5
Wavelengths:	Excitation: 485/14, Emission: 528/14
Optics Position:	Top
Gain:	Top: Auto, Scale to High Wells, C1, 50000
Read Speed:	Normal
Delay after plate movement:	100 msec
Measurements per data point:	50
Lamp Energy:	Low (faster)
Dynamic Range:	Standard
Read Height (for top optics):	5.00 mm

Synergy Neo 2_TRF-N2.prt (requires Gen5 v. 3.06 or higher)

To qualify Synergy Neo 2 “T” models. Dedicated filter cube required: Cube #126 (337/620; LUM) PN 1035126

Parameter	Default Setting
Detection Method:	Time-resolved fluorescence -laser
Read Type:	Endpoint
Plate Type:	Costar 96-well white opaque
Delay Step:	3 minutes
Read Step 1:	
Kinetic:	Run Time: 0:00:30 Intervale: 0:00:02 Reads: 16
Step Label:	Sensitivity Read

Parameter	Default Setting
Read Well:	A8
Filter Sets:	Single PMT 337/620 Optics position: Top Gain: Auto, Scale to High Wells, A8, 50000
Read Speed:	Normal
Delay after plate movement:	100 msec
Measurements per data point:	50
Read Height:	4.50 mm
Light source:	TRF Laser
Dynamic Range:	Standard
Delay:	250 μ sec
Data collection time:	1000 μ sec
Read Step 2:	
Kinetic:	Run Time: 0:01:15 Interval: 0:00:05 Reads: 16
Step Label:	Sensitivity Read Buffer
Read Wells:	A6–C6
Filter Sets:	Single PMT 360/620 Optics position: Top Gain: Auto, Use first filter set gain from FIRST Read Step
Read Speed:	Normal
Delay after plate movement:	100 msec
Measurements per data point:	50
Read Height:	4.50 mm
Light source:	TRF Laser

Parameter	Default Setting
Dynamic Range:	Standard
Delay:	250 μ sec
Data collection time:	1000 μ sec

Alternate/Supplemental Tests Using Methylumbelliferone (MUB)

As an alternative to using Sodium Fluorescein, Methylumbelliferone (“MUB”) can be used to test the top and bottom optics for filter-based fluorescence systems and the top optics for monochromator-based fluorescence systems.

Required Materials

Microplates should be perfectly clean and free from dust or bottom scratches. Use new microplates from sealed packages.

Manufacturer part numbers are subject to change over time.

BioTek offers a liquid test kit (PN 7160012) containing the microplates and solutions used in the MUB fluorescence liquid tests for the Top optics. If also testing the Bottom optics (filter-based system), you will need an additional microplate (see below).

- Methylumbelliferone (“MUB”) (10-mg vial, BioTek PN 98156)
- Carbonate-Bicarbonate buffer (“CBB”) capsules (BioTek PN 98158)
- 100% methanol (BioTek PN 98161)
- **Top optics:** A new, clean 96-well solid black plate microplate, such as Corning Costar #3915 or equivalent. The same plate is used to test both filter- and monochromator-based systems.
- **Bottom optics:** A new, clean 96-well glass-bottom Greiner SensoPlate (Mfr. #655892), a clean Hellma Quartz 96-well titration plate (Mfr. #730.009.QG), or equivalent
- Filter cube 107 or equivalent (360/460)
- Deionized or distilled water
- Various beakers, graduated cylinders, and pipettes
- 95% ethanol (for cleaning clear-bottom plates)

- Aluminum foil
- (Optional, but recommended) 0.45-micron filter
- (Optional) Black polyethylene bag(s) to temporarily store plate(s)
- Gen5 protocols described on page 132:
 - **Synergy Neo 2_FI_T_MUB.prt** tests the top filter-based fluorescence system
 - **Synergy Neo 2_FI_B_MUB.prt** tests the bottom filter-based fluorescence system
 - **Synergy Neo 2_M_FI_T_MUB.prt** tests the top optics of the monochromator-based fluorescence system

Test Solutions

Filter solutions to remove particulates that could cause erroneous readings. Do not allow dust to settle on the surface of the solution; use microplate covers or seals when not reading the plate.

Wrap the vial containing the MUB stock solution in foil to prevent exposure to light.

Discard any open, unused solutions after seven days.

1. Prepare the buffer (CBB) solution:
 - (Optional, but recommended) Using a 0.45-micron filter, filter 200 mL of deionized or distilled water.
 - Open and dissolve the contents of two CBB capsules (do not dissolve the outer gelatin capsule) into 200 mL of the water.
 - Stir the solution (preferably using a stir table) until the CBB is completely dissolved.
2. Prepare the MUB stock solution:
 - Add 1 mL of 100% methanol to the 10 mg vial of MUB.
 - Make sure all of the dye has completely dissolved and is well mixed. This yields a 10 mg/mL stock solution.

- Wrap the solution in aluminum foil to prevent exposure to light.
3. Prepare the dilutions. Label each with “MUB” and the concentration.

Mix This MUB Solution:	With:	To Make:
0.5 mL of 10 mg/mL stock solution	4.5 mL of 100% methanol	1 mg/mL
0.88 mL of 1 mg/mL solution	4.12 mL of CBB	176 µg/mL
0.1 mL of 176 µg /mL solution	9.9 mL of CBB	1.76 µg /mL
0.5 mL of 1.76 µg /mL solution	4.5 mL of CBB	176 ng/mL
1 mL of 176 ng/mL solution	9 mL of CBB	17.6 ng/mL (100 nM)

Procedure

1. If you have not already done so, create the Gen5 protocols as described on page 132.
2. If you have not already done so, prepare the test solutions. See page 127.
3. Perform the Sensitivity and Linearity tests using the Top optics of the filter-based fluorescence system:
 - Refer to the pipette map in the next section and pipette the solutions into a clean, 96-well plate.
 - Create an experiment based on the **Synergy Neo 2_FI_T_MUB.prt** protocol.
 - Perform the steps under **Determine the Optimal Read Height** to determine and set the optimal read height for the filter cube/fluid height combination.
 - Read the plate, and then save the experiment.
 - If you intend to run this test in the future **with the same filter cube**, select **File > Save Protocol As**, and save the updated protocol with the adjusted read height. You must, however, rerun these steps if you change the filter cube used with this test.
4. Perform the Sensitivity and Linearity tests for the monochromator-based fluorescence system:
 - Create an experiment based on the **Synergy Neo 2_M_FI_T_MUB.prt** protocol.
 - Read the plate from step 3, and then save the experiment.

5. If applicable, perform the Sensitivity and Linearity tests using the Bottom optics of the filter-based fluorescence system:
 - Refer to the pipette map in the next section, and pipette the solutions into a clean, 96-well clear-bottom plate.
 - Create an experiment based on the **Synergy Neo 2_FI_B_MUB.prt** protocol.
 - Perform the steps under Determine the Optimal Read Height to determine and set the optimal read height for the filter cube/fluid height combination.
 - Read the plate, and then save the experiment.
 - (Optional) Save the updated protocol with the adjusted read height.
6. Calculate and evaluate the results as described under Results Analysis, starting on page 131.

Determine the Optimal Read Height

1. Select **Protocol > Procedure**, and open the "Sensitivity Read" step.
2. Click **Options**, clear **Automatic Gain Adjustment**, and click **OK**.
3. Click **Show advanced options** to display the Read Height information.
4. Click **Auto-Adjust**. Set the Test Well to **D7**, and click **Start Calibration**. When prompted, place the plate on the carrier and click **OK**.
5. When calibration is complete, a graph appears. Click **Select** to use the Optimal Height.
6. Click **Options** again, and select **Automatic Gain Adjustment**.
7. Record the optimal read height, then click **OK** to close the read step and return to the Procedure dialog.
8. Open the "Sensitivity Read Buffer" step. Set the read height to match the recorded value from step 7 above, then return to the Procedure dialog.
9. Open the "Linearity Read" step. Repeat steps 2–6 above; set the Test Well to **C1** in step 4.
10. Return to the Procedure dialog, then click **OK** to close the Procedure and return to the experiment.

Pipette Map

Seal the plate with foil or store it in a black polyethylene bag until use. When using a clear-bottom plate, if the base of the plate is touched, clean the entire base with alcohol (95% ethanol) and then wipe with a lint-free cloth. Before placing a plate in the instrument, blow the bottom of the plate with an aerosol duster.


For the Sensitivity test:

- Pipette 200 μL of 100 nM MUB solution into well D7.
- Pipette 200 μL of buffer into wells C9, D9, and E9.

For the Linearity test:

Using a multi-channel pipette with just four tips installed to process rows C–F:

- Pipette 150 μL of buffer into columns 2–5 (**not column 1**). Discard the tips.
- Pipette 150 μL of the 17.6 ng/mL (100 nM) solution into column 1. Discard the tips.
- Pipette 150 μL of the 17.6 ng/mL (100 nM) solution into column 2. Do not discard the tips.
- Aspirate 150 μL from column 2 and dispense it into column 3. Mix the wells using the pipette. Do not discard the tips.
- Aspirate 150 μL from column 3 and dispense it into column 4. Mix the wells using the pipette. Do not discard the tips.
- Aspirate 150 μL from column 4 and dispense it into column 5. Mix the wells using the pipette. Do not discard the tips.
- Aspirate 150 μL from column 5. Discard the tips.

	1	2	3	4	5	6	7	8	9	10	11	12
A												
B												
C	100nM_150	50nM_150	25nM_150	12_5nM_150	6_25nM_150				BUF_200			
D	100nM_150	50nM_150	25nM_150	12_5nM_150	6_25nM_150		100nM_200		BUF_200			
E	100nM_150	50nM_150	25nM_150	12_5nM_150	6_25nM_150				BUF_200			
F	100nM_150	50nM_150	25nM_150	12_5nM_150	6_25nM_150							
G												
H												

Results Analysis

Sensitivity Test

1. Calculate the Mean and Standard Deviation of the 16 reads for each of the buffer wells (C9, D9, E9).
2. Among the three buffer wells, find the Median Standard Deviation and corresponding Mean.
3. Calculate the Mean for the 16 reads of the MUB concentration well (D7).
4. Calculate the Signal-to-Noise Ratio (SNR) using the Mean MUB Concentration and Buffer Media STD with its corresponding Buffer Mean:

$$(\text{Mean MUB} - \text{Buffer Mean}) / (3 * \text{Buffer STD})$$
5. Calculate the Detection Limit, in ng/mL, using the known concentration value of MUB and the calculated SNR: $17.6/\text{SNR}$

Filter-Based Fluorescence System:

Filter Cube	Optic Position	To pass, the Detection Limit must be:
Cube 107 (or equivalent)—360/40, 460/40, 400 nm dichroic mirror	Top/Bottom	≤ 0.16 ng/mL (0.91 nM)

Monochromator-Based Fluorescence System:

Wavelengths	Optic Position	To pass, the Detection Limit must be:
Excitation—360 nm Emission—460 nm	Top	≤ 0.16 ng/mL (0.91 nM)

Linearity Test

1. Calculate the Mean of the four wells for each concentration in columns 1–5.
2. Perform linear regression using these values as inputs:

x	y
100	Mean of the 100 nM wells
50	Mean of the 50 nM wells
25	Mean of the 25 nM wells
12.5	Mean of the 12.5 nM wells
6.25	Mean of the 6.25 nM wells

3. Calculate the R-Square value; it must be ≥ 0.9500 to pass.

Gen5 Protocol Reading Parameters

The information in the following table represents the recommended reading parameters. It is possible that your test will require modifications to some of these parameters, such as the Plate Type (see Troubleshooting on page 115).

The Plate Type setting in the Gen5 protocol should match the plate you are actually using.

Neo 2_FI_T_MUB.prt and Neo 2_FI_B_MUB.prt

Parameter	Default Setting
Plate Type:	Top Optics: Costar 96-well black opaque (#3915) Bottom Optics: Greiner Sensoplate
Detection Method:	Fluorescence
Read Type:	Endpoint
Read Step 1:	
Kinetic	Run Time: 0:00:45 Interval: 0:00:03 Reads: 16
Step Label:	Sensitivity Read
Read Well:	D7
Filter Sets:	Single PMT 360/460 Optics position: Top/Bottom Gain: Auto, Scale to High Wells, D7, 80000
Read Speed:	Normal Delay after plate movement: 350 msec Measurements per data point: 50
Read Height:	4.50 mm (top read); 8.75 mm (bottom read)
Dynamic Range:	Standard
Light Source:	Xenon Flash
Lamp Energy:	Low (faster)
Read Step 2:	
Kinetic:	Run Time: 0:01:35 Interval: 0:00:06 Reads: 16
Step Label:	Sensitivity Read Buffer
Read Well:	C9, D9, E9

Parameter	Default Setting
Filter Sets:	Single PMT 360/460 Optic position: Top/Bottom Gain: Auto, Use first filter set gain from FIRST Read Step
Read Speed:	Normal Delay after plate movement: 350 msec Measurements per data point: 50
Read Height:	4.50 mm (top read); 8.75 mm (bottom read)
Dynamic Range:	Standard
Light Source:	Xenon Flash
Lamp Energy:	Low (faster)
Read Step 3:	
Step Label:	Linearity Read
Read Well:	C1–F5
Filter Sets:	Single PMT 360/460 Optics position: Top/Bottom Gain: Auto, Scale to High Wells, C1, 80000
Read Speed:	Normal Delay after plate movement: 350 msec Measurements per data point: 50
Read Height:	4.50 mm (top read); 8.75 mm (bottom read)
Dynamic Range:	Standard
Light Source:	Xenon Flash
Lamp Energy:	Low (faster)

Neo 2_M_FI_T_MUB.prt

Parameter	Default Setting
Plate Type:	Costar 96-well black opaque (#3915)

Parameter	Default Setting
Detection Method:	Fluorescence
Read Type:	Endpoint
Read Step 1:	
Kinetic:	Run Time: 0:00:45 Interval: 0:00:03 Reads: 16
Step Label:	Sensitivity Read
Read Wells:	D7
Wavelengths:	Excitation: 360/14, Emission: 460/14
Optics Position:	Top
Gain:	Auto, Scale to High Wells, D7, 80000
Read Speed:	Normal
Delay after plate movement:	100 msec
Measurements per data point:	50
Lamp Energy:	Low (faster)
Dynamic Range:	Standard
Read Height:	4.50 mm
Read Step 2:	
Kinetic:	Run Time: 0:01:35 Interval: 0:00:03 Reads: 16
Step Label:	Sensitivity Read Buffer
Read Wells:	C9, D9, E9
Wavelengths:	Excitation: 360/14, Emission: 460/14
Optics Position:	Top
Gain:	Auto, Use first filter set gain from FIRST Read Step
Read Speed:	Normal

Parameter	Default Setting
Delay after plate movement:	100 msec
Measurements per data point:	50
Lamp Energy:	Low (faster)
Dynamic Range:	Standard
Read Height:	4.50
Read Step 3:	
Step Label:	Linearity Read
Read Wells:	C1–F5
Wavelengths:	Excitation: 360/14, Emission: 460/14
Optics Position:	Top
Gain:	Auto, Scale to High Wells, C1, 80000
Read Speed:	Normal
Delay after plate movement:	100 msec
Measurements per data point:	50
Lamp Energy:	Low (faster)
Dynamic Range:	Standard
Read Height:	4.50 nm

Luminescence Test

Luminescence performance is verified using the Harta™ Luminometer Reference Microplate, which is an LED-based test plate available for purchase.

Depending on the instrument model, one or two luminescence tests should be performed. Filter-based luminescence is available in most instrument models (except NEO2MON-N). “M”-models with a monochromator also include a broadband fiber for luminescence. BioTek recommends performing both luminescence tests at initial setup. On a monthly or quarterly basis, your organization’s instrument qualification practices and the frequency of performing luminescence assays with each method should dictate the necessity of testing performance with both methods.

Required Materials

- Harta Luminometer Reference Microplate, PN 8030015
- Gen5 protocol (see page 139)
- Filter-based test: LUM Upper and Lower Top filter cubes (e.g., #3 and #112, or equivalent)

Procedure

1. Turn on the Harta reference plate using the I/O switch on the back of the plate.
2. Check the plate’s battery by pressing the test button on the back of the plate and ensuring that the test light turns on.

The test light may be difficult to see in bright light. Change your angle of view or move to a darker environment if you cannot see it.

3. Place the Harta plate adapter on the reader’s carrier, then place the test plate on top of the adapter (PN 8042263).
4. Create an experiment based on the **Synergy Neo2 LumTest_Harta.prt** or **Synergy Neo2 M-LumTest_Harta.prt** protocol and read the plate.
5. Calculate and evaluate results as described under **Results Analysis** below.

Plate Layout

Filter LUM test: Synergy Neo2 LumTest_Harta.prt

	1	2	3	4	5	6	7	8	9	10	11	12
A		REF					LED7	LED8				
B												
C												
D	Buffer	Buffer	Buffer	Buffer								
E	Buffer	Buffer	Buffer	Buffer								
F	Buffer	Buffer	Buffer	Buffer								
G	Buffer	Buffer	Buffer	Buffer								
H												

Monochromator broadband LUM: Synergy Neo2 M-LumTest_Harta.prt

	1	2	3	4	5	6	7	8	9	10	11	12
A		REF					LED7	LED8				
B												
C												
D												
E												
F	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer
G	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer	Buffer
H												

Results Analysis

Through a manual correlation process, it was found that the system requires approximately 35 photons per attomole of ATP, thus a conversion factor of 0.02884 attomole/photon was applied to determine ATP concentration from the NIST data in photons/s.

1. On the Harta plate's Calibration Certificate, locate the NIST measurement for the A2 position and convert it to attomoles:
(A2 NIST measurement * 0.02884)
2. Determine if the plate's battery is functioning properly:
 - If $A8 > (0.2 * A7)$, the battery is good. Otherwise, it requires replacement.

A replacement battery is included with each Harta plate. A new spare battery will be supplied when the plate is recertified.

3. Calculate the signal-to-noise ratio:
(A2 - Mean of the buffer cells) / (3 * Standard deviation of buffer cells)
4. Calculate the detection limit:
A2 NIST measurement in attomoles / signal-to-noise ratio
 - If the reader is equipped with the low-noise PMT, the detection limit must be:
 - filter-based luminescence: **≤ 50 amol** to pass.
 - broadband-fiber luminescence: **≤ 75 amol** to pass.
 - If the reader is equipped with the red-shifted PMT (monochromator or filter), the detection limit must be **≤ 500 amol** to pass with either method.

Gen5 Protocol Reading Parameters

The information in the following tables represents the recommended reading parameters.

Synergy Neo 2_LumTest_Harta.prt

Parameter	Default Setting
Plate Type:	8030015 Harta - w/o 8032028 adapter

Parameter	Default Setting
Delay Step:	3 minutes
Read Step 1:	
Detection Method:	Luminescence
Read Type:	Endpoint
Optics Type:	Filters
Step Label:	Reference well A2
Read Wells:	A2
Light path 1:	LUM (Single PMT)
Optics position:	Top
Gain:	200
Integration Time:	0:10.00 MM:SS.ss
Delay After Plate Movement:	0 msec
Read Height:	10.00 mm
Read Step 2:	
Detection Method:	Luminescence
Read Type:	Endpoint
Step Label:	Background
Read Wells:	D1..G4
Light path 1:	LUM (Single PMT)
Gain:	200
Integration Time:	0:10.00 MM:SS.ss
Delay After Plate Movement:	0 msec
Read Height:	10.00 mm
Read Step 3:	
Detection Method:	Luminescence

Parameter	Default Setting
Read Type:	Endpoint
Read Wells:	A7–A8
Step Label:	Battery Check
Light path 1:	LUM (Single PMT)
Optics Position:	Top
Gain:	60
Integration Time:	0:01.00 MM:SS.ss
Delay After Plate Movement:	0 msec
Read Height:	10.00 mm

Synergy Neo 2 M-LumTest_Harta.prt (requires Gen5 v. 3.06 or higher)

Parameter	Default Setting
Plate Type:	8030015 Harta - w/o 8032028 adapter
Delay Step:	3 minutes
Read Step 1:	
Detection Method:	Luminescence
Read Type:	Endpoint
Optics Type:	Luminescence fiber
Step Label:	Reference well A2
Read Wells:	A2
Light path 1:	LUM (Single PMT)
Gain:	200
Integration Time:	0:10.00 MM:SS.ss
Delay After Plate Movement:	0 msec
Read Height:	1.00 mm

Parameter	Default Setting
Read Step 2:	
Detection Method:	Luminescence
Read Type:	Endpoint
Step Label:	Background
Read Wells:	F1..G12
Light path 1:	LUM (Single PMT)
Gain:	200
Integration Time:	0:10.00 MM:SS.ss
Delay After Plate Movement:	0 msec
Read Height:	1.00 mm
Read Step 3:	
Detection Method:	Luminescence
Read Type:	Endpoint
Read Wells:	A7–A8
Step Label:	Battery Check
Light path 1:	LUM (Single PMT)
Gain:	60
Integration Time:	0:01.00 MM:SS.ss
Delay After Plate Movement:	0 msec
Read Height:	1.00 mm

Luminescence Troubleshooting

If the luminescence test fails, try the following suggestions. If it continues to fail, print the results and contact Technical Support.

- Ensure that the reading is performed through a hole in the filter cube, not through a glass filter.
- Verify that the filter cube settings in Gen5 match the physical cube.

- If the test continues to fail, the optical probe(s) may need to be cleaned. Contact Technical Support for instructions.

Dispense Module Tests

This section applies only to models with the dispenser.

BioTek has developed a set of tests that you can perform to ensure that the dispense module performs to specification initially and over time. We recommend that you perform these tests before first use (e.g., during the Initial OQ), and then every three months.

- The **Accuracy Test** is a measure of the mean volume per well for multiple dispenses. The actual weight of the dispensed fluid is compared to the expected weight and must be within a certain percentage to pass. Pass/Fail criteria depends on the per-well volume dispensed: 2.0% for 80 μL , 5.0% for 20 μL , and 20.0% for 5 μL . It is assumed that one gram is equal to one milliliter.

The test uses a single green dye test solution and a 96-well microplate (per injector) to test the three different volumes. The balance is tared with the empty plate, and then the 80 μL dispense is performed for columns 1–4. The fluid is weighed and the balance is tared again (with the plate on the balance). This process is repeated for the 20 μL and 5 μL dispenses. It is assumed that the solutions used are at room temperature. A precision balance (three-place) is used to weigh the plate.

- The **Precision Test** is a measure of the variation among volumes dispensed to multiple wells. For each volume dispensed (80 μL , 20 μL , and 5 μL) to four columns, the %CV (coefficient of variation) of 32 absorbance readings is calculated. Pass/Fail criteria depends on the per-well volume dispensed: 2.0% for 80 μL , 7.0% for 20 μL , and 10.0% for 5 μL . The plate is read in an absorbance reader at 405/750 nm for columns 1–4 and at 630/750 nm for columns 5–12.

The two tests are performed simultaneously and use the same plate.

Failures

If any tests fail, prime the fluid lines and rerun the test(s).

If the test(s) fail again, the injectors may require cleaning; see **Periodic Maintenance**.

If tests continue to fail, contact Technical Support.

Required Materials

Manufacturer part numbers are subject to change over time.

- Absorbance reader with capability of reading at 405, 630, and 750 nm. The reader must have an accuracy specification of $\pm 1.0\% \pm 0.010$ OD or better and a repeatability specification of $\pm 1.0\% \pm 0.005$ OD or better.

The Synergy Neo 2 may be used if it has passed the Absorbance Plate Test or Absorbance Liquid Test 2, and Absorbance Liquid Test 1.

- Microplate shaker (if the absorbance reader does not support shaking)
- Precision balance with capacity of 100 g minimum and readability of 0.001 g
- 50–200 μ L hand pipette and disposable tips
- Deionized water
- Supply bottles
- 250-mL beaker
- New 96-well, clear, flat-bottom microplates
- BioTek's Green Test Dye Solution (PN 7773003) undiluted, or one of the alternate test solutions listed in the next section
- 100-mL graduated cylinder and 10-mL pipettes (if not using BioTek's Green Test Dye Solution)
- Gen5 software installed on the host PC
- Gen5 protocols as defined by the procedure on page 147

Test Solution Formulas

80 μ L of test solution with 150 μ L of deionized water should read between 1.300 and 1.700 OD at 405/750 nm. The solutions should be at room temperature.

If you do not have the Green Test Dye Solution (PN 7773003), prepare a dye solution using one of the following methods:

Using the Blue and Yellow Concentrate Dye Solutions:

Material	Quantity
Concentrate Blue Dye Solution (PN 7773001, 125 mL)	4.0 mL

Material	Quantity
QC (Yellow) Solution (PN 7120782, 125 mL)	5.0 mL
Deionized water	90.0 mL

Using FD&C Blue and Yellow Dye Powder:

Material	Quantity
FD&C Blue No. 1	0.200 grams
FD&C Yellow No. 5	0.092 grams
Tween 20	1.0 mL
Sodium Azide N ₃ Na	0.100 gram
Deionized water	Make to 1 liter

Procedure

If you have not already done so, create Gen5 protocols **Synergy Neo2_Disp 1 Test.prt** and **Synergy Neo2_Disp 2 Test.prt**. Instructions begin on page 147.

1. Prime both dispensers with 4000 µL of deionized or distilled water.
2. Remove the inlet tubes from the supply bottles. Prime both dispensers with the Volume set to 2000 µL. This prevents the water from diluting the dye.
3. Fill a beaker with at least 20 mL of the green dye solution. Prime both dispensers with 2000 µL of the solution. When finished, remove the priming plate from the carrier.
4. In Gen5, create an experiment based on **Synergy Neo2_Disp 1 Test.prt**.
5. Place a new 96-well microplate on the balance and tare the balance.
6. Place the plate on the microplate carrier.

Running a dispense procedure without placing a plate in the reader will result in contamination of the reader from spilled liquid.

When each dispense step is finished, you will weigh the plate, record the weight, tare the balance with the plate on it, and then place the plate back on the carrier for the next step.

7. Initiate a plate read. Gen5 prompts you to empty the tip priming trough.
8. When ready, click **OK** at the Load Plate dialog to begin the experiment. The sequence is as follows:
 - a. Dispense 80 μL /well to columns 1–4.
 - b. Remove the plate and weigh it. Record the weight and tare the balance.
 - c. Place the plate on the carrier and dispense 20 μL /well to columns 5–8.
 - d. Remove the plate and weigh it. Record the weight and tare the balance.
 - e. Place the plate on the carrier and dispense 5 μL /well to columns 9–12.
 - f. Remove the plate and weigh it. Record the weight.
 - g. Manually pipette 150 μL of deionized or distilled water into all 12 columns, on top of the green test dye solution.
 - h. Place the plate on the carrier for a 30-second shake, the “80 μL ” read at 405/750 nm, and the “20 and 5 μL ” read at 630/750 nm.
9. When finished, select **File > Save As**. Enter an identifying file name and click **Save**.
10. Remove the plate from the carrier and set it aside.
11. Repeat steps 4–9 using **Synergy Neo2_Disp 2 Test.prt**.
12. See page "Results Analysis" below for instructions on analyzing the results.

When all tests are complete, prime both dispensers with at least 5000 μL of deionized water to flush out the green dye solution.

Results Analysis

For your convenience, worksheets are included at the end of this chapter for recording the dispense weights, Delta OD values, calculations, and pass/fail.

The pass/fail criteria for each set of 32 wells with the same dispense volume is based on the calculated coefficient of variation (% CV) and Accuracy % Error.

For each volume dispensed (80 µL, 20 µL, 5 µL), for each dispenser (1, 2):

- Calculate the Standard Deviation of the 32 wells
- Calculate the Mean of the 32 wells
- Calculate the %CV: $(\text{Standard Deviation} / \text{Mean}) \times 100$
- Calculate the Accuracy % Error:
 $((\text{Actual Weight} - \text{Expected Weight}) / \text{Expected Weight}) * 100$

Expected Weights for 32 wells: 80 µL (2.560 g), 20 µL (0.640 g), 5 µL (0.160 g). It is assumed that one gram is equal to one milliliter.

Dispense Volume	To pass, %CV must be	To pass, Accuracy % Error must be:
80 µL	≤ 2.0%	≤ 2.0%
20 µL	≤ 7.0%	≤ 5.0%
5 µL	≤ 10.0%	≤ 20.0%

Failures

If any tests fail, prime the fluid lines and rerun the test(s).

- If the test(s) fail again, the injectors may require cleaning (see **Periodic Maintenance**).

If tests continue to fail, contact Technical Support.

Gen5 Test Protocols

This section contains instructions for creating two Gen5 protocols specifically for performing the Synergy Neo 2 Dispense Precision and Accuracy test.

1. Select **System > Instrument Configuration**, and add/configure the Synergy Neo 2 (if it is not already there).
2. Create a new protocol.

To edit a protocol category, double-click its “branch” in the tree.

3. Perform the steps in the following three sections to define the Procedure, customize the Plate Layout, and add Data Reduction steps, to test Dispenser #1.
4. When finished, select **File > Save As** and save the file as **Synergy Neo2_Disp 1 Test.prt**.
5. Repeat steps 2–4 above to create **Synergy Neo2_Disp 2 Test.prt** to test Dispenser 2.

Define the Procedure

In brief, the protocol's procedure follows the sequence below. After each Dispense step, the plate is ejected to allow the operator to weigh it and then tare the balance.

- Dispense 80 μ L dye to columns 1–4
- Dispense 20 μ L dye to columns 5–8
- Dispense 5 μ L dye to columns 9–12
- Shake the plate for 30 seconds
- Read columns 1–4 at 405/750 nm and calculate the Delta OD
- Read columns 5–12 at 630/750 nm and calculate the Delta OD

The detailed procedure is described on the next page. To add a step to the procedure, click the appropriate button on the left side of the Procedure dialog and define the required parameters.

The comments suggested for use with the Plate Out/In steps are optional, but they may be useful for the person running the test. When the Plate Out/In step is executed, Gen5 displays its comment in a message box.

Gen5 Procedure Steps	
Step Type	Details
1. Dispense	Dispenser <select 1 or 2, depending on the protocol> Dispense to wells A1..H4 Tip prime before this dispense step, 20 μ L Dispense 80 μ L at rate 275 μ L/sec
2. Plate Out,In	Suggested comment: Weigh the plate (80 μ L test). RECORD the weight, TARE the balance. Place the plate back on the carrier. Click OK to continue.

Gen5 Procedure Steps	
Step Type	Details
3. Dispense	<p>Dispenser <select 1 or 2, depending on the protocol></p> <p>Dispense to wells A5..H8</p> <p>Tip prime before this dispense step, 20 μL</p> <p>Dispense 20 μL at rate 250 μL/sec</p>
4. Plate Out,In	<p>Suggested comment: Weigh the plate (20 μL test). RECORD the weight and TARE the balance. Place the plate back on the carrier. Click OK to continue.</p>
5. Dispense	<p>Dispenser <select 1 or 2, depending on the protocol></p> <p>Dispense to wells A9..H12</p> <p>Tip prime before this dispense step, 5 μL</p> <p>Dispense 5 μL at rate 225 μL/sec</p>
6. Plate Out,In	<p>Suggested comment: Weigh the plate (5 μL test). RECORD the weight. PIPETTE 150 μL/well of DI water into all 12 columns. Place the plate back on the carrier. Click OK to perform the Read steps.</p>
7. Shake	<p>Orbital at 425 cpm (3 mm) for 30 seconds.</p>
8. Read	<p>Step label: "80 ul Read_Dispatch 1" (or _Disp 2)</p> <p>Wells: A1..H4</p> <p>Detection Method: Absorbance</p> <p>Read Type: Endpoint</p> <p>Read Speed: Normal</p> <p>Two Wavelengths: 405 and 750 nm</p>

Gen5 Procedure Steps	
Step Type	Details
9. Read	Step label: "20 and 5 ul Read_Displ 1" (or _Disp 2) Wells: A5..H12 Detection Method: Absorbance Read Type: Endpoint Read Speed: Normal Two Wavelengths: 630 and 750 nm

Customize the Plate Layout (Optional)

The results analysis worksheet at the end of this chapter requires the calculation of the Standard Deviation, Mean, and % CV of the ODs read for each dispense volume in each plate (six sets of calculations). By identifying the wells by their dispense volumes in the Plate Layout, Gen5 will calculate these values for you.

1. In the protocol, open the Plate Layout dialog.
2. Set the Type set to **Assay Control** and define three control types: Disp_80, Disp_20, and Disp_5.
3. In the Plate Layout, select **Disp_80** and highlight wells A1 to H4.
4. Select **Disp_20** and highlight wells A5 to H8.
5. Select **Disp_5** and highlight wells A9 to H12.
6. Click **OK** to save the changes and close the dialog.

After running the experiment, view the Statistics for each Delta OD Data Set to view the calculations

Add Data Reduction Steps

Each Read step is performed using two wavelengths, so you will create two data reduction steps to calculate the Delta OD values.

1. In the protocol, open the Data Reduction dialog and click **Custom. Transformation.**
2. Click **Select multiple data sets** and then select **DS2.**
3. Set the Data In for DS1 to the 80 μ L Read step at 405 nm.
4. Set the Data In for DS2 to the 80 μ L Read step at 750 nm.
5. Click **OK** to return to the Transformation dialog.
6. In the New Data Set Name field, type an identifying name such as **Delta OD 80 uL_Dispatch 1.**
7. Clear **Use single formula for all wells.**
8. In the Current Formula field, type DS1–DS2 and then highlight wells A1 to H4 to assign the formula.
9. Click **OK** to add the transformation to the Data Reduction list.
10. Create another Transformation similar to the above, with these characteristics:
 - DS1 set to the 20 and 5 μ L Read step at 630 nm
 - DS2 set to the 20 and 5 μ L Read step at 750 nm
 - New Data Set Name resembling **Delta OD 20 and 5 uL_Dispatch 1**
 - Remember to clear **Use single formula for all wells**
 - Formula DS1–DS2 applied to wells A5 to H12
11. When you are finished, the Data Reduction Steps list shows two Delta OD transformations:
12. Click **OK** to close the Data Reduction dialog.

Alpha Detection Test

This section applies only to models with the Alpha module.

The alpha laser has been factory-calibrated to meet specification. BioTek has developed a test protocol that can be used with AlphaScreen TruHits Kit, available from PerkinElmer (Mfg. #6760627), to verify the functionality of the alpha laser system.

The **Crosstalk** test is a measure of how well the optical system can distinguish the signals emitted from the well being read from those of any adjacent well. This test also determines the signal-to-noise ratio (SNR) of the test plate and verifies that the signal is at an acceptable level for the sample material used. The test is designed for use with AlphaScreen TruHits, and it is assumed that 96-well plates are used with 100 μ L well volumes.

After pipetting the solutions into the microplate, allow the plate to sit at room temperature for 15 minutes before starting the experiment. Cover the plate to prevent any light from hitting the beads.

Required Materials

Manufacturer part numbers are subject to change.

- Recommended test solution, AlphaScreen TruHits Kit, available from PerkinElmer (Mfg. #6760627)
- Buffer: Phosphate-Buffered Saline (PBS), pH 7.2–7.6 (e.g., Sigma tablets, Mfg. #P4417 or equivalent)
- Clean 96-well solid white microplate
- 1.5-mL conical-bottom, polypropylene sample tube
- Alpha filter cubes
 - Dual PMT system: 2 or equivalent (Alpha Top)
 - 1 or equivalent (Alpha Bottom)
- Gen5 protocol: "Gen5 Protocol Reading Parameters" on page 155

Test Solutions

AlphaScreen beads are light sensitive. All tests should be performed under subdued laboratory lighting of less than 100 lux.

1. Prepare the PBS buffer solution:
 - a. (Optional, but recommended) Use a 0.45-micron filter to filter 200 mL of deionized or distilled water.
 - b. Follow the manufacturer's instructions on the PBS packaging to create 200 mL, dissolving the necessary amount of PBS into the filtered water.
 - c. Stir the solution (preferably using a stir table) until the PBS is completely dissolved.
 - d. Check the pH; it should be between 7.2 and 7.6 at 25°C.
2. Prepare the TruHits bead suspension:

- a. Shake the container of TruHits bead suspension vigorously for 30 seconds prior to pipetting. Alternatively, sonicate or vortex the container.
- b. Pipette the Donor and Acceptor beads into the bottom of a 1.5 mL conical-bottom, polypropylene sample tube
 - Donor beads 8 μ L
 - Acceptor beads 8 μ L
- c. Replace the tube cap and cover the tube with foil.
- d. Incubate for 15 minutes at room temperature.
- e. Add diluent to a final volume of 1 mL.
 - PBS 984 μ L
- f. Replace the tube cap and cover the tube with foil.
- g. Incubate for 15 minutes at room temperature.
- h. This yields a TruHits bead suspension that is 20 μ g/mL with respect to the Acceptor beads.
- i. Discard the remaining bead solution.
- j. If extra beads have been reacted, they can be stored for later use. For storage, refrigerate any unused portions of the TruHits bead suspension. The temperature must be between +2°C and +6°C.

Do **NOT** use sodium azide as a preservative as it will affect the AlphaScreen signal. Procline™ 300 at 0.03% is suggested.

Pipette Map

	1	2	3	4	5	6	7	8	9	10	11	12
A	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF
B	BUF	OMB1	BUF	BUF	OMB1	BUF	BUF	OMB1	BUF	BUF	OMB1	BUF
C	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF
D												
E												
F												
G	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF
H	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF	BUF

Procedure

Crosstalk Test

1. Pipette 100 μ L of PBS in wells A1–A12, B1, B3–B4, B6–B7, B9–B10, B12, C1–C12, and G1–H12 (see pipette map, “BUF” wells).
2. Pipette 100 μ L of 20 μ g/mL AlphaScreen beads suspension into wells B2, B5, B8, and B11 (see pipette map, “20 μ g/mL OMB” wells).

Allow the plate to sit at room temperature for approximately 15 minutes prior to use.

3. Create an experiment based on the **Synergy Neo2 AlphaTest_Crosstalk.prt** protocol. Read the plate, and then save the experiment.
4. Open the Plate menu and export the data to the embedded Power Export spreadsheet. The spreadsheet reports pass or fail for the test performed. See **Results Analysis**, next, for descriptions of the calculations and troubleshooting tips.
5. Print the spreadsheets and store them with your test records.

Results Analysis

1. Calculate the crosstalk for each of the four wells of AlphaScreen beads solution by dividing the background-subtracted Mean value of the surrounding adjacent wells by the background-subtracted AlphaScreen beads suspension well.
2. Average the % crosstalk of the four test wells to determine level of crosstalk (Crosstalk Mean).
3. Verify that the % crosstalk is less than or equal to 0.1%
4. Calculate the signal-to-noise ratio (SNR) by using the following equation:
$$\text{SNR} = (\text{signal mean} - \text{background mean}) / (\text{SQRT}(\text{signal STD}^2 + \text{background STD}^2))$$
5. Verify that SNR is greater than or equal to **10**.

Troubleshooting Alpha Tests

If the test fails, please try the following suggestions. If the test(s) continue to fail, print the results and contact Technical Support.

- Are the solutions fresh?
- Have the solutions been stored properly (between +2°C and +6°C)?
- Has the kit been exposed to excessive light (in excess of 100 lux)?
- If the Crosstalk test continues to fail, the laser may not be firing. Contact Technical Support.

Gen5 Protocol Reading Parameters

The information in the following table represents the recommended reading parameters. It is possible that your tests will require modifications to some of these parameters, such as Plate Type or Gain value (see **Troubleshooting Tips** above).

The Plate Type setting in the Gen5 protocol must match the plate you are actually using.

Synergy Neo 2 AlphaTest_Crosstalk.prt

This procedure contains one read step and calculates crosstalk based on the full plate data.

Parameter	Default Setting
Plate Type	Costar 96-well white opaque
Detection Method	Alpha
Read Type	Endpoint
Read Wells	Full plate
Gain	120
Delay after plate movement	0 msec
Excitation time	100 msec
Delay after excitation	120 msec
Integration time	100 msec
Read height	7.00 mm

Specifications

This appendix contains the published specifications for the Synergy Neo 2.

General Specifications	158
Dispense/Read Specifications	159
Absorbance Specifications	161
Fluorescence Specifications (Mono-Based)	163
Fluorescence Specifications (Filter-Based)	164
Luminescence Specifications	166
Alpha Detection Specifications	166

General Specifications

Microplates	
The Synergy Neo 2 accommodates standard 6-, 12-, 24-, 48-, 96-, 384-, and 1536-well microplates with 128 x 86 mm geometry and the Take3 and Take3 Trio Micro-Volume Plates.	
Maximum Plate Height: 0.89"	
Hardware and Environmental	
Light Source	
Absorbance, Fluorescence (FI), monochromator-based:	Xenon flash light source, 20W maximum average power
Fluorescence (FI/FP), filter-based:	Xenon flash light source, 5W maximum average power
TRF xenon flash filter-based:	Xenon flash light source, 5W maximum average power
TRF laser filter-based:	Nitrogen (N ₂) laser light source
With all modules installed, without power supply or dispense module attached:	
Weight:	approximately 100 lbs. (45.36 kg)
Dimensions	16.25" H x 15" W x 24.5" D 46.4 cm H x 38.1 cm W x 62.3 cm D
Environment:	Power-up temperature 18° to 30°C Operational temperature 18° to 40°C Storage temperature -25°C to 50°C
Humidity:	Operational: 10% to 85% relative humidity (non-condensing) Storage: 10% to 80% relative humidity (non-condensing)
Power:	Powered from an external 250W (minimum), 24VDC power supply compatible with 100-240 volts AC @50-60Hz

Microplates	
Incubation:	<p>AlphaScreen only: Temperature control range from 3°C over ambient to 30°C.</p> <p>All other modes: Temperature control range from 4°C over ambient to 65°C (70°C for models with 70°C support, indicated in the instrument system test report). For incubation setpoints > 65°C, ambient room temperature must be at least 20°C.</p> <p>Temperature variation $\pm 0.5^\circ\text{C}$ across the plate @ 37°C, tested with Innovative Instruments, Inc. temperature test plate</p>

Dispense/Read Specifications

Maximum Delay between End of Dispense and Beginning of Read 96/384-Well Plates, Default Probe Heights	
Bottom Mono Fluorescence	T < 1.0 second
Luminescence (filter or broadband)	T < 1.0 second
Top Filter Fluorescence	T < 1.0 second
Top Mono Fluorescence	T < 1.0 second
Absorbance	T < 1.0 second
Bottom Filter Fluorescence	T < 1.0 second

Dispense/Read, for Models with the Dual-Reagent Dispense Module	
Plate Type	Both injectors dispense to standard height 6-, 12-, 24-, 48-, 96-, and 384-well microplates.
Detection Method	Absorbance, Fluorescence (FI, FP, TRF), Luminescence
Volume Range	5–1000 μL with a 5–20 μL tip prime
Reagent Dead Volume	< 1100 μL , with dead volume recovery function (back flush)
Injection Speeds	225, 250, 275, 300 $\mu\text{L}/\text{second}$
Accuracy	$\pm 1\mu\text{L}$ or 2.0%, whichever is greater

Dispense/Read, for Models with the Dual-Reagent Dispense Module	
Precision	< 2.0% for volumes of 50–200 μL < 4.0% for volumes of 25–49 μL < 7.0% for volumes of 10–24 μL < 10.0% for volumes of 5–9 μL

Absorbance Specifications

Note: For the performance specifications described in this section, the gain on the optics test should be ≤ 8 .

Optics	
Wavelength Range	230 to 999 nm
Wavelength Bandpass	< 4 nm (230–285 nm), < 8 nm (> 285 nm)
Measurement Range	0.000 to 4.000 OD
Resolution	0.0001 OD
Increment	1 nm
Wavelength Accuracy	± 2 nm
Wavelength Precision	± 0.2 nm
Minimum kinetic interval (450 nm)	< 12 seconds, sweep mode, 384-well microplate

Plate In/Plate Out Speed

< 20 seconds, 450 nm, sweep mode, 384-well microplate

Accuracy, Linearity, Repeatability

Specifications apply from 250–999 nm, 200 μ L (96-well microplates) and 80 μ L (384-well microplates).

Accuracy (tested with certified neutral density glass)

96-well plate, normal read speed

0–2 OD: $\pm 1\%$ ± 0.010 OD, delay after plate movement = 100 ms

2–2.5 OD: $\pm 3\%$ ± 0.010 OD, delay after plate movement = 100 ms

384-well plate, normal read speed

0–2 OD: $\pm 2\%$ ± 0.010 OD, delay after plate movement = 100 ms

2–2.5 OD: $\pm 5\%$ ± 0.010 OD, delay after plate movement = 100 ms

Accuracy, Linearity, Repeatability	
96-well and 384-well plate, sweep read speed 0–1 OD: +/-2% +/-0.010 OD	
Linearity (by liquid dilution)	
96-well plate, normal read speed 0–2 OD: +/-1% +/-0.010 OD, delay after plate movement = 100 ms 2–2.5 OD: +/-3% +/-0.010 OD, delay after plate movement = 100 ms	
384-well plate, normal read speed 0–2 OD: +/-2% +/-0.010 OD, delay after plate movement = 100 ms 2–2.5 OD: +/-5% +/-0.010 OD, delay after plate movement = 100 ms	
96-well and 384-well plate, sweep read speed 0–1 OD: +/-2% +/-0.010 OD	
Repeatability (tested with certified neutral density glass/measured by one standard deviation: 8 measurements per data point)	
96-well and 384-well plate, normal read speed 0–2 OD: +/-1% +/-0.005 OD, delay after plate movement = 100 ms 2–2.5 OD: +/-3% +/-0.005 OD, delay after plate movement = 100 ms	
96-well and 384-well plate, sweep read speed 0–1 OD: +/-2% +/-0.010 OD	
Take3 Plate	
Detection Limit, 260 nm dsDNA	< 5 ng/ μ L

Fluorescence Specifications (Mono-Based)

The Synergy Neo 2 measures fluorescence with monochromators from the top and bottom of 6- to 1536-well plates.

Monochromator-Based Fluorescence	
Excitation range	250-700 nm with low-noise PMT 250–850 nm with red-shifted PMT
Emission range	250-700 nm with low-noise PMT 250–850 nm with red-shifted PMT; 300 nm minimum with EM scan.
Bandpass	Excitation and Emission: Variable, from 3 nm–50 nm (wavelength dependent) in 1 nm increments

Optical Probes	
Top position	2.00 mm diameter fixed
Bottom position	2.00 mm diameter fixed

Plate In/Plate Out Speed
< 20 seconds for filter set, sweep mode, 384-well microplate

Sensitivity	
96-/384-well plates	<p><i>Sodium Fluorescein in phosphate buffered saline (PBS)</i> DL <= 20 pM top or bottom read Excitation 485 nm, Emission 528 nm</p> <p><i>Methylumbelliferone (MUB) in carbonate-bicarbonate buffer (CBB)</i> DL <= 0.16 ng/mL (0.91 nM) top read Excitation 360 nm, Emission 460 nm</p> <p><i>Propidium Iodide (PI) in PBS</i> DL <= 62.5 ng/mL bottom read Excitation 485 nm, Emission 645 nm</p>

Fluorescence Specifications (Filter-Based)

The Synergy Neo 2 measures fluorescence with filters from the top and bottom of 6- to 1536-well plates. The xenon flash bulb is used for TRF flash reads except when the TRF laser is selected in T-model instruments. The T-models use an N2 laser to perform TRF using filters and top probe.

Optical Probes	
Top position	2.6-mm diameter fixed
Bottom position	2.6-mm diameter fixed

Top Single PMT

Plate In/Plate Out Speed
< 20 seconds for filter set, sweep mode, 384-well microplate

Sensitivity	
96-/384-well plates	DL < 5 pM solution of Sodium Fluorescein in PBS Excitation 485/20, Emission 528/20, 510 nm mirror DL ≤ 0.16 ng/mL (0.91 nM) solution of Methylumbelliferone in CBB, Excitation 360/40, Emission 460/40, 400 nm mirror

Time-Resolved Fluorescence Flash	
96-/384-well plates	DL Europium ≤ 250 fM Excitation 360/40 nm, Emission 620/40 nm 400 nm mirror
Integration Interval	20 to 2000 μs
Delay	0 to 2000 μs
Granularity	1-μs step

T-Models with TRF Laser only: Gen5 version 3.06 or higher.

Time-Resolved Fluorescence Laser	
96-/384-well plates	DL Europium ≤ 20 fM (50 measurements/well - 384-well plate, 100 μ L/well); DL Europium ≤ 40 fM (50 measurements/well - 96-well plate, 200 μ L/well)
Delay	250 μ sec
Collection time	1000 μ sec
HTRF (Cisbio®)	HTRF Assays successfully performed with white and black microplates

Top Double PMT

Plate In/Plate Out Speed	
< 20 seconds for filter set, sweep mode, 384-well microplate	

Fluorescence Polarization	
96-/384-well plate	< 5 mP standard deviation at 1 nM Sodium Fluorescein Excitation 485/20 nm, Emission 528/20 nm, 510 nm mirror

Bottom Optics

Sensitivity	
96-/384-well plates	DL < 5 pM solution of Sodium Fluorescein in PBS Excitation 485/20, Emission 528/20, 510 nm mirror DL ≤ 0.16 ng/mL (0.91 nM) solution of Methylumbelliferone in CBB, Excitation 360/40, Emission 460/40, 400 nm mirror

Luminescence Specifications

The Synergy Neo 2 measures luminescence from the top of 6- to 384-well plates. The following requirements apply to 96-well plates with 200 μ L/well, at room temperature. Product testing is performed using a Harta plate.

Luminescence - Filter-based detection	
96-well plate	≤ 50 amol/well—standard low-noise PMT ≤ 500 amol/well—red-shifted PMT
Integration Time	10 seconds
Blank Wells	16
Cross-talk	$< 0.1\%$ in 96-well plate $< 0.25\%$ in low-volume 384-well plate

Synergy Neo 2 models with this capability require Gen5 version 3.06 or higher:

Luminescence - Broadband fiber detection	
96-well plate	≤ 75 amol/well —standard low-noise PMT ≤ 500 amol/well—red-shifted PMT
Integration Time	10 seconds
Blank Wells	16

Alpha Detection Specifications

Alpha Laser	
Laser Output Power	100 mW \pm 10%

Appendix B

Error Conditions

This appendix lists and describes Synergy Neo 2 error codes that may appear in Gen5.

Error Codes Overview	168
Error Codes	168

Error Codes Overview

When a problem occurs during operation with the Synergy Neo 2, an error message appears in Gen5. Error messages typically contain four characters, such as “4168,” and in most cases are accompanied by descriptive text, such as “PMT overload error.” With many errors, the instrument will beep repeatedly; press the carrier eject button to stop this alarm.

Some problems can be solved easily, such as “2B0A: Priming plate not detected” (place a priming plate on the carrier). Some problems can be solved only by trained Technical Support personnel. This appendix lists the most common and easily resolved error codes that you may encounter.

Error codes beginning with “A” (e.g., A100) indicate conditions that require immediate attention. If this type of code appears, turn the instrument off and on. If the System Test does not conclude successfully, record the error code, note the instrument serial number and contact Technical Support.

If an error code appears in Gen5, you may want to run a System Test for diagnostic purposes. In Gen5, select **System > Diagnostics > Run System Test**.

If an error message appears while an experiment is in process and after having received measurement data, it is your responsibility to determine if the data is valid.

Technical Support

Use this appendix to diagnose problems and solve them if possible. If you need further assistance, contact Technical Support; see page vi.

For errors that are displayed during operation of the Synergy Neo 2 with the stacker or the BioStack 4, refer to the stacker's or Biostack's user manuals.

Error Codes

This table lists the most common and easily resolved error codes that you may encounter. If an error code appears in Gen5, look for it here. If you find the code, follow the suggestions provided for solving the problem. If you cannot find the code or if you are unable to solve the problem, please contact Technical Support. The Gen5 Help system also provides troubleshooting tips.

Code	Description and Possible Remedy
2353	<p>Expected plug/hole/filter not found in filter cube</p> <p>This error indicates that the filter cube is not installed and is required for the read. Install the filter cube or check that it is installed correctly.</p>
2700	<p>Error attempting to run the barcode scanner SET command</p> <p>A response returned from the scanner is invalid.</p>
2701	<p>Error attempting to run the barcode scanner SET command</p> <p>The command message is calling out an invalid barcode location. Valid numbers are 1–4.</p>
2702	<p>Error attempting to get barcode scanner information regarding one of the barcode types</p> <p>The barcode type returned by the scanner is not one of those expected.</p>
2703	<p>Barcode type is not supported</p> <p>One of the four barcode types is not supported by the scanner.</p>
2704	<p>Error disabling start/stop character transmission</p> <p>While attempting to tell the scanner to disable the transmission of start and stop characters along with the barcode value for the Codebar barcode type, an error occurred.</p>
2B0A	<p>Priming plate not detected</p>
2B0x	<p>Dispenser syringe 1 or 2 (respectively) did not home</p> <p>x=1–3</p> <p>Generally, this error indicates the syringe was not properly installed. Make sure the syringe's thumbscrews are properly threaded. (Refer to the Installation chapter for instructions.) Restart the reader.</p>
2B04	<p>Dispenser syringe 1 or 2 (respectively) failed position verify</p> <p>Generally, this error indicates the syringe was not properly installed. Make sure the syringe's thumbscrews are properly threaded. Restart the reader. (Refer to the Installation chapter for instructions.)</p>

Code	Description and Possible Remedy
2D15	<p>Kinetic interval is out of range.</p> <p>Valid range is 9 to 9999 seconds.</p>
37x0/47x0 38x0/48x0 39xy/49xy	<p>Noise Test Errors Offset Test Errors Dark Range Errors</p> <p>x=0, 1; y=0–6</p> <p>This series of System Test errors may indicate too much light inside the chamber. Make sure the plate carrier door and the front hinged door are properly closed. For models with the dispense module, if the dispense tubes are not connected to the reader, re-install the light shield that shipped with the instrument (or cover the hole with black tape). Restart the reader.</p>
4Fxy	<p>Fluorescence signal out of range.</p> <p>x=0, 1; y=0–6</p> <p>Too low of a reading indicates a light signal problem. Ensure that Gen5 Fluor/Lum wavelengths table matches the actual filter installed in the filter cube.</p>
40xx	<p>PMT overload well error at <well #xxx></p> <p>This error typically means that the fluid in a well has oversaturated the PMT (i.e., the well is too bright). Try lowering the sensitivity value in the read step.</p> <p>To identify the well:</p> <p>Wells are counted starting at A1, moving left-to-right, row-by-row. The row and column of the well can be extracted from the well number code by applying the following formula (example uses 8 x 12 geometry, 96-well plate):</p> <ol style="list-style-type: none"> 1. Convert the ASCII hex string to a decimal equivalent. Ex: “057” indicates 57 hex, yielding a well code of 87 decimal. 2. Row = (well code) / (columns in plate), rounded up to a whole number. Ex: 87/12 = 7.25, indicating row 8 (or H). 3. Column = (well code) - ((row-1) * (columns in plate)). Ex: 87 - ((8 - 1) * 12) = column 3.

Code	Description and Possible Remedy																													
4Exy	<p>Detector saturated (too much light). Relative Fluorescing Units (RFU) reached (99999).</p> <p>X = Fluorescence channel</p> <p>Fluorescence/Luminescence Channel Codes</p> <table border="0" data-bbox="548 541 1224 856"> <tr> <td>Mono optics reference channel</td> <td>0</td> </tr> <tr> <td>Mono optics PMT channel</td> <td>1</td> </tr> <tr> <td>Top filter optics reference channel</td> <td>2</td> </tr> <tr> <td>Top filter optics, top PMT channel</td> <td>3</td> </tr> <tr> <td>Top filter optics, side PMT channel</td> <td>4</td> </tr> <tr> <td>Bottom filter optics reference channel</td> <td>5</td> </tr> <tr> <td>Bottom filter optics PMT channel</td> <td>6</td> </tr> </table> <p>Y = PMT Test Type Code</p> <table border="0" data-bbox="391 968 1317 1188"> <tr> <td>Connection Test</td> <td>0</td> <td>PMT not connected</td> </tr> <tr> <td>High Voltage Test</td> <td>1</td> <td>Failure during test at higher voltage</td> </tr> <tr> <td>Low Voltage Test</td> <td>2</td> <td>Failure during test at lower voltage</td> </tr> <tr> <td>Well Overload Test</td> <td>5</td> <td>Failure during test at well</td> </tr> <tr> <td>Background Overload Test</td> <td>8</td> <td>Failure during background overload test</td> </tr> </table> <p>OR Y = filter readset</p> <p>This error can indicate one of several scenarios. It is possibly due to incorrect chemistry, e.g., the fluorescence standards dispensed to the plate exceed expectations.</p> <p>For models with the dispense module, the internal chamber may require cleaning (contact Technical Support).</p>	Mono optics reference channel	0	Mono optics PMT channel	1	Top filter optics reference channel	2	Top filter optics, top PMT channel	3	Top filter optics, side PMT channel	4	Bottom filter optics reference channel	5	Bottom filter optics PMT channel	6	Connection Test	0	PMT not connected	High Voltage Test	1	Failure during test at higher voltage	Low Voltage Test	2	Failure during test at lower voltage	Well Overload Test	5	Failure during test at well	Background Overload Test	8	Failure during background overload test
Mono optics reference channel	0																													
Mono optics PMT channel	1																													
Top filter optics reference channel	2																													
Top filter optics, top PMT channel	3																													
Top filter optics, side PMT channel	4																													
Bottom filter optics reference channel	5																													
Bottom filter optics PMT channel	6																													
Connection Test	0	PMT not connected																												
High Voltage Test	1	Failure during test at higher voltage																												
Low Voltage Test	2	Failure during test at lower voltage																												
Well Overload Test	5	Failure during test at well																												
Background Overload Test	8	Failure during background overload test																												
2D46	<p>Fluorescence wavelength not found in table</p> <p>This error indicates that the wavelength specified in the procedure is not detected in the instrument's filter table. In Gen5, verify the Fluorescence filter table has the wavelengths loaded into the reader. Compare the contents of the table with the filters installed in the filter cube (see the Gen5 Help system for more information). Restart the reader.</p>																													

Code	Description and Possible Remedy						
50xx 510x	<p data-bbox="386 296 646 327">Axis failed to home</p> <table data-bbox="418 365 1110 495"> <tr> <td data-bbox="418 365 1068 396">Top filter optics, lower filter/mirror slider</td> <td data-bbox="1068 365 1110 396">03</td> </tr> <tr> <td data-bbox="418 411 1068 443">Top filter optics, upper filter/mirror slider</td> <td data-bbox="1068 411 1110 443">04</td> </tr> <tr> <td data-bbox="418 457 1068 489">Bottom filter optics, filter/mirror slider</td> <td data-bbox="1068 457 1110 489">06</td> </tr> </table> <p data-bbox="386 531 1336 646">Generally, this error indicates the filter cube is not seated properly in the reader. Remove it, ensure each filter or plug is properly positioned and reinstall it securely. Restart the reader.</p>	Top filter optics, lower filter/mirror slider	03	Top filter optics, upper filter/mirror slider	04	Bottom filter optics, filter/mirror slider	06
Top filter optics, lower filter/mirror slider	03						
Top filter optics, upper filter/mirror slider	04						
Bottom filter optics, filter/mirror slider	06						
540x	<p data-bbox="386 680 834 711">Filter cube failed positional verify</p> <table data-bbox="418 749 1110 879"> <tr> <td data-bbox="418 749 1068 781">Top filter optics, lower filter/mirror slider</td> <td data-bbox="1068 749 1110 781">03</td> </tr> <tr> <td data-bbox="418 795 1068 827">Top filter optics, upper filter/mirror slider</td> <td data-bbox="1068 795 1110 827">04</td> </tr> <tr> <td data-bbox="418 842 1068 873">Bottom filter optics, filter/mirror slider</td> <td data-bbox="1068 842 1110 873">06</td> </tr> </table> <p data-bbox="386 915 1336 1031">Generally, this error indicates the filter cube is not seated properly in the reader. Remove it, ensure each filter or plug is properly positioned and reinstall it securely. Restart the reader.</p>	Top filter optics, lower filter/mirror slider	03	Top filter optics, upper filter/mirror slider	04	Bottom filter optics, filter/mirror slider	06
Top filter optics, lower filter/mirror slider	03						
Top filter optics, upper filter/mirror slider	04						
Bottom filter optics, filter/mirror slider	06						

Code	Description and Possible Remedy																																				
55xy	<p><Motor> not homed successfully</p> <p>xy=axis</p> <p>Axis Codes</p> <table data-bbox="444 499 1305 1339"> <tbody> <tr><td>Carrier X</td><td>00</td></tr> <tr><td>Carrier Y</td><td>01</td></tr> <tr><td>Top probe Z</td><td>02</td></tr> <tr><td>Top filter optics, lower filter/mirror slider</td><td>03</td></tr> <tr><td>Top filter optics, upper filter/mirror slider</td><td>04</td></tr> <tr><td>Bottom probe Z</td><td>05</td></tr> <tr><td>Bottom filter optics, filter/mirror slider</td><td>06</td></tr> <tr><td>Mono optics, excitation filter sector</td><td>07</td></tr> <tr><td>Mono optics, emission filter sector</td><td>08</td></tr> <tr><td>Mono optics, excitation slit wheel</td><td>09</td></tr> <tr><td>Mono optics, emission slit wheel</td><td>0A</td></tr> <tr><td>Mono optics, excitation probe changer</td><td>0B</td></tr> <tr><td>Mono optics, emission probe changer</td><td>0C</td></tr> <tr><td>Excitation monochromator grating</td><td>0D</td></tr> <tr><td>Emission monochromator grating</td><td>0E</td></tr> <tr><td>Dispenser syringe 1</td><td>0F</td></tr> <tr><td>Dispenser syringe 2</td><td>10</td></tr> <tr><td>N2/XE selector</td><td>11</td></tr> </tbody> </table> <p>This error indicates that an axis failed a previous verify function and now needs to be homed. Check for any obstructions that may prevent the carrier, syringes, or filter cube from moving normally. Restart the reader.</p>	Carrier X	00	Carrier Y	01	Top probe Z	02	Top filter optics, lower filter/mirror slider	03	Top filter optics, upper filter/mirror slider	04	Bottom probe Z	05	Bottom filter optics, filter/mirror slider	06	Mono optics, excitation filter sector	07	Mono optics, emission filter sector	08	Mono optics, excitation slit wheel	09	Mono optics, emission slit wheel	0A	Mono optics, excitation probe changer	0B	Mono optics, emission probe changer	0C	Excitation monochromator grating	0D	Emission monochromator grating	0E	Dispenser syringe 1	0F	Dispenser syringe 2	10	N2/XE selector	11
Carrier X	00																																				
Carrier Y	01																																				
Top probe Z	02																																				
Top filter optics, lower filter/mirror slider	03																																				
Top filter optics, upper filter/mirror slider	04																																				
Bottom probe Z	05																																				
Bottom filter optics, filter/mirror slider	06																																				
Mono optics, excitation filter sector	07																																				
Mono optics, emission filter sector	08																																				
Mono optics, excitation slit wheel	09																																				
Mono optics, emission slit wheel	0A																																				
Mono optics, excitation probe changer	0B																																				
Mono optics, emission probe changer	0C																																				
Excitation monochromator grating	0D																																				
Emission monochromator grating	0E																																				
Dispenser syringe 1	0F																																				
Dispenser syringe 2	10																																				
N2/XE selector	11																																				
570x	<p>Area Scan matrix too large for perimeter wells</p> <p>x=0, 1</p> <p>For some plate type and read probe combinations, it might not be possible to define the entire area scan matrix offered by Gen5 for some perimeter wells, due to the physical limitations of carrier travel. Redefine the area scan to include a smaller matrix or select wells in a different row or column.</p>																																				

Code	Description and Possible Remedy
5A0x	<p>Plate could not be moved inside</p> <p>x=0, 1</p> <p>Make sure the Plate Type defined in the Gen5 Protocol matches the plate you are using. This error can also occur if the plate type is correct but the lid was left on the plate. If you wish to read the plate with a lid on it, create a new plate type in Gen5 with the correct Plate Height. This error can also indicate that the plate is not seated in the carrier correctly. Remove the plate and replace it securely in the carrier.</p> <p>Verify the tip priming trough has not become dislodged. Remove any plates from the carrier and power cycle the instrument to see if the error is resolved.</p>
5B00	<p>Required carrier in when expected to be outside</p> <p>Carrier is inside the chamber when it should be outside. This may occur if the read was aborted and "home all axes" not performed.</p> <p>This error can also occur if the carrier is inside and the newly defined plate height is different from the most recently specified plate height.</p> <p>To resolve the error, eject the carrier prior to running the experiment.</p>
70xx	<p>Motor synchronization error</p> <p>If at any time it is discovered that a motor is currently on the wrong microstep boundary for a move of a specified microstep size, a motor synch error is flagged.</p> <ul style="list-style-type: none"> • Verify there is no binding of axis movement. • Turn the instrument off and on to see if the error is repeatable. • Reload the basecode.

Microplate Barcode Scanner

This appendix contains instructions for installing the optional microplate barcode scanner and specifications for microplate barcode label format and positioning. It also tells you how to set up Gen5 to view and use the scanned barcodes.

Overview	176
Install the Microplate Barcode Scanner	176
About Barcode Labels	179
Placing a Barcode Label	180
Microplate Scanner Diagnostic Utility	182

Overview

The external barcode scanner allows you to automate the reading of microplate barcode labels in a robotic system. If a valid barcode is read, the barcode value is automatically passed to Gen5 for storage in the Gen5 experiment.

WARNING



Laser Beam. Serious eye injury may occur if you stare directly into the laser beam of the barcode scanner during operation of the scanner. This hazard is noted by the symbol shown here. Do not look directly into the laser beam during operation of the scanner.

The barcode scanner can be installed on either side of or in front of the reader, depending on which side of the microplate the barcode label is located.

Package Contents

PN 1030008 contains:

- Barcode scanner assembly: scanner (housed in a protective cover), scanner cable, and mirror attached to a removable mounting bracket
- Scanner alignment bracket




Install the Microplate Barcode Scanner

Be very careful not to scratch or damage the scanner's mirror when unpacking or installing the microplate barcode scanner! If you are using a stacker, install the stacker with the interfacing instrument before installing the barcode scanner. See Chapter 2, Installation.

Make sure the reader is powered off during installation of the barcode scanner.


The barcode scanner can be placed to the left, the right, or in front of the reader's carrier. Select the mounting position based on the location of the barcode labels on your microplates.

1. Lift the front of the Synergy Neo 2, place the mounting bracket under the locating pins, and carefully lower the instrument into the mounting bracket.

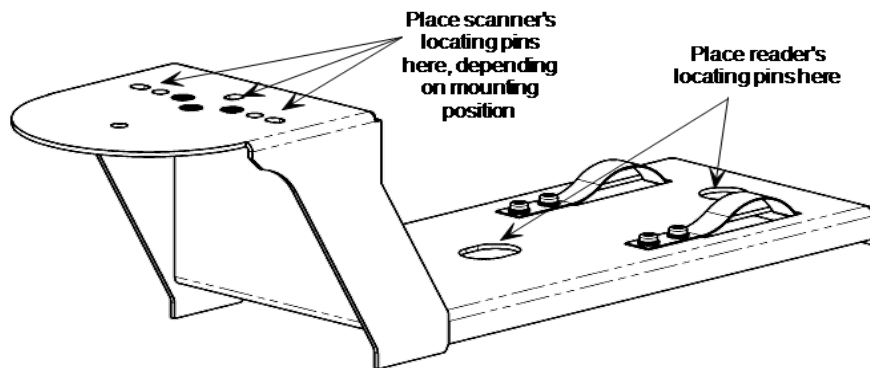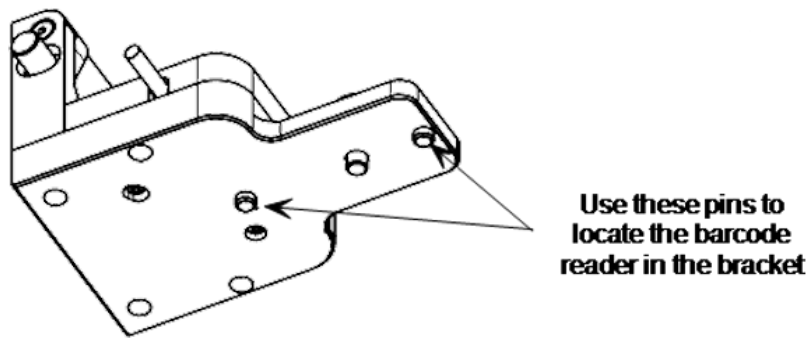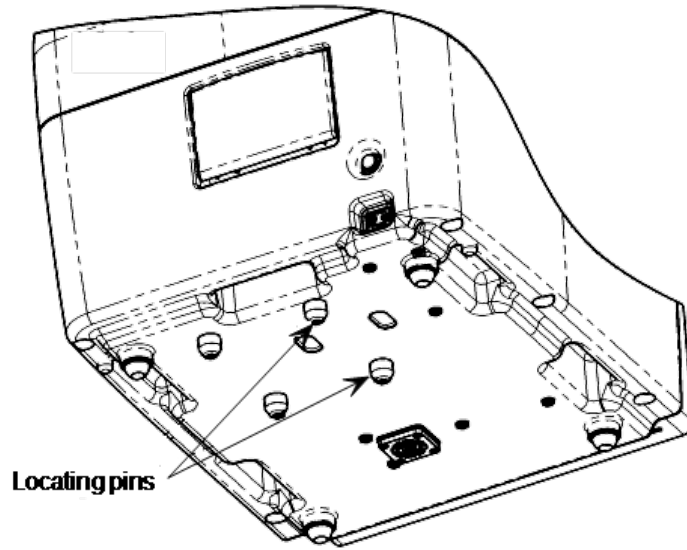


2. Determine the appropriate mounting position for the barcode scanner: left, front, or right.

3. Mount the barcode scanner in the bracket, then tighten the thumbscrew.

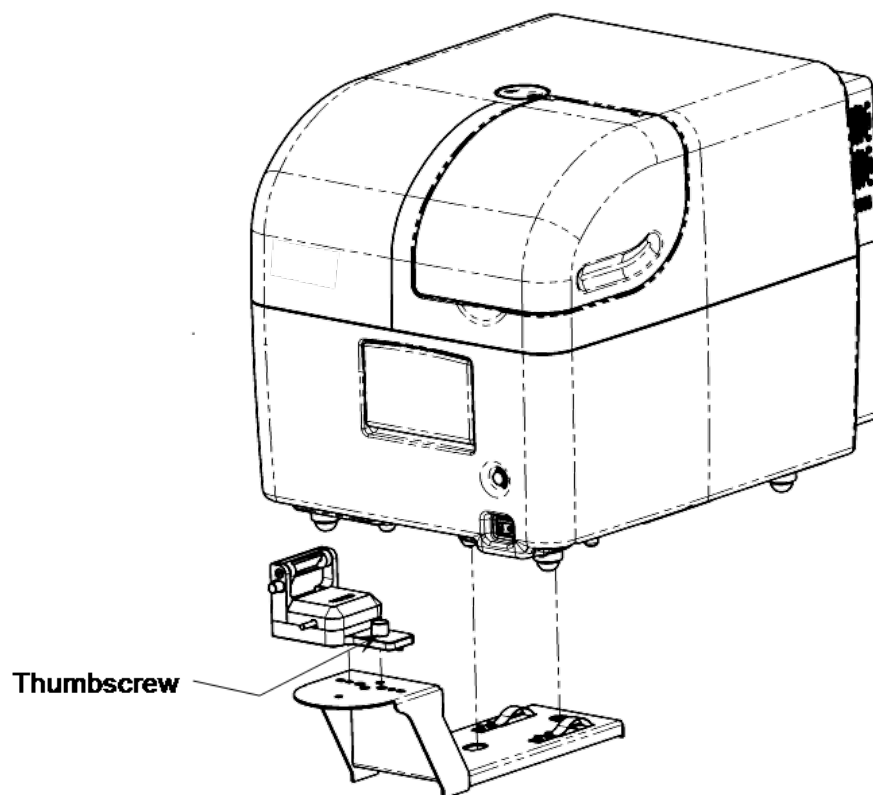

4. Plug one end of the communications cable into the barcode reader and connect the other end to the Synergy Neo 2.

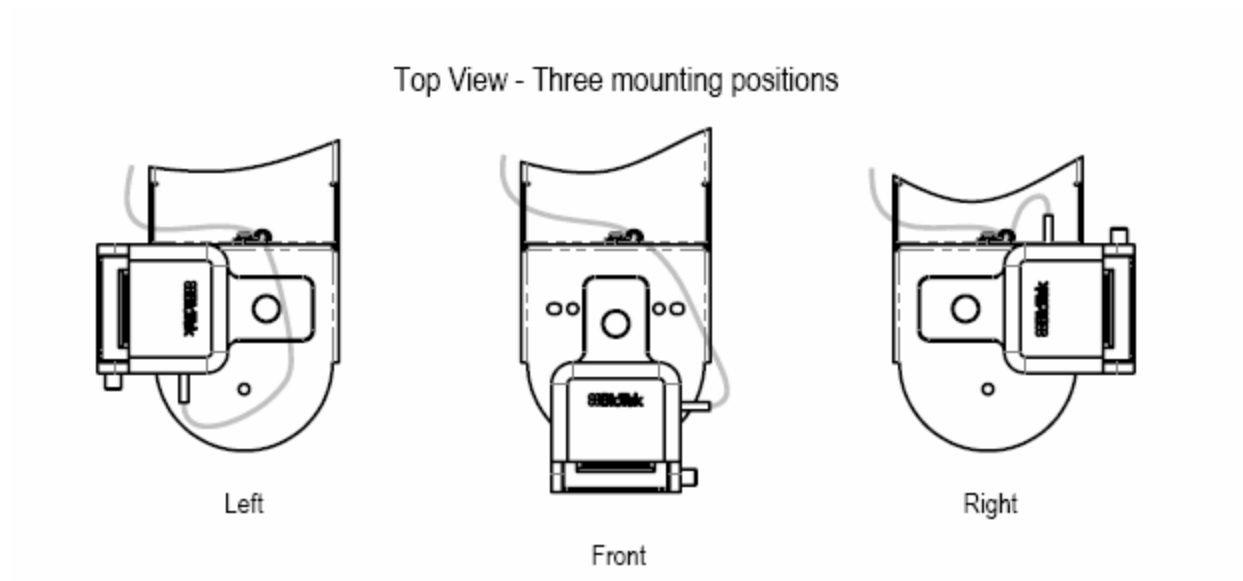



About Barcode Labels

Barcode labels are available for purchase, or they may be created using barcode software and label products that meet the following specifications:

- **Label format (“symbology”):** Required formats are:
 - Code 128
 - Code 39
 - Datamatrix
 - PDF 417
 - Codebar
 - UPS
 - QR
 - GS1 Databar
- **Industry Standards:** The labels should be made in accordance with Automation Identification Manufacturers (AIM) uniform symbol specification for all codes supported. Label decodability should meet ANSI Specification X2, 182-199 “Bar-Code Print Quality Guideline” for a rating of A/05/880.
- **Label quality:** The labels should be printed on good-quality copier machines or laserjet printers.
- **Label position on plate:** The labels should not extend above or below the edges of the plate, because the plates may stick to one another when they are stacked.


See "**Placing a Barcode Label**" on the next page on page 180 for more information about where to place the label on the microplate.

The following figures contain barcode label artwork that may be submitted to a print vendor for creating labels. The areas marked as quiet zones are critical for the reader to locate the barcode image.

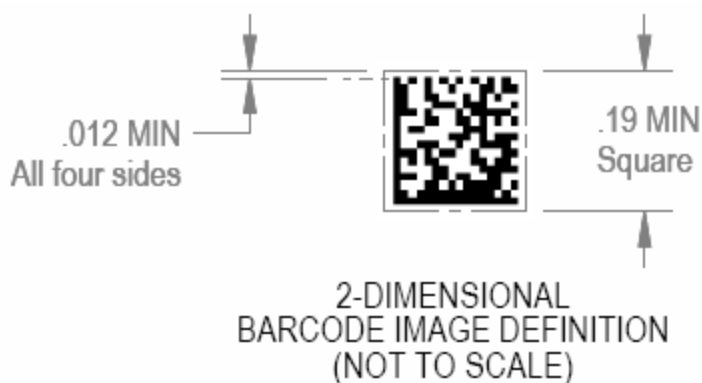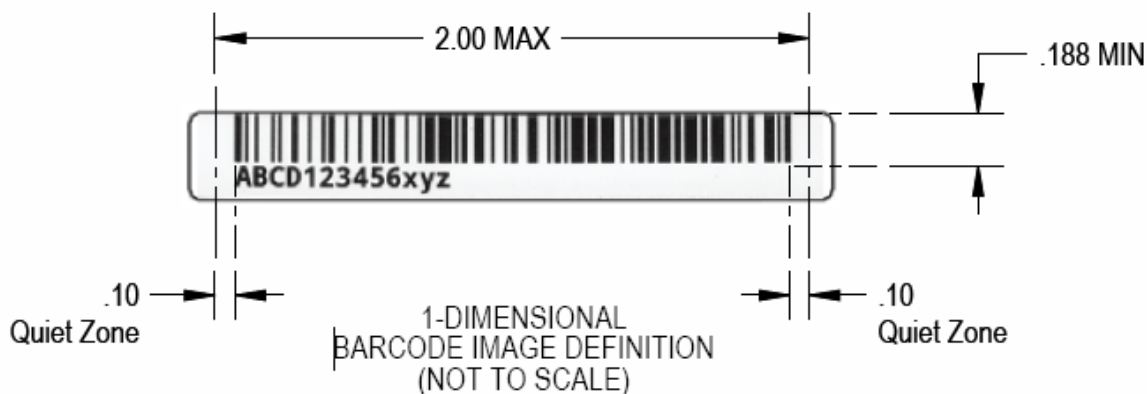

Placing a Barcode Label

Barcode labels can be located on either the long sides of a microplate or the short side opposite the well A1 location. The next figure shows the location on the microplate where the label must be placed in order to be read properly.

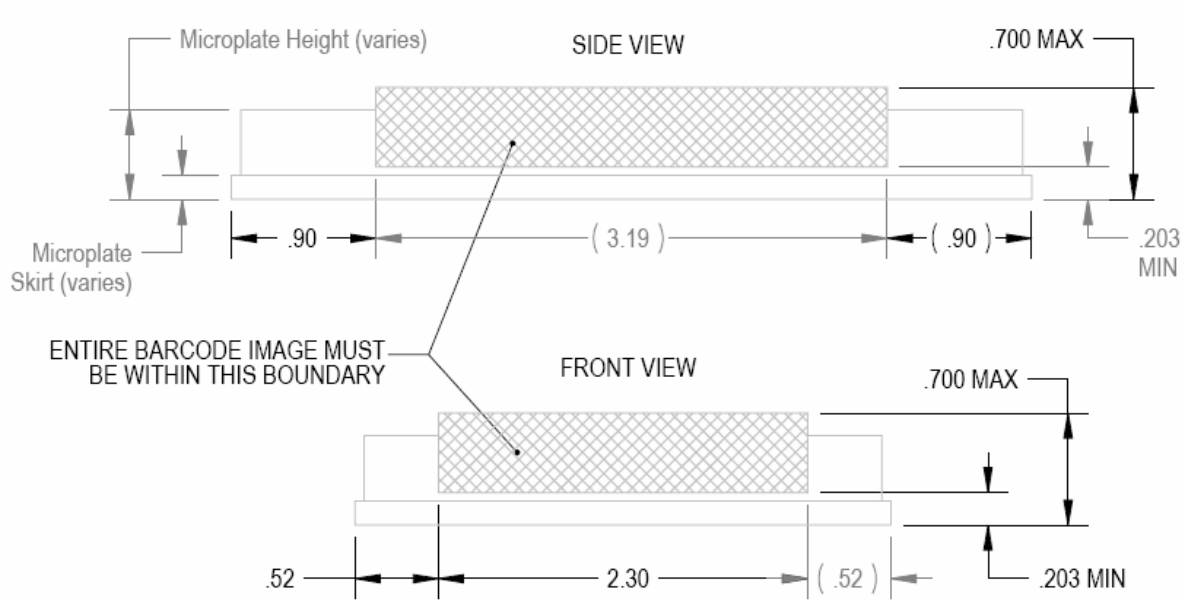


Label Placement Examples


Review the following examples of good label placement:



These examples demonstrate poor placement of the barcode labels:

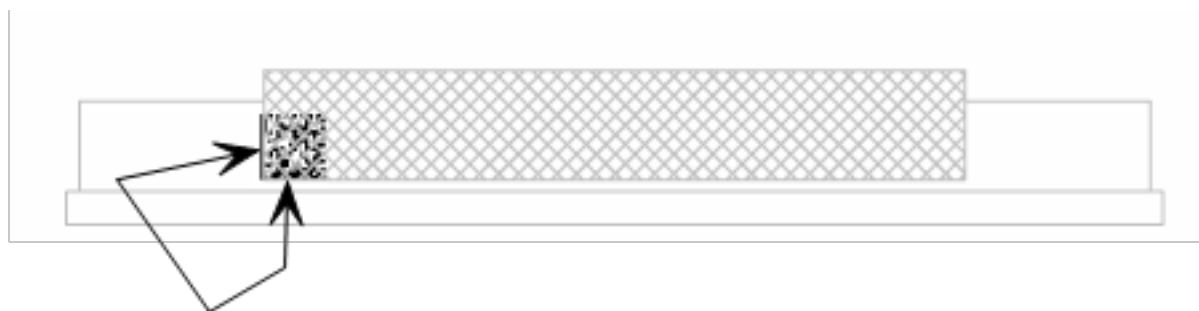

Quiet zone outside of boundary

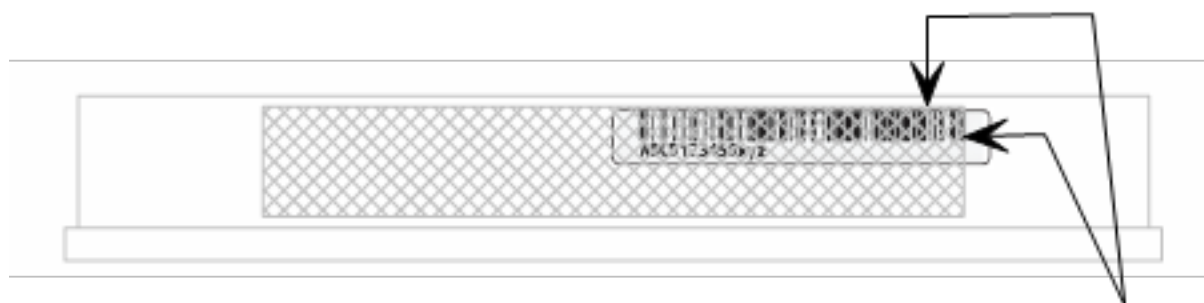


Quiet zone outside of boundary

Microplate Scanner Diagnostic Utility

To ensure that the microplate barcode scanner mirror is correctly positioned to scan the barcodes, use the Barcode Scanner Diagnostic utility in Gen5.

1. Place a microplate with a barcode label on the Synergy Neo 2 carrier.
2. From the main Gen5 screen, click **System > Instrument Configuration**, select the **Synergy Neo 2**, and click **View/Modify > Setup**.
3. On the Barcode Scanner Diagnostic tab, click **Start Diagnostic**. You can see the reflected red light from the laser shining toward the microplate carrier.
4. If necessary, click **Jog Carrier In** to gradually move the carrier so that the barcode label is within the range of the laser light. When the barcode scanner can read the label, the results are displayed in the Reader Response area.

5. If the scanner is unable to read the label, adjust the scanner's mirror position using the mirror adjustment thumbscrew and run the utility again.

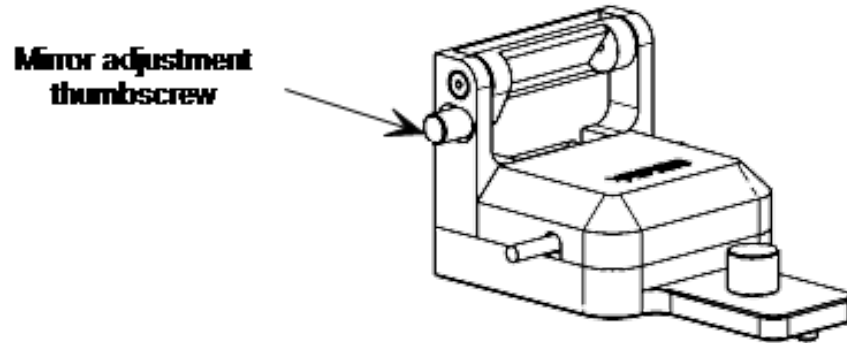

6. When the mirror is set correctly, click **Stop Diagnostic**. The barcode scanner is now ready to use.

Appendix D

Safety Information

Veiligheidsmededelingen

Avis de sécurité

Sicherheitshinweise

Avvisi di sicurezza

Avisos de seguridad

This appendix contains safety information for the 800 TS, translated into Dutch, French, German, Italian, and Spanish.

Safety Notices

Veiligheidsmededelingen

Avis de sécurité

Sicherheitshinweise

Avvisi di sicurezza

Avisos de seguridad

Pay special attention to the following safety notices in all product documentation.

Let vooral op de volgende veiligheidsmededelingen in alle productdocumentatie.

Portez une attention particulière aux avis de sécurité suivants dans l'ensemble de la documentation du produit.

Achten Sie besonders auf die folgenden Sicherheitshinweise in allen Produktdokumentationen.

Prestare particolare attenzione agli avvisi di sicurezza presenti in tutta la documentazione del prodotto.

Preste especial atención a los siguientes avisos de seguridad en toda la documentación del producto.

WARNING

A WARNING notice denotes a hazard. It calls attention to an operating procedure, practice, or the like that, if not correctly performed or adhered to, could result in personal injury or death. Do not proceed beyond a WARNING notice until the indicated conditions are fully understood and met.

De aanduiding WAARSCHUWING duidt op een gevaar. Deze vestigt de aandacht op een bedieningsprocedure, praktijk of iets dergelijks die, indien niet correct uitgevoerd of nageleefd, persoonlijk letsel of de dood tot gevolg kan hebben. Ga niet verder bij een aanduiding WAARSCHUWING voordat de aangegeven voorwaarden volledig begrepen zijn en eraan voldaan is.

Un AVERTISSEMENT signale un danger. Il attire l'attention sur une procédure d'utilisation, une pratique ou autre qui, si elle n'est pas correctement exécutée ou respectée, peut entraîner des dommages corporels, voire un décès. Ne passez pas outre l'AVERTISSEMENT uniquement si les conditions indiquées sont entièrement comprises et remplies.

Ein WARNHINWEIS weist auf eine Gefahr hin. Er weist auf ein Betriebsverfahren, eine Vorgehensweise oder ähnliches hin, deren falsche Ausführung oder Nichtbeachtung zu Verletzungen oder zum Tod führen können. Fahren Sie bei einem WARNHINWEIS erst dann mit Ihrer Arbeit fort, wenn die angegebenen Bedingungen vollständig verstanden und erfüllt sind.

Un avviso di AVVERTENZA indica un pericolo. Richiama l'attenzione su procedure operative, pratiche o azioni simili che, se non rispettate o eseguite correttamente, potrebbero causare lesioni personali o decesso. Non procedere ignorando un avviso di AVVERTENZA fino a quando le condizioni indicate non sono state completamente comprese e soddisfatte.

Un aviso de ADVERTENCIA indica un peligro. Destaca la importancia de un procedimiento operativo, una práctica o un proceso similar que, si no se realiza o se sigue correctamente, podría provocar lesiones o la muerte. No siga adelante sin antes comprender y cumplir plenamente los requisitos indicados en el aviso de ADVERTENCIA.

CAUTION

A CAUTION notice denotes a hazard. It calls attention to an operating procedure, practice, or the like that, if not correctly performed or adhered to, could result in damage to the product or loss of important data. Do not proceed beyond a CAUTION notice until the indicated conditions are fully understood and met.

De aanduiding VOORZICHTIG duidt op een gevaar. Deze vestigt de aandacht op een bedieningsprocedure, praktijk of iets dergelijks die, indien niet correct uitgevoerd of nageleefd, schade aan het product of verlies van belangrijke gegevens tot gevolg kan hebben. Ga niet verder bij een aanduiding VOORZICHTIG voordat de aangegeven voorwaarden volledig begrepen zijn en eraan voldaan is.

Une MISE EN GARDE signale un danger. Elle attire l'attention sur une procédure d'utilisation, une pratique ou autre qui, si elle n'est pas correctement exécutée ou respectée, peut endommager le produit ou entraîner la perte de données importantes. Ne passez pas outre la MISE EN GARDE uniquement si les conditions indiquées sont entièrement comprises et remplies.

Ein VORSICHTSHINWEIS weist auf eine Gefahr hin. Er weist auf ein Betriebsverfahren, eine Vorgehensweise oder ähnliches hin, deren falsche Ausführung oder Nichtbeachtung zu einer Beschädigung des Produkts oder zum Verlust wichtiger Daten führen kann. Fahren Sie bei einem VORSICHTSHINWEIS erst dann mit Ihrer Arbeit fort, wenn die angegebenen Bedingungen vollständig verstanden und erfüllt sind.

Un avviso di ATTENZIONE indica un pericolo. Richiama l'attenzione su procedure operative, pratiche o azioni simili che, se non rispettate o eseguite correttamente, potrebbero causare danni al prodotto o perdita di dati importanti. Non procedere ignorando un avviso di ATTENZIONE fino a quando le condizioni indicate non sono state completamente comprese e soddisfatte.

Un aviso de PRECAUCIÓN indica un peligro. Destaca la importancia de un procedimiento operativo, una práctica o un proceso similar que, si no se realiza o no se sigue correctamente, podrían provocar daños en el producto o la pérdida de datos importantes. No siga adelante sin antes comprender y cumplir plenamente los requisitos indicados en el aviso de PRECAUCIÓN.

Warnings and Precautions

Electrical Hazards

Elektrische gevaren

Risques électriques

Elektrische Gefahren

Rischi elettrici

Peligros eléctricos

WARNING **Internal Voltage.** Always turn off the power switch and unplug the power supply before cleaning the outer surface of the instrument.

Interne spanning. Zet altijd de stroomschakelaar uit en haal de stekker uit het stopcontact voordat de buitenkant van het instrument wordt gereinigd.

Tension interne. Désactivez toujours l'interrupteur d'alimentation électrique et débranchez l'alimentation avant de nettoyer la surface extérieure de l'instrument.

Spannung im Geräteinneren. Vor dem Reinigen der Außenfläche des Geräts grundsätzlich den Stromschalter ausschalten und das Stromkabel aus der Steckdose ziehen.

Tensione interna. Spegner sempre l'interruttore dell'alimentazione e scollegare l'alimentazione prima di pulire le superfici esterne dello strumento.

Tensión interna. Siempre apague el interruptor y desconecte la fuente de alimentación antes de limpiar la superficie exterior del instrumento.

WARNING

Power Rating. The instrument's power supply or power cord must be connected to a power receptacle that provides voltage and current within the specified rating for the system. Use of an incompatible power receptacle may produce electrical shock and fire hazards.

Vermogensklasse. De voeding of het netsnoer van het instrument moet worden aangesloten op een stopcontact dat spanning en stroom levert binnen de gespecificeerde nominale waarden voor het systeem. Gebruik van een niet-compatibel stopcontact kan leiden tot elektrische schokken en brandgevaar.

Puissance électrique nominale. L'alimentation ou le cordon d'alimentation de l'instrument doit être raccordé(e) à une prise de courant qui fournit la tension et le courant correspondants à la puissance spécifiée du système. L'emploi d'une prise de courant incompatible peut entraîner un choc électrique et un risque d'incendie.

Leistungsbemessung. Die Stromversorgung des Geräts bzw. das Anschlusskabel muss mit einer Steckdose verbunden werden, deren Spannungs- und Stromwerte innerhalb der für das System vorgeschriebenen Nennwerte liegen. Die Verwendung einer nicht kompatiblen Steckdose kann zu einem elektrischen Schlag und Brandgefahr führen.

Potenza nominale. L'alimentazione o il cavo di alimentazione dello strumento devono essere collegati a una presa di corrente che fornisca tensione e corrente comprese entro il valore nominale previsto per il sistema. L'uso di una presa di alimentazione non compatibile può causare scosse elettriche e rischi di incendio.

Potencia nominal. La fuente de alimentación o el cable de alimentación del instrumento tienen que conectarse a un receptáculo que suministre tensión y corriente dentro de la potencia especificada para el sistema. El uso de un receptáculo incompatible puede producir descargas eléctricas y riesgo de incendio.

WARNING

Electrical Grounding. Never use a plug adapter to connect primary power to the external power supply. Use of an adapter disconnects the utility ground, creating a severe shock hazard. Always connect the power cord directly to an appropriate receptacle with a functional ground.

Elektrische aarding. Gebruik nooit een stekkeradapter om de primaire stroom aan te sluiten op de externe voeding. Het gebruik van een adapter verbreekt de verbinding met de aarding van het elektriciteitsnet, waardoor een ernstige schok kan ontstaan. Sluit het netsnoer altijd rechtstreeks aan op een geschikt stopcontact met werkende aarding.

Mise à la terre électrique. N'utilisez jamais d'adaptateur de prise pour raccorder l'alimentation principale à l'alimentation électrique extérieure. L'utilisation d'un adaptateur déconnecte la terre du secteur, créant un risque important de choc. Raccordez toujours le cordon d'alimentation directement à une prise appropriée dotée d'une mise à la terre fonctionnelle.

Elektrische Erdung. Verwenden Sie niemals einen Steckeradapter zum Anschließen der Primärstromversorgung an die externe Stromversorgung. Bei Verwendung eines Adapters wird die Verbindung zur Gebäudeerde unterbrochen, sodass ein erhebliches Stromschlagrisiko besteht. Das Stromkabel ist immer direkt an eine geeignete Steckdose mit Funktionserdung anzuschließen.

Messa a terra elettrica. Non usare mai un adattatore per collegare l'alimentazione principale all'alimentazione esterna. Se si usa un adattatore, si scollega la messa a terra della rete elettrica creando un grave pericolo di scosse elettriche. Collegare sempre il cavo di alimentazione direttamente a una presa idonea dotata di messa a terra funzionale.

Conexión a tierra. Nunca use un adaptador de enchufe para conectar la corriente principal a la fuente de alimentación externa. El uso de un adaptador desconecta la tierra del servicio y crea un riesgo de descarga grave. Conecte siempre el cable de alimentación directamente a un receptáculo adecuado con una toma de tierra funcional.

WARNING

Service. Only qualified technical personnel should perform service procedures on internal components.

Service. Alleen gekwalificeerd technisch personeel mag serviceprocedures aan interne onderdelen uitvoeren.

Entretien. L'exécution des procédures d'entretien des composants internes doit être réservée au personnel technique qualifié.

Wartung. Wartungsarbeiten an Komponenten im Geräteinneren sollten nur von qualifizierten Servicetechnikern durchgeführt werden.

Manutenzione. Le procedure di manutenzione sui componenti interni devono essere eseguite esclusivamente da personale tecnico qualificato.

Revisión. Solo puede realizar procedimientos de revisión de los componentes internos el personal técnico cualificado.

CAUTION

Power Supply. Use only the power supply shipped with the instrument, and operate it within the range of line voltages listed on it.

Voeding. Gebruik alleen de voeding die bij het instrument is geleverd en gebruik deze binnen het bereik van de netspanningen die op de voeding staan vermeld.

Alimentation électrique. Utilisez exclusivement l'alimentation électrique fournie avec l'instrument dans la plage de tension de ligne indiquée dessus.

Stromversorgung. Verwenden Sie nur die im Lieferumfang des Geräts enthaltene Stromversorgung und betreiben Sie diese innerhalb des darauf angegebenen Netzspannungsbereichs.

Alimentazione. Usare esclusivamente l'alimentatore fornito con lo strumento, utilizzando quest'ultimo entro l'intervallo delle tensioni di linea indicato sull'unità.

Fuente de alimentación. Use únicamente la fuente de alimentación incluida con el instrumento y úsela en el rango de tensiones de línea indicado en ella.

Chemical/Environmental

Chemisch/Milieu

Substances chimiques/Environnement

Chemie/Umwelt

Rischi chimici/ambientali

Riesgos químicos y medioambientales

WARNING



Potential Biohazards. Some assays or specimens may pose a biohazard. Adequate safety precautions should be taken as outlined in the assay's package insert. Always wear safety glasses and appropriate protective equipment, such as chemical-resistant rubber gloves and apron.

Potentiële biologische gevaren. Sommige tests of specimens kunnen een biologisch gevaar inhouden. Er moeten adequate veiligheidsmaatregelen worden getroffen zoals aangegeven in de bijsluiters van de test. Draag altijd een veiligheidsbril en geschikte beschermingsmiddelen, zoals chemicaliënbestendige rubberen handschoenen en een schort.

Risques biologiques potentiels. Certains tests ou échantillons peuvent présenter un risque biologique. Des précautions de sécurité adéquates doivent être prises, comme indiqué dans la notice de l'emballage du test. Portez toujours des lunettes de sécurité et un équipement de protection approprié, comme des gants en caoutchouc résistant aux substances chimiques et un tablier.

Potenzielle Biogefahren. Manche Assays oder Proben stellen eine Biogefahr dar. Es sollten angemessene Sicherheitsvorkehrungen entsprechend der Packungsbeilage des Assays ergriffen werden. Tragen Sie immer eine Schutzbrille und eine geeignete Schutzausrüstung, wie chemikalienbeständige Gummihandschuhe und Schürze.

Potenziali rischi biologici. Alcuni test o campioni potrebbero comportare un rischio biologico. Implementare misure di sicurezza adeguate secondo quanto delineato nel foglietto della confezione del test. Indossare sempre occhiali di sicurezza e dispositivi di protezione appropriati, ad esempio guanti e grembiule in gomma resistenti alle sostanze chimiche.

Riesgos biológicos potenciales. Algunos ensayos y especímenes pueden constituir un riesgo biológico. Se han de tomar precauciones de seguridad suficientes tal como se indica en el folleto del paquete del ensayo. Use siempre gafas de seguridad y equipos protectores adecuados, como guantes de caucho resistentes a productos químicos y un delantal.

WARNING

Liquids. Avoid spilling liquids on the instrument; fluid seepage into internal components creates a potential for shock hazard or instrument damage. If a spill occurs while a program is running, stop the program and turn off the instrument. Wipe up all spills immediately. Do not operate the instrument if internal components have been exposed to fluid.

Vloeistoffen. Voorkom dat vloeistoffen op het instrument worden gemorst; het doorsijpelen van vloeistoffen in interne onderdelen kan leiden tot schokgevaar of beschadiging van het instrument. Als een lekkage optreedt terwijl een programma loopt, stopt u het programma en schakelt u het instrument uit. Veeg alle gemorste vloeistof onmiddellijk op. Gebruik het instrument niet als interne onderdelen aan vloeistof zijn blootgesteld.

Liquides. Évitez de renverser des liquides sur l'instrument ; les infiltrations de liquide dans les composants internes créent un risque potentiel de choc ou de détérioration de l'instrument. En cas de déversement de liquide alors qu'un programme est en cours d'exécution, arrêtez le programme et mettez l'instrument hors tension. Essayez immédiatement tout liquide renversé. N'utilisez pas l'instrument si les composants internes ont été exposés à du liquide.

Flüssigkeiten. Keine Flüssigkeiten auf dem Gerät verschütten! In die Bauteile im Geräteinneren bilden einsickernde Flüssigkeiten ein Potenzial für die Gefahr von Stromschlägen oder Schäden am Gerät. Bei Verschütten von Flüssigkeiten während ein Programm läuft, ist dieses zu stoppen und das Gerät auszuschalten. Verschüttete Flüssigkeiten sind unverzüglich abzuwischen. Das Gerät darf nicht betrieben werden, wenn Komponenten im Geräteinneren Flüssigkeiten ausgesetzt waren.

Liquidì. Evitare di versare liquidi sullo strumento; l'infiltrazione di fluidi nei componenti interni crea rischi di scosse elettriche o danni allo strumento. Se si verifica un versamento durante l'esecuzione di un programma, arrestare il programma e spegnere lo strumento. Ripulire immediatamente tutti i versamenti. Non utilizzare lo strumento se i componenti interni sono stati esposti a fluidi.

Líquidos. Procure no derramar líquidos sobre el instrumento, ya que si se filtran fluidos en los componentes internos se puede producir un riesgo de descarga o de deterioro del instrumento. Si se produce un derramamiento mientras se está ejecutando un programa, detenga el programa y apague el instrumento. Limpie el derrame inmediatamente. No utilice el instrumento si los componentes internos han estado expuestos a fluidos.

CAUTION

Liquids. Do not immerse the instrument, spray it with liquid, or use a dripping-wet cloth on it. Do not allow water or other cleaning solution to run into the interior of the instrument. If this happens, contact Technical Support.

Vloeistoffen. Dompel het instrument niet onder, bespuit het niet met vloeistof en gebruik er geen druijpnatte doek op. Zorg ervoor dat er geen water of andere schoonmaakmiddelen in het inwendige van het instrument terechtkomen. Als dit gebeurt, neem dan contact op met de afdeling Technische Ondersteuning.

Liquides. N'immergez pas l'instrument, ne le vaporisez pas de liquide et n'utilisez pas de chiffon non essoré dessus. Ne laissez pas d'eau ou autre solution de nettoyage pénétrer à l'intérieur de l'instrument. Le cas échéant, contactez l'assistance technique.

Flüssigkeiten. Das Gerät nicht in Flüssigkeit eintauchen oder damit einsprühen und keine tropfnassen Tücher verwenden. Kein Wasser oder andere Reinigungslösung in das Geräteinnere eindringen lassen. Sollte dies vorkommen, setzen Sie sich mit dem technischen Kundendienst in Verbindung.

Liquidì. Non immergere lo strumento, nebulizzarlo con liquidi né usare un panno che non sia stato strizzato bene. Evitare che acqua o soluzioni detergenti penetrino all'interno dello strumento. Se si verifica un'infiltrazione, contattare il supporto tecnico.

Líquidos. No sumerja el instrumento, no lo pulverice con líquidos y no use un paño mojado que gotee sobre él. No permita que entre agua ni otra solución de limpieza en el interior del instrumento. Si esto sucediera, póngase en contacto con el servicio de soporte técnico.

CAUTION

Environmental Conditions. Do not expose the instrument to temperature extremes. For proper operation, temperature near the instrument should remain within the range in the *Specifications* section of this document. Performance may be adversely affected if temperatures fluctuate above or below this range.

Omgevingsvoorwaarden. Stel het instrument niet bloot aan extreme temperaturen. Voor een goede werking moet de temperatuur in de buurt van het instrument binnen het bereik blijven zoals aangegeven in het gedeelte Specificaties van dit document. De prestaties kunnen nadelig worden beïnvloed als de temperatuur boven of onder dit bereik schommelt.

Conditions environnementales. N'exposez pas l'instrument à des températures extrêmes. Pour assurer un bon fonctionnement, la température à proximité de l'instrument doit demeurer dans la plage indiquée sous la rubrique Spécifications du présent document. La performance peut être affectée négativement si les températures fluctuent au-dessus ou au-dessous de cette plage.

Umgebungsbedingungen. Das Gerät darf keinen Extremtemperaturen ausgesetzt werden. Für den ordnungsgemäßen Betrieb müssen die Temperaturen in Gerätenähe in den im Abschnitt Spezifikationen dieses Dokuments angegebenen Grenzen bleiben. Temperaturschwankungen über diese Grenzwerte hinaus können die Geräteleistung beeinträchtigen.

Condizioni ambientali. Non esporre lo strumento a temperature estreme. Per il corretto funzionamento, la temperatura nei pressi dello strumento deve restare nell'intervallo indicato nella sezione Specifiche di questo documento. Fluttuazioni delle temperature al di sopra o al di sotto di questo intervallo possono compromettere le prestazioni dello strumento.

Condiciones ambientales. No exponga el instrumento a temperaturas extremas. Para su correcto funcionamiento, la temperatura que rodee al instrumento deberá estar dentro del rango indicado en la sección Especificaciones de este documento. Si las temperaturas fluctúan por encima o por debajo de este rango, el rendimiento puede verse afectado negativamente.

CAUTION

Sodium Hypochlorite. Do not expose any part of the instrument to the recommended diluted sodium hypochlorite solution for more than 20 minutes. Prolonged contact may damage the instrument surfaces. Be certain to rinse and thoroughly wipe all surfaces.

Natriumhypochloriet. Stel geen enkel deel van het instrument langer dan 20 minuten bloot aan de aanbevolen verdunde natriumhypochlorietoplossing. Langdurig contact kan de oppervlakken van het instrument beschadigen. Zorg ervoor dat alle oppervlakken goed worden afgespoeld en schoongeveegd.

Hypochlorite de sodium. N'exposez aucune pièce de l'instrument à la solution d'hypochlorite de sodium diluée comme recommandé pendant plus de 20 minutes. Un contact prolongé peut endommager les surfaces de l'instrument. Veillez à rincer et essuyer soigneusement toutes les surfaces.

Natriumhypochlorit. Kein Teil des Geräts darf der empfohlenen verdünnten Natriumhypochloritlösung länger als 20 Minuten lang ausgesetzt werden. Bei längerem Kontakt drohen Beschädigungen an den Geräteoberflächen. Alle Oberflächen unbedingt abspülen und gründlich abwischen.

Ipoclorito di sodio. Non esporre nessun componente dello strumento alla soluzione di ipoclorito di sodio diluita raccomandata per più di 20 minuti. Un contatto prolungato potrebbe danneggiare le superfici dello strumento. Accertarsi di sciacquare e ripulire accuratamente tutte le superfici.

Hipoclorito sódico. No exponga ninguna parte del instrumento a la solución de hipoclorito sódico diluido recomendada durante más de 20 minutos. Un contacto demasiado prolongado puede dañar las superficies del instrumento. Asegúrese de aclarar y secar concienzudamente todas las superficies.

CAUTION

Lubricants. Do not apply lubricants to moving parts. Lubricant on components in the carrier compartment will attract dust and other particles, which may cause the instrument to produce an error.

Smeermiddelen. Breng geen smeermiddelen aan op bewegende delen. Smeermiddel op onderdelen in het draagcompartiment zal stof en andere deeltjes aantrekken, waardoor het instrument een fout kan produceren.

Lubrifiants. N'appliquez pas de lubrifiants sur les pièces mobiles. La présence de lubrifiant sur les composants dans le compartiment du portoir attire la poussière et autres particules, ce qui peut provoquer une erreur de l'instrument.

Schmierstoffe. Keine Schmierstoffe auf bewegliche Teile auftragen. Schmierstoffe auf Komponenten im Trägerfach ziehen Staub und andere Teilchen an, die zu einem Gerätefehler führen können.

Lubrificanti. Non applicare lubrificanti alle parti in movimento. La presenza di lubrificante sui componenti del vano portapietra attira polvere e altre particelle che potrebbero causare errori dello strumento.

Lubricantes. No aplique lubricantes en las piezas móviles. El lubricante en los componentes del compartimento del portador atraerá polvo y otras partículas que pueden hacer que el instrumento muestre un error.

Components

Onderdelen

Composants

Komponenten

Componenti

Componentes

WARNING



Laser Beam. Serious eye injury may occur if you stare directly into the laser beam of the barcode scanner during operation of the scanner. This hazard is noted by the symbol shown here. Do not look directly into the laser beam during operation of the scanner.

Laserstraal. Er kan ernstig oogletsel ontstaan als u tijdens het gebruik van de scanner rechtstreeks in de laserstraal van de barcodescanner kijkt. Dit gevaar wordt aangegeven met het hier afgebeelde symbool. Kijk tijdens het gebruik van de scanner niet rechtstreeks in de laserstraal.

Faisceau laser. De graves lésions oculaires peuvent survenir si vous regardez directement dans le faisceau laser du lecteur de code-barres pendant le fonctionnement du scanner. Ce danger est signalé par le symbole illustré ci-contre. Ne regardez pas directement dans le faisceau laser pendant le fonctionnement du scanner.

Laserstrahl. Schwere Augenverletzungen können auftreten, wenn Sie während des Betriebs des Scanners direkt in den Laserstrahl des Barcode-Scanners blicken. Diese Gefahr wird durch das hier abgebildete Symbol angezeigt. Während des Betriebs des Scanners nicht direkt in den Laserstrahl blicken.

Raggio laser. Per evitare gravi lesioni oculari, non fissare direttamente il raggio laser dello scanner di codici a barre durante il funzionamento dello scanner. Questo pericolo è segnalato dal simbolo mostrato qui. Non fissare direttamente il raggio laser durante il funzionamento dello scanner.

Haz del láser. Mirar directamente al haz del láser del escáner de código de barras cuando el escáner está en funcionamiento puede provocar lesiones oculares graves. El símbolo que se muestra aquí advierte sobre este peligro. No mire directamente al haz del láser cuando el escáner esté funcionando.

WARNING



Class 1 Laser Product

Class 1 Laser Product. Alpha laser "A" and TRF laser "T" models.

See "Report on laser safety" on page 208.

Klasse 1 laserproduct. Alpha Laser "A" en TRF laser "T" - modellen.

Zie "Rapport over laserveiligheid" op pagina 208.

Laser de classe 1. Alpha Laser "A" et laser TRF modèles "T".

Voir « Rapport sur la sécurité laser » à la page 209.

Klasse 1 Laserprodukt. Alpha Laser "A" und TRF Laser "T" - Modellen.

Siehe "Bericht zur Lasersicherheit" auf Seite 210.

Classe di prodotti laser 1. Alpha Laser "A" e laser modelli "T" della Fondazione.

Vedere "Rapporto sulla sicurezza laser" a pagina 211.

Producto láser de clase 1. Modelos de láser alfa "A" y láser TRF "T".

Ver "Informe sobre seguridad láser" en la página 212.

WARNING

Pinch Hazard. Some areas of the external dispense module can present pinch hazards when the instrument is operating. Keep hands and fingers clear of these areas when the instrument is operating.

Beknellingsgevaar. Sommige delen van de externe uitgiftemodule kunnen beknellingsgevaar opleveren wanneer het instrument in bedrijf is. Houd handen en vingers uit de buurt van deze gebieden wanneer het instrument in bedrijf is.

Risque de pincement. Certaines zones du module de dispense externe peuvent présenter des risques de pincement lors du fonctionnement de l'instrument. Gardez vos mains et vos doigts à l'écart de ces zones lors du fonctionnement de l'instrument.

Quetschgefahr. In einigen Bereichen des externen Dispenser-Moduls können beim Betrieb des Geräts Quetschgefahren auftreten. Hände und Finger von diesen Bereichen fernhalten, wenn das Gerät in Betrieb ist.

Rischio di pizzicamento. Alcune aree del modulo di erogazione esterno possono presentare rischi di pizzicamento quando lo strumento è in funzione. Tenere le mani e le dita lontane da queste aree quando lo strumento è in funzione.

Peligro de atrapamiento. Algunas áreas del módulo dispensador externo pueden presentar riesgos de atrapamiento cuando el instrumento está en funcionamiento. Mantenga las manos y los dedos alejados de estas áreas cuando el instrumento esté en funcionamiento.

WARNING

Hot Surface. The lamp assembly is hot when the instrument is turned on. Turn off the reader and allow the bulb to cool for at least 15 minutes before attempting to replace it.

Heet oppervlak. De lamp is heet wanneer het instrument wordt ingeschakeld. Zet het leesapparaat uit en laat de lamp ten minste 15 minuten afkoelen alvorens te proberen deze te vervangen.

Surface chaude. La lampe est chaude lorsque l'instrument est allumé. Éteignez le lecteur et laissez l'ampoule refroidir pendant 15 minutes au moins avant de la remplacer.

Heiße Oberflächen. Bei eingeschaltetem Gerät ist die Lampe heiß. Schalten Sie den Reader aus und lassen Sie die Glühbirne mindestens 15 Minuten lang abkühlen, bevor Sie versuchen, sie auszutauschen.

Superficie molto calda. Il gruppo lampada diventa molto caldo quando lo strumento è acceso. Prima di tentare di sostituirlo, spegnere il lettore e lasciare raffreddare la lampadina per almeno 15 minuti.

Superficie caliente. El conjunto de piezas de la lámpara está caliente cuando el instrumento está encendido. Apague el lector y deje que la bombilla se enfríe durante al menos 15 minutos antes de proceder a cambiarla.

WARNING

Two-person lift. The instrument should be lifted by two people. The instrument with all available modules weighs up to 45.5.kg.

Tillen door twee personen. Het instrument moet door twee personen worden opgetild. Het instrument met alle beschikbare modules weegt maximaal 45,5 kg.

Charge à soulever par deux personnes. L'instrument doit être soulevé par deux personnes. L'instrument avec tous les modules disponibles pèse jusqu'à 345,5 kg

Anheben durch zwei Personen. Das Gerät sollte von zwei Personen angehoben werden. Das Gerät mit allen verfügbaren Modulen wiegt bis zu 45,5 kg.

Due persone per il sollevamento. Lo strumento deve essere sollevato da due persone. Lo strumento con tutti i moduli disponibili pesa fino a 45,5 kg.

Levantamiento por dos personas. Es necesario que dos personas levanten el instrumento. El instrumento con todos los módulos disponibles puede pesar hasta 45,5 kg.

WARNING

Accessories. Only accessories that meet the manufacturer's specifications shall be used with the instrument.

Accessoires. Bij het instrument mogen alleen accessoires worden gebruikt die voldoen aan de specificaties van de fabrikant.

Accessoires. L'instrument doit être utilisé exclusivement avec des accessoires correspondant aux spécifications du fabricant.

Zubehör. In Verbindung mit dem Gerät dürfen nur Zubehörkomponenten verwendet werden, die den Spezifikationen des Herstellers entsprechen.

Accessori. Utilizzare esclusivamente accessori dello strumento che rispettano le specifiche del fabbricante.

Accesorios. Solamente aquellos accesorios que cumplan las especificaciones del fabricante deberán usarse con el instrumento.

CAUTION

Shipping Hardware. All shipping hardware must be removed before operating the instrument and reinstalled before repackaging the instrument for shipment.

Verzendingshardware. Alle verzendingshardware moet worden verwijderd voordat het instrument wordt gebruikt en opnieuw worden geïnstalleerd voordat het instrument opnieuw wordt verpakt voor verzending.

Matériel d'expédition. Tout le matériel d'expédition doit être retiré avant d'utiliser l'instrument et réinstallé avant de remballer l'équipement pour expédition.

Festes Versandmaterial. Alle festen Versandmaterialien müssen vor der Inbetriebnahme des Geräts entfernt und vor der Wiederverpackung des Geräts zum Versand neu angebracht werden.

Minuteria di spedizione. Prima di utilizzare lo strumento, rimuovere tutta la minuteria di spedizione, che dovrà essere reinstallata prima di reimballare lo strumento per la spedizione.

Equipo de envío. Antes de utilizar el instrumento es necesario retirar todo el equipo de envío y, del mismo modo, habrá que volver a colocárselo cuando el instrumento se vaya a enviar.

CAUTION

Spare Parts. Only approved spare parts should be used for maintenance. The use of unapproved spare parts and accessories may result in a loss of warranty and potentially impair instrument performance or cause damage to the instrument.

Reserveonderdelen. Voor onderhoud mogen alleen goedgekeurde reserveonderdelen worden gebruikt. Het gebruik van niet-goedgekeurde onderdelen en accessoires kan tot gevolg hebben dat de garantie vervalt en mogelijk de prestaties van het instrument nadelig beïnvloeden of het instrument beschadigen.

Pièces de rechange. Utilisez exclusivement des pièces de rechange approuvées pour l'entretien. L'utilisation de pièces de rechange et accessoires non approuvés peut entraîner l'annulation de la garantie et potentiellement nuire à la performance de l'instrument ou l'endommager.

Ersatzteile. Für die Wartung sollten nur genehmigte Ersatzteile verwendet werden. Die Verwendung nicht genehmigter Ersatzteile und Zubehörkomponenten kann zum Verlust der Garantie führen und möglicherweise die Geräteleistung beeinträchtigen oder Schäden am Gerät verursachen.

Parti di ricambio. Per la manutenzione, usare esclusivamente parti di ricambio approvate. L'uso di parti di ricambio e accessori non approvati potrebbe dare luogo all'annullamento della garanzia e ripercuotersi negativamente sulle prestazioni o causare danni allo strumento.

Repuestos. Durante el mantenimiento, solo deben emplearse repuestos originales. El uso de repuestos y accesorios no autorizados puede producir la pérdida de la garantía y reducir el funcionamiento del instrumento o provocar daños en él.

CAUTION

Service. Only qualified technical personnel should perform service procedures on internal components.

Service. Alleen gekwalificeerd technisch personeel mag serviceprocedures aan interne onderdelen uitvoeren.

Entretien. L'exécution des procédures d'entretien des composants internes doit être réservée au personnel technique qualifié.

Wartung. Wartungsarbeiten an Komponenten im Geräteinneren sollten nur von qualifizierten Servicetechnikern durchgeführt werden.

Manutenzione. Le procedure di manutenzione sui componenti interni devono essere eseguite esclusivamente da personale tecnico qualificato.

Revisión. Solo puede realizar procedimientos de revisión de los componentes internos el personal técnico cualificado.

Intended Product Use

Beoogd productgebruik

Utilisation prévue du produit

Vorgesehene Produktverwendung

Uso previsto del prodotto

Uso previsto del producto

WARNING **Software Quality Control.** The operator must follow the manufacturer's assay package insert when modifying software parameters and establishing reading methods. It is considered good laboratory practice to run laboratory samples according to instructions and specific recommendations included in the assay package insert for the test to be conducted. Failure to conduct quality control checks could result in erroneous test data.

Softwarekwaliteitscontrole. Bij het wijzigen van de softwareparameters en het vaststellen van afleesmethoden moet de operator de bijsluiter van de test van de fabrikant volgen. Het wordt beschouwd als een goede laboratoriumpraktijk om laboratoriummonsters te onderzoeken volgens de instructies en specifieke aanbevelingen die zijn opgenomen in de bijsluiter van de verpakking van de uit te voeren test. Het niet uitvoeren van kwaliteitscontroles kan leiden tot foutieve testgegevens.

Contrôle de qualité du logiciel. L'opérateur doit respecter la notice présente dans l'emballage du test lorsqu'il modifie les paramètres du logiciel et établit les méthodes de lecture. L'exécution d'échantillons de laboratoire conformément aux instructions et aux recommandations spécifiques présentées dans la notice de l'emballage du test à réaliser est considérée comme une bonne pratique de laboratoire. Ne pas exécuter les vérifications de contrôle de qualité peut produire des données de test erronées.

Report on laser safety

Synergy Neo2 T-models with the N2 laser

Synergy Neo2 "T" models include a nitrogen laser MNL-100 manufactured by LTB. This laser is classified to class 3B and therefore the safety aspects are made according to laser class 3B.

Synergy Neo2 A-models with the Alpha Laser

Synergy Neo2 "A" models include an Alpha Laser, K6880E09FN-0.800W manufactured by BWT Beijing. This laser is classified to class 3B and therefore the safety aspects are made according to laser class.

Concerning class 3B lasers, the International Standard IEC 60825-1 on the safety of laser products states the following: "Class 3B: Lasers are normally hazardous when direct intrabeam exposure occurs (i.e., within the NOHD). Viewing diffuse reflections is normally safe."

Regarding laser safety, it is extremely important to prevent the user from being exposed to laser beam, either directly or through reflection from a mirror surface. Although the laser used in Synergy Neo2 is a class-3B laser, the Synergy Neo2 instrument is a class 1 laser product. This is possible, because Synergy Neo2 has adequate laser shielding which prevents the user being exposed to the laser radiation.

The Synergy Neo2's interlock system prevents user exposure to the laser beam by disabling the laser whenever the top-filter door, plate carrier, or side (bottom-filter) hatch of the instrument is open. The state of the interlock is made visible to users via LEDs on the laser unit attached to the rear of the instrument. Red LEDs indicate interlock activity; green LEDs indicate safe laser usage.

When closed the three doors form a light-tight system. Light tightness is essential not only for laser safety, but also for the functional performance of the instrument. The laser itself has been installed on the instrument under a protective casing, and the laser beam is led by an optical fiber via the electronics compartment to the light-tight instrument compartment.

Rapport over laserveiligheid

Synergy Neo2 T-modellen met de N2-laser

Synergy Neo2 "T"-modellen bevatten een stikstoflaser MNL-100 vervaardigd door LTB. Deze laser is geclassificeerd in klasse 3B en daarom zijn de veiligheidsaspecten gemaakt volgens laserklasse 3B.

Synergy Neo2 A-modellen met de Alpha Laser

Synergy Neo2 "A"-modellen bevatten een Alpha Laser, K6880E09FN-0.800W, vervaardigd door BWT Beijing. Deze laser is geclassificeerd in klasse 3B en daarom zijn de veiligheidsaspecten gemaakt volgens laserklasse.

Met betrekking tot lasers van klasse 3B stelt de internationale norm IEC 60825-1 over de veiligheid van laserproducten het volgende: "Klasse 3B: lasers zijn normaal gesproken gevaarlijk wanneer directe blootstelling binnen de bundel plaatsvindt (d.w.z. binnen de NOHD). Het bekijken van diffuse reflecties is normaal gesproken veilig."

Met betrekking tot laserveiligheid is het uiterst belangrijk om te voorkomen dat de gebruiker wordt blootgesteld aan laserstralen, hetzij direct, hetzij door reflectie van een spiegeland oppervlak. Hoewel de laser die in Synergy Neo2 wordt gebruikt een laser van klasse 3B is, is het Synergy Neo2-instrument een laserproduct van klasse 1. Dit is mogelijk omdat Synergy Neo2 een adequate laserafscherming heeft die voorkomt dat de gebruiker wordt blootgesteld aan de laserstraling.

Het interlocksysteem van de Synergy Neo2 voorkomt blootstelling van de gebruiker aan de laserstraal door de laser uit te schakelen wanneer de bovenste filterdeur, plaatdrager of zijklep (onderfilter) van het instrument open is. De status van de vergrendeling wordt voor gebruikers zichtbaar gemaakt via LED's op de lasereenheid die aan de achterkant van het instrument is bevestigd. Rode LED's geven interlock-activiteit aan; groene LED's geven veilig lasergebruik aan.

In gesloten toestand vormen de drie deuren een lichtdicht systeem. Lichtdichtheid is niet alleen essentieel voor de laserveiligheid, maar ook voor de functionele prestaties van het instrument. De laser zelf is onder een beschermende behuizing op het instrument geïnstalleerd en de laserstraal wordt door een optische vezel via het elektronikacompartiment naar het lichtdichte instrumentencompartiment geleid.

Rapport sur la sécurité laser

Modèles Synergy Neo2 T avec le laser N2

Les modèles Synergy Neo2 « T » comprennent un laser à azote MNL-100 fabriqué par LTB. Ce laser est classé dans la classe 3B et donc les aspects de sécurité sont réalisés selon la classe laser 3B.

Modèles Synergy Neo2 A avec le laser Alpha

Les modèles Synergy Neo2 "A" incluent un laser Alpha, K6880E09FN-0.800W fabriqué par BWT Beijing. Ce laser est classé dans la classe 3B et donc les aspects de sécurité sont faits selon la classe laser.

Concernant les lasers de classe 3B, la norme internationale IEC 60825-1 sur la sécurité des produits laser stipule ce qui suit : « Classe 3B : les lasers sont normalement dangereux lorsqu'une exposition directe à l'intérieur du faisceau se produit (c'est-à-dire dans le NOHD). La visualisation des réflexions diffuses est normalement sans danger.

En ce qui concerne la sécurité laser, il est extrêmement important d'empêcher l'utilisateur d'être exposé au faisceau laser, soit directement, soit par réflexion sur une surface de miroir. Bien que le laser utilisé dans Synergy Neo2 soit un laser de classe 3B, l'instrument Synergy Neo2 est un produit laser de classe 1. Ceci est possible, car Synergy Neo2 dispose d'un blindage laser adéquat qui empêche l'utilisateur d'être exposé au rayonnement laser.

Le système de verrouillage du Synergy Neo2 empêche l'exposition de l'utilisateur au faisceau laser en désactivant le laser chaque fois que la porte du filtre supérieur, le support de plaque ou la trappe latérale (filtre inférieur) de l'instrument est ouvert. L'état du verrouillage est rendu visible aux utilisateurs via des LED sur l'unité laser fixée à l'arrière de l'instrument. Les LED rouges indiquent l'activité de verrouillage ; les LED vertes indiquent une utilisation sûre du laser.

Lorsqu'elles sont fermées, les trois portes forment un système étanche à la lumière. L'étanchéité à la lumière est essentielle non seulement pour la sécurité laser, mais aussi pour les performances fonctionnelles de l'instrument. Le laser lui-même a été installé sur l'instrument sous un boîtier de protection, et le faisceau laser est conduit par une fibre optique via le compartiment électronique vers le compartiment de l'instrument étanche à la lumière.

Bericht zur Lasersicherheit

Synergy Neo2 T-Modelle mit dem N2-Laser

Die Synergy Neo2 „T“-Modelle enthalten einen Stickstofflaser MNL-100, hergestellt von LTB. Dieser Laser ist in die Klasse 3B eingestuft und daher werden die Sicherheitsaspekte nach der Laserklasse 3B ausgeführt.

Synergy Neo2 A-Modelle mit dem Alpha Laser

Die Synergy Neo2 „A“-Modelle umfassen einen Alpha-Laser, K6880E09FN-0.800W, hergestellt von BWT Beijing. Dieser Laser ist in die Klasse 3B eingestuft und daher werden die Sicherheitsaspekte entsprechend der Laserklasse vorgenommen.

In Bezug auf Laser der Klasse 3B heißt es in der Internationalen Norm IEC 60825-1 zur Sicherheit von Laserprodukten: „Klasse 3B: Laser sind normalerweise gefährlich, wenn eine direkte Strahlenexposition (d. h. innerhalb des NOHD) auftritt. Das Betrachten von diffusen Reflexionen ist normalerweise sicher.“

In Bezug auf die Lasersicherheit ist es äußerst wichtig zu verhindern, dass der Benutzer direkt oder durch Reflexion von einer Spiegeloberfläche dem Laserstrahl ausgesetzt wird. Obwohl der in Synergy Neo2 verwendete Laser ein Laser der Klasse 3B ist, ist das Synergy Neo2-Gerät ein Laserprodukt der Klasse 1. Dies ist möglich, da Synergy Neo2 über eine ausreichende Laserabschirmung verfügt, die verhindert, dass der Benutzer der Laserstrahlung ausgesetzt wird.

Das Interlock-System des Synergy Neo2 verhindert, dass der Benutzer dem Laserstrahl ausgesetzt wird, indem der Laser deaktiviert wird, wenn die obere Filtertür, der Plattenträger oder die seitliche (untere Filter-) Klappe des Instruments geöffnet ist. Der

Zustand der Verriegelung wird dem Benutzer über LEDs an der an der Geräterückseite angebrachten Lasereinheit sichtbar gemacht. Rote LEDs zeigen Interlock-Aktivität an; grüne LEDs zeigen eine sichere Lasernutzung an.

Im geschlossenen Zustand bilden die drei Türen ein lichtdichtes System. Lichtdichtheit ist nicht nur für die Lasersicherheit, sondern auch für die Funktionstüchtigkeit des Instruments unerlässlich. Der Laser selbst ist unter einem Schutzgehäuse am Gerät installiert und der Laserstrahl wird über einen Lichtwellenleiter über den Elektronikraum in den lichtdichten Geräteraum geführt.

Rapporto sulla sicurezza laser

Synergy Neo2 T-modelli con il laser N2

I modelli Synergy Neo2 "T" includono un laser ad azoto MNL-100 prodotto da LTB. Questo laser è classificato alla classe 3B e quindi gli aspetti di sicurezza sono realizzati secondo la classe laser 3B.

Synergy Neo2 modelli A con Alpha Laser

I modelli Synergy Neo2 "A" includono un laser Alpha, K6880E09FN-0.800W prodotto da BWT Beijing. Questo laser è classificato in classe 3B e quindi gli aspetti di sicurezza sono realizzati in base alla classe laser.

Per quanto riguarda i laser di classe 3B, lo standard internazionale IEC 60825-1 sulla sicurezza dei prodotti laser afferma quanto segue: "Classe 3B: i laser sono normalmente pericolosi quando si verifica un'esposizione diretta all'interno del raggio (cioè, all'interno del NOHD). La visualizzazione di riflessi diffusi è normalmente sicura."

Per quanto riguarda la sicurezza del laser, è estremamente importante evitare che l'utente venga esposto al raggio laser, direttamente o attraverso la riflessione da una superficie a specchio. Sebbene il laser utilizzato in Synergy Neo2 sia un laser di classe 3B, lo strumento Synergy Neo2 è un prodotto laser di classe 1. Ciò è possibile, poiché Synergy Neo2 dispone di un'adeguata schermatura laser che impedisce all'utente di essere esposto alla radiazione laser.

Il sistema di interblocco di Synergy Neo2 impedisce l'esposizione dell'utente al raggio laser disabilitando il laser ogni volta che lo sportello del filtro superiore, il portapiastre o il portello laterale (filtro inferiore) dello strumento sono aperti. Lo stato dell'interblocco è reso visibile agli utenti tramite LED sull'unità laser fissata sul retro dello strumento. I LED rossi indicano l'attività dell'interblocco; i LED verdi indicano un utilizzo sicuro del laser.

Quando sono chiuse le tre ante formano un sistema a tenuta di luce. La tenuta alla luce è essenziale non solo per la sicurezza del laser, ma anche per le prestazioni funzionali dello strumento. Il laser stesso è stato installato sullo strumento sotto un involucro protettivo e il raggio laser è guidato da una fibra ottica attraverso il vano dell'elettronica al vano dello strumento a tenuta di luce.

Informe sobre seguridad láser

Modelos Synergy Neo2 T con el láser N2

Los modelos Synergy Neo2 "T" incluyen un láser de nitrógeno MNL-100 fabricado por LTB. Este láser está clasificado en la clase 3B y, por lo tanto, los aspectos de seguridad se realizan de acuerdo con la clase de láser 3B.

Modelos Synergy Neo2 A con Alpha Laser

Los modelos Synergy Neo2 "A" incluyen un Alpha Laser, K6880E09FN-0.800W fabricado por BWT Beijing. Este láser está clasificado en la clase 3B y, por lo tanto, los aspectos de seguridad se realizan de acuerdo con la clase de láser.

Con respecto a los láseres de clase 3B, la norma internacional IEC 60825-1 sobre la seguridad de los productos láser establece lo siguiente: "Clase 3B: los láseres son normalmente peligrosos cuando se produce una exposición directa al haz (es decir, dentro del NOHD). Normalmente, ver reflejos difusos es seguro".

Con respecto a la seguridad del láser, es extremadamente importante evitar que el usuario se exponga al rayo láser, ya sea directamente o mediante el reflejo de la superficie de un espejo. Aunque el láser utilizado en Synergy Neo2 es un láser de clase 3B, el instrumento Synergy Neo2 es un producto láser de clase 1. Esto es posible, porque Synergy Neo2 tiene un blindaje láser adecuado que evita que el usuario se exponga a la radiación láser.

El sistema de enclavamiento del Synergy Neo2 evita la exposición del usuario al rayo láser al desactivar el láser siempre que la puerta del filtro superior, el portador de la placa o la trampilla lateral (filtro inferior) del instrumento estén abiertos. El estado del enclavamiento se hace visible para los usuarios a través de LED en la unidad láser adjunta a la parte posterior del instrumento. Los LED rojos indican actividad de interbloqueo; Los LED verdes indican un uso seguro del láser.

Cuando están cerradas, las tres puertas forman un sistema hermético a la luz. La estanqueidad a la luz es esencial no solo para la seguridad del láser, sino también para el rendimiento funcional del instrumento. El láser en sí se ha instalado en el instrumento debajo de una carcasa protectora, y el rayo láser es conducido por una fibra óptica a través del compartimiento de la electrónica al compartimiento del instrumento a prueba de luz.

In This Book

This document contains installation, operation, maintenance, and qualification information for all models of the Synergy Neo 2.

Document Revision History			
Part Number	Revision	Date	Modifications
1351000N	A	December 2021	First issue.

Original Language – EN



Manufactured by Agilent Technologies, Inc.

5301 Stevens Creek Blvd.

Santa Clara, CA 95051

